



OCCASION

This publication has been made available to the public on the occasion of the 50th anniversary of the United Nations Industrial Development Organisation.



DISCLAIMER

This document has been produced without formal United Nations editing. The designations employed and the presentation of the material in this document do not imply the expression of any opinion whatsoever on the part of the Secretariat of the United Nations Industrial Development Organization (UNIDO) concerning the legal status of any country, territory, city or area or of its authorities, or concerning the delimitation of its frontiers or boundaries, or its economic system or degree of development. Designations such as "developed", "industrialized" and "developing" are intended for statistical convenience and do not necessarily express a judgment about the stage reached by a particular country or area in the development process. Mention of firm names or commercial products does not constitute an endorsement by UNIDO.

FAIR USE POLICY

Any part of this publication may be quoted and referenced for educational and research purposes without additional permission from UNIDO. However, those who make use of quoting and referencing this publication are requested to follow the Fair Use Policy of giving due credit to UNIDO.

CONTACT

Please contact <u>publications@unido.org</u> for further information concerning UNIDO publications.

For more information about UNIDO, please visit us at www.unido.org

2/990

A Handbook of Chinese Pharmaceutical Machinery & Equipment for Export

中国制药装备出口产品手册

MAY.1998 一九九八年五月

A Handbook of Chinese Pharmaceutical Machinery & Equipment for Export

中国制药装备出口产品手册

Department of International Cooperation,
The State Pharmaceutical Administration of China

国家医药管理局 国际合作司

The China Association for Pharmaceutical Equipment

中国制药装备行业协会

Copyright ©1998
Copyright in the material included in this book is reserved by the individual organizations

版权所有,违者必究

Preface

With the support of UNIDO,"A Handbook of Chinese Pharmaceutical Machinery & Equipment for Export " Compiled by The State Pharmaceutical Administration of China is now published. I am very pleased to extend my warm congratulations on the coming out of this book.

After the founding of new China, particularly since China's economic reform and opening up to the outside world, the Chinese pharmaceutical industry has witnessed tremendous development. Its industry output value increased at an average annual rate of 21% from 1990 to 1995. In 1996, the output value amounted to RMB \(\frac{1}{2}\) 120 billion, 15\% more than that of 1995. In 1997 the industry reached an output value of RMB ¥ 140 billion, while the commercial sales value was amounted to RMB ¥ 102 billion. With the development of pharmaceutical industry, the pharmaceutical machinery sector has grown rapidly. At pesent, there are around 400 pharmaceutical machinery enterprises. A pharmaceutical machinery industrial system consisting of scientific research, design, production and distribution has been formed with its products more than 1100 varieties and specifications covering eight large categories - machinery for bulk pharmaceuticals and preparations, pharmaceutical machinery, pulverizing machinery, hydraulic equipment, packaging machinery, drug testing instruments and other related pharmaceutical machinery.

"A Handbook of Chinese Pharmaceutical Machinery & Equipment for Export" introduces the products, specifications and the models of the main pharmaceutical machinery enterprises in China with more than 85 pictures. The publishing of this book will make contribution to understanding of Chinese pharmaceutical machinery and enterprises, and to promoting cooperation among developing countries in the field of pharmaceutical industry.

Zheng Xiaoyu

Director-General

The State Pharmaceutical Administration of China

前言

在 UNIDO 的支持下,由国家医药管理局编辑的《中国制药装备出口产品手册》现已出版发行。我对该书的出版表示热烈的祝贺。

建国之后,特别是改革开放以来,我国医药工业发展迅速。1990 到 1995 年,医药工业总产值以年均 21%的速度递增。1996 年医药工业总产值 1200 亿元,比 1995 年递增 15%。1997 年医药工业总产值 1400 亿元,医药商业销售额 1021 亿元。随着医药工业的发展,制药机械工业也得到了迅速发展。目前,有近 400 家制药机械企业,形成了具有科研、设计、生产及销售组成的制药机械工业体系,可生产原料药设备及机械、制剂机械、药用分碎机械、饮片机械、制药用水设备、药品包装机械、药物检测设备及其他制药机械及设备等八大类、1100 余个品种、规格的制药机械产品。

《中国制药装备出口产品手册》介绍了我国主要制药设备企业的产品、规格、型号,书中的插页收载了85张产品图片。这本书的发行,对于了解中国制药设备及企业情况,为发展中国家的制药工业作出贡献,为人类医药健康事业的发展,都是很有意义的。

国家医药管理局局长 郑筱萸

The progress of the pharmaceutical industry in China is indeed very impressive during the past ten years. We note with interest to see the active role played by the renowned international companies in the country through joint research and production. China certainly possesses the capability to provide technical and institutional support to other developing countries in the pharmaceutical sector.

We welcome the publication of this compendium, which will be a useful tool for the diffusion of information and promotion of mutual cooperation.

YoMaruno

Managing Director

Investment and Technology Promotion Division

在过去十年,中国制药工业的发展确实令世人瞩目, 我们很有兴趣地看到,著名的国际大公司通过合作研究 和合作生产在这个国家所起的积极角色。在制药领域,中 国有向其它发展中国家提供技术和援助的能力。

我们欢迎这本手册的出版,它在信息传播和增进相互合作方面将成为一个有用的工具。

联合国工业发展署 投资技术促进司司长 尤. 玛儒诺

CONTENTS

MACHINERY AND EQUIPMENT FOR BULK PHARMACEUTICALS

原料药设备及机械

| DOUBLE—CIRCULATION AIR | LIFT FERMENT REACTOR |
|-------------------------------------|----------------------|
| (10) | |
| 双环流气升式发酵反应器 | 10 2 |
| K-TYPE GLASS LINED REAC | TOR (K 2000) 3 |
| 搪玻璃开式反应罐 | K2000 4 |
| GLASS LINED REACTOR (F 30 | 000) 5 |
| 搪玻璃闭式反应 罐 | F3000 6 |
| HORIZONTAL CLOSED PEELE | R CENTRIFUGE |
| (WFG 1050) ······ | •••••• 7 |
| 卧式封闭刮刀离心机 | WFG1050 8 |
| THREE-COLUMN MANUAL | TOP DISCHARGE |
| CENTRIFUGE (SS800) ······· | 9 |
| 三足式上部卸料离心机 | SS800 10 |
| TUBULAR SEPARATOR (GQ 1 | 05 A) 11 |
| 管式分离机 | GQ105A 12 |
| STARCH CENTRIFUGE (LLC 3 | 350) 13 |
| 淀粉离心机 | LLC350 14 |
| DRY DISC SEPARATOR (DRY | 400) 15 |
| 碟片分离机 | DRY400 16 |
| THREE-PHASE HORIZONTAI | SCREW DISCHARGE |
| CENTRIFUGE (LWS 320×1280) |) 17 |
| 三相卧式螺旋卸料离心机 | LWS320×1280 18 |
| QUICK-OPEN FILTER PRESS | (KY-500A) 19 |
| 快开式压滤机 | KY-500A 20 |
| PLATE-FRAME STERILE FIL | TER (LY-2, 5B) 21 |

| 无菌过滤器 | LY-2.5B | 22 |
|--------------------------|--------------------------|----|
| PLATE-FRAME PRESSURE F | TILTER (BAY 20/635-25) | |
| | •••••••••••• | 23 |
| 板框压滤机 | BAY20/635~25 ······ | 24 |
| PRESSURIZED REVERSE FLO | W EXTRACTING | |
| MACHINE SET (DT. V1 m³) | | 25 |
| 加压逆流提取机组 | DT • V1(m ³) | 26 |
| MINIATURE VERSATILE EXT | RACTING AND | |
| CONCENTRATING MACHINES | S SET [DT. V 100(L)] | 27 |
| 小型立式提取・浓缩机组 | DT0 • V100(L) | 28 |
| PRODUCTION LINE FOR MAK | ING TRADITIONAL | |
| CHINESE MEDICINE BY DYNA | AMIC WATER | |
| LEACHING (50) | | 29 |
| 中药动态水提取生产线 | 50 | 30 |
| MINIATUR MULTI-PURPOSE | E EXTRACING CON- | |
| CENTRATING & RECOVERING | G DEVICE (DT 300) ······ | 31 |
| 小型多能提取浓缩回收装置 | DT300 | 32 |
| FAST SETTING ALCOHOL RE | COVERY TOWER | |
| (JH-400) | | 33 |
| 快装酒精回收塔 | JH-400 | 34 |
| TRI EFFECT ENERGY SAVING | G EVAPORATING & | |
| CONCENTRATING DEVICE (S | JZN-1000) | 35 |
| 三效节能蒸发浓缩器 | SJZN—1000 | 36 |
| VACUUM REDUCED PRESSUF | RE CONCENTRA— | |
| TING TANK (500L) | | 37 |
| 真空减压浓缩罐 | 500L | 38 |
| SCRAPER FILM EVAPORATO | $R (GZ-2m^2) \cdots$ | 39 |
| 刮板薄膜蒸发器 | $GZ-2m^2$ | 40 |
| MULTI-FUNCTIONAL EXTR. | ACTING AND VA- | |

| CUUM CONCENTRATE-RECI | LAIMING SET (TN 50) ···· | •• |
|----------------------------|--------------------------|----|
| ************************** | | 41 |
| 多能提取-真空浓缩机组 | TN50 | 42 |
| HEAT CYCLE OVEN (RXH-14 | | |
| 热风循环烘箱 | RXH—14—B | 44 |
| HIGH-SPEED CENTRIFUGAL | | |
| CHINESE HERBAL MEDICINE | (ZLPG-40) | 45 |
| 中药用离心喷雾干燥机 | ZLPG40 | 46 |
| FLUID-BED CYCLONE CATC | HING DRYER (FGX-120) | |
| •••••••• | | 47 |
| 流化床旋风捕集式干燥机 | FGX-120 | 48 |
| VACUUM FREEZE DRYER (TI | L-LYOEPIC 40) | 49 |
| 真空冷冻干燥机 | TL-LYOEPIC40 ······· | 50 |
| FREEZE VACUUM DRYER (GI | | |
| 冷冻真空干燥机 | GDA-3 | 52 |
| VIBRATING SCREENING MAC | HINE (ZS-350) | 53 |
| 振荡筛 | ZS-350 | 54 |
| PULSE VACUUM STERILIZER | (XGI. P) | 55 |
| 脉动真空灭菌器 | XGI. P | 56 |
| ELECTRONIC STERILIZER FO | R TRADITIONAL | |
| CHINESE MEDICINE (GM20-1 | | |
| 食品中药灭菌机 | GM20—1 ······ | 58 |
| PREPARATION MACHINERY | | |
| 制剂机械 | | |
| MIXER (CH200-2) | •••••• | 59 |
| 槽式混合机 | CH200—2 ······ | 60 |
| MULTI-DIRECTION MOVEMI | | |
| •••••• | | |
| 多向运动混合机 | HD-400 | 62 |

| OSCILLATING GRANULATOR | R (YK 160) | 63 |
|---|-------------------------|-----|
| 摇摆式颗粒机 | YK160 | 64 |
| EXTRUSION GRANULATOR (| (JZL 60) | 65 |
| 挤压造粒机 | JZL60 | 66 |
| HIGH SPEED MIXER/GRANUI | CATOR (SHK-220A) ······ | .67 |
| 湿法混合颗粒机 | SHK-220A | 68 |
| FLUID-BED DRYING GRANU | JLATOR (FL-120A) ······ | 69 |
| 沸腾干燥制粒机 | FL-120A | 70 |
| FLUID-BED GRANULATOR | (FL 60 B) | 71 |
| 沸腾制粒机 | FL60B | 72 |
| FLUID—BED ROTARY GRANI | ULATING/COATING | |
| MACHINE (FXLY-60) ······· | | |
| 流化床旋转制粒/包衣机 | FXLY-60 | 74 |
| GRINDING AND GRANULATI | | |
| ••••• | | 75 |
| 粉碎整粒机 | FZ—150 ··············· | 76 |
| ROTARY TABLET PRESS (ZP | W 21A) | 77 |
| 旋转式压片机 | ZPW21A | 78 |
| ROTARY TABLET PRESS (ZP | 33) | 79 |
| 旋转式压片机 | ZP33 | 80 |
| HIGH SPEED ROTARY TABLE | ET PRESS (GZPK 37A) ··· | 81 |
| 高速旋转式 压片机 | GZPK37A | 82 |
| FULLY AUTOMATIC HIGH-S | SPEED TABLET PRESS | |
| (GZP-28) | | 83 |
| 全自动高速压片机 | GZP-28 | 84 |
| STEPLESS VARYING SPEED S | STAINLESS SUGAR- | |
| COATING MACHINE (TNW 80 | | |
| 无级调速不锈钢糖衣机 | TNW800 | 86 |
| HIGH EFFICIENCY COATING | MACHINE (BGB-150) ···· | •• |
| ,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,, | | 87 |

| 高效包衣机 | BGB150 | •••••• | 88 |
|--------------------------|---|---|-----|
| HIGH-EFFICIENCY SUGAR C | OATING M. | ACHINE | |
| (BGB-75B) | •••••• | • | 89 |
| 高效包衣机 | BGB-75B | | 90 |
| TABLET DEDUSTING UNIT (| CFQ-I) | | 91 |
| 除粉器 | | • | |
| AMPOULE CLEANING - DRYI | NG AND FIL | LING — | |
| SEALING MACHINE (BXSZ1/2 | 0) | • | 93 |
| 安瓿洗、烘、灌、封联动机组 | BXSZ1/20 | | 94 |
| AMPOULE WASHING, DRYIN | G & SEALIN | IG PRO- | |
| DUCTION LINE (ACSD 1-20) | •••••• | • | 95 |
| 安瓿洗、烘、灌、封联动机 | ACSD1-20 | •••••• | 96 |
| COMPUTERIZED AUGER FILL | | | |
| (DLF 5PC) | *** *** *** *** | | 97 |
| 微机自控螺杆分装机 | DLF-5PC | •••• | 98 |
| COMPUTER-CONTROLLED S | CREW TYP | E POW- | |
| DER BOTTLING MACHINE (F | ZL2-120) ·· | | 99 |
| 微机数控螺杆粉剂灌装机 | FZL2-120 | ••••• | 100 |
| ANTIBIOTIC POWDER PRODU | CTION LIN | E | |
| (HFZ-300) | | ••••• | 101 |
| 抗生素粉针生产联动线 | HFZ-300 | | 102 |
| RUBBER STOPPER WASHING | STERILIZA | TION— | |
| DRYING MACHINE (JS-140) | ************* | | 103 |
| 胶塞洗涤一灭菌一干燥联合机组 | JS-140 ···· | | 104 |
| PROGRAM CONTROL RUBBEI | R STOPPER | WASHER | |
| (MJK-80) | • | | 105 |
| 程控胶塞清洗机 | MJK80 | | 106 |
| MULTI FUNCTIONAL FULLY | -AUTOMA | TIC UL- | |
| TRASONIC RUBBER-CORK W | ASHER (DO | Q X-1B) | 107 |
| 多功能全自动超生波胶塞清洗机 | DQX—1B | | 108 |

| WATER BATH TYPE STERILIZER (SR 2540/90) ······· | 109 |
|---|-----|
| 水浴式灭菌柜 SR2540/90 ·············· | 110 |
| QUICK COOLING STERILIZER (KLM-B) | 111 |
| 快速冷却灭菌器 KLM-B ······· | 112 |
| BOTTLE INFUSION PRODUCTION LINE (PSY 100/500E | 3) |
| | 113 |
| 玻瓶输液生产联动线 PSY100/500B | 114 |
| SEMI-AUTOMATIC CAPSULE FILLING MACHINE | |
| (JTJ-II) | 115 |
| 半自动胶囊充填机 | 116 |
| SEMI-AUTOMATIC CAPSULE FILLING MACHINE | |
| (NCB-200B) | 117 |
| 半自动胶囊充填机 NCB200B | 118 |
| SEMI AUTOMATIC CAPSULE FILLING MACHINE | |
| (JTJ-III) | 119 |
| 半自动胶囊充填机 JTJ—III ·································· | 120 |
| HARD CAPSULE FILLING MACHINE (YJF 15) | 121 |
| 硬胶囊分装机 YJF 15 ··································· | 122 |
| FULLY AUTOMATIC CAPSULE FILLING MACHINE | |
| (ZTJ-400) | 123 |
| 全自动胶囊充填机 ZTJ-400 | 124 |
| AUTOMATIC HARD CAPSULE FILLING MACHINE | |
| (CFM 600) | 125 |
| 全自动胶囊充填机 | 126 |
| CAPSULE POLISHING MACHINE (JMJ-1) | 127 |
| 胶囊磨光机 | 128 |
| SOFT GELATIN ENCAPSULATOR (RJNJ -2) | 129 |
| 软胶囊机 | 130 |
| AUTOMATIC SEAMLESS ENCAPSULATOR (DWJ-1) | ••• |
| •••••••••••••••••••••••• | 131 |

| 全自动无缝滴丸机 | OWJ-I | 132 |
|----------------------------|------------------------|-----|
| GRANULATING AND COATING | MACHINE | |
| (BZJ-1000F) | | 133 |
| 包衣造粒机 | BZJ—1000F | 134 |
| AUTOMATIC PILL MAKING MA | ACHINE (ZW-120) ··· | 135 |
| 中药自动制丸机 2 | ZW—120 ······ | 136 |
| THREE ROLLER HONEY BOLUS | S MACHINE | |
| (ZTM80-3) | | 137 |
| 三辊蜜丸机 2 | ZTM80—3 ······ | 138 |
| OINTMENT FILLING MACHINE | C (RHG1-20B) | 139 |
| 金属软管灌封机 | RHG1—20B | 140 |
| AUTOMATIC ULTRASONIC TU | BE FILLER (ZGC 40) ··· | ••• |
| •••••• | | 141 |
| 全自动超生波软膏灌装机 2 | | 142 |
| SUPPOSITORY DRUG FILLING | | |
| 拴剂灌装机组 5 | SJ—II | 144 |
| AUTOMATIC LIQUID-FILLING | | |
| 液体灌装自动线 | YZ25/500 ······ | 146 |
| FULLY-AUTOMATIC PRODUC | CTION LINE FOR | |
| FILLING ORAL LIQUID SEALIN | | |
| ON CAPS (YLG-3) | | 147 |
| 易拉盖口服液全自动灌装生产线 | YLG—3 | 148 |
| LYING ULTRASONIC BOTTLE O | CLEANING | |
| MACHINE (CXP 5-10) | | 149 |
| 卧式超生波洗瓶机 (| CXP5—10 | 150 |
| PHARMACEUTICAL MILI | LING MACHINERY | |
| 药用粉碎 | 机械 | |
| UNIVERSAL COMMINUTOR (30 | OB II) | 151 |
| 万能粉碎机 3 | BOBII | 152 |

| HIGH EFFICIENT MII | LL (GFSJ-16A) 153 | |
|---|----------------------------------|--|
| 高效粉碎机 | GFSJ—16A 154 | |
| ULTRAFINE POWDER | R PULVERIZER (CWF-600) 155 | |
| 超微粒粉碎机组 | CWF—600 ····· 156 | |
| SELECTING PARTICI | LE MILL ASSEMBLY (TF-700) | |
| *** *** *** *** *** *** *** *** *** | 157 | |
| 分粒式粉碎机组 | TF-700 158 | |
| ULTRA-CRYOGENIO | C MILL (DFJ-2) 159 | |
| 超低温粉碎机组 | DFJ—2 ····· 160 | |
| | | |
| PHARMACEUTIC | CAL PUIVERIZING MACHINERY | |
| | 饮片机械 | |
| SPRAY AND RINSE M | MACHINE WITH HIGHPRES- | |
| SURE WATER (BGX | 500A) ······ 161 | |
| 高压水喷射式洗净机 | BGX-500A 162 | |
| | CING MACHINE (WQY 240-2) ····· | |
| *** | 163 | |
| 往复式切药机 | WQY240-2 164 | |
| ROTATING TYPE SL | ICING MACHINE (QY 120-4) ··· 165 | |
| 旋转式切药机 | QY120—4 ····· 166 | |
| | | |
| | ATMENT EQUIPMENT FOR | |
| PHARMACEUTICAL USE | | |
| 制药用水设备 | | |
| MULTI - EFFECT T | TUBULAR WATER DISTILLER (LD | |
| - | 167 | |
| 列管式多效蒸馏水机 | LD2000 168 | |
| | ION DISTILLER (RY 1000) 169 | |
| 热压式蒸馏水机 | RY1000 170 | |

| ELECTRODIALYSIS-MIXED I | BED PURE WATER | |
|-----------------------------|---|-----|
| MAKING EQUIPMENT (10T). | | 171 |
| 电渗析一混床纯水制备生产线 | 10T | 172 |
| | | |
| PHARMACEUTICAL PAC | KAGING MACHINERY | |
| 药品包装 | | |
| ELECTROMAGNETIC INDUCT | TION SEALING CAP | |
| MACHINE (DSPF-I) | | 173 |
| 电磁感应式铝箔封口机 | DSPF-1 | 174 |
| CAPSULE PRINTER, TABLET | | 175 |
| 胶囊. 药片印字机 | AP-III | 176 |
| DRUM LABEL-STICKING MA | | 177 |
| 转鼓贴标机 | ZT20/1000 | 178 |
| FULLY-AUTOMATIC LABEL | ER FOR ANTIBIOTIC | |
| BOTTLES (TQ-3) | | 179 |
| 抗生素瓶全自动标签粘贴机 | TQ-3 | 180 |
| GRAIN AUTOMATIC PACKAG | SING MACHINE (SB-80) | ••• |
| •••••• | | 181 |
| 颗粒自动包装机 | SB—80 | 182 |
| BLISTER PACKING MACHINE | | 183 |
| 泡罩包装机 | DPA250 | 184 |
| AUTOMATIC ALUMINUM-P | LASTICS PACKER | |
| (DPT-140) | | |
| 自动铝塑包装机 | DPT140 | 186 |
| PLATE AL-PLASTIC & AL-A | AL PACKAGING | |
| MACHINE (DPP 170-S) ······· | • | 187 |
| 平板铝塑、铝一铝包装机 | DPP170—S | 188 |
| MULTI-FUNCTION AND PLA | TE TYPE AUTO | |
| ALPLASTIC BUBBLE CAP PAC | KING MACHINE | |
| (DPP-200K) | | 180 |

| 多功能辊板式自动铝塑包装机] | DPP-200K | ••••• | 190 |
|-----------------------------|----------------|---|---------|
| AUTOMATIC HIGH - SPEED D | | | |
| PACKER (ZGL-160)······ | •••••••• | • | 191 |
| 自动高速双铝箔包装机 2 | ZGL-160 ··· | • | 192 |
| PLASTIC BOTTLE FILLING, CO | OUNTING & | SEALING | |
| ASSEMBLY LINE (SBX-2)······ | ••••• | • | 193 |
| 塑料瓶包装生产线 | SBX-2 ······ | • | 194 |
| PLASTIC BOTTLE ASSEMBLY I | LINE (BXSP | 15/100) | • • • • |
| | •••••• | | 195 |
| 塑料瓶瓶装自动线 | BXSP15/100 | ••••• | 196 |
| AUTOMATIC STRAPPER (SK- | 1A) | • | 197 |
| 自动捆扎机 | SK-1A | | 198 |
| GLASS TUBE DRAWING PRODU | JCTION LIN | E (BLGX 10 |)) |
| | | | 199 |
| 玻璃拉管生产线 | BLGX10 ······ | • | 200 |
| VERTICAL AMPOULE PRODUC | | | |
| | | | |
| 立式安瓿生产线 | | | |
| GLASS TUBE VIAL MAKING LI | | | |
| 管制瓶生产线 (| | | |
| PRODUCTION LINE FOR HARD | | | |
| | | | |
| 硬胶囊生产自动线 | JSZ—III ······ | • | 206 |
| MEDICINE DETECTIN | G INSTRUM | ENT | |
| 药物检测 | 设备 | | |
| TABLET HARDNESS TESTER (| YD-II) | ••••• | 207 |
| 片剂硬度计 | YD—II ······ | | 208 |
| INTELLECTUAL PHARMACY D | DISSOLUTIO | N TESTER | |
| RC-610(100) | | ••••• | 209 |

| 智能型药物溶出度测定仪 | RC-610(100) | 210 |
|---|----------------------|---------|
| INTELLECTUAL DISINTEGRA | TION TESTER (BJ I) | • • • • |
| *************************************** | •••••• | 211 |
| 智能崩解仪 | BJI | 212 |
| INTELLIGENT DISINTEGRAT | | ••• |
| •••••• | | 213 |
| 智能崩解试验仪 | ZBS-6D | 214 |
| GEL STRENGTH TESTER (JS- | | |
| 冻力测试仪 | JS—I | 216 |
| INTELLIGENT DRUG DISSOLU | JTION TESTER (D-800) |) |
| *************************************** | | 217 |
| 智能药物溶出仪 | D-800 | 218 |
| INTELLIGENT PARTICLE DET | | 219 |
| 智能微粒检测仪 | WJ-2 | 220 |
| BACTERIAL ENDOTOXINS TE | STER (BET-16) | 221 |
| 细菌内毒素测定仪 | BET-16 | 222 |
| APPENDIX:Manufacturers' Dire | ctory ······ | 223 |
| 附录:通讯录 | | |

COMMODITY: DOUBLE-CIRCULATION AIR LIFT FER-MENT REACTOR

MODEL: 10

ŧ

APPLICATION:

It is suited for biological ferment engineering in pharmaceutical, chemical, foodstuff and light industries.

FEATURES:

It is based on the principle of air lift circulation, with which liquid circulating flow and turbulent motion are formed to get air, liquid and solid mixing evenly in tank, thus providing nutrition for microbe. This reactor can be enlarged in size, and installed outdoors.

TECHNICAL DATA:

Whole volume (m³):

10

Parameter comparison for airlift tank and mechanical stirring tank

| Acid rate or ferment unit(%): | $+5.7\sim12.1$ |
|-------------------------------|----------------|
| Capacity or ferment index(%): | 0~11 |
| Energy consumption rate(%). | $-37 \sim -56$ |

POWER(kW):

OVERALL DIMENSIONS (mm):

OTHER SPECIFICATIONS:

20, 30, 50, 100, 300, 500

产品名称:双环流气升式发酵反应器

型号规格:10

主要用途:用于医药、化工、轻工、食品等行业生物发酵工程。

工作原理及结构特点:应用气体提升原理,在气升罐内形成液体循环流动及湍动,达到气、液、固相均匀混合供给微生物营养,促进微生物健康成长,发育。气升罐易于大型化,并可安装在室外。

主要技术参数:

全容积(m³):

10

气升罐与机械搅拌罐参数对比:

产酸率或发酵单位(%):

 $+5.7\sim12.1$

生产能力或发酵指数(%):

 $0 \sim 11$

能耗比(%):

 $-37 \sim -56$

配套电机(kW):

外形尺寸(mm):

其它规格:20,30,50,100,300,500

COMMODITY: K-TYPE GLASS LINED REACTOR

MODEL: K 2000

APPLICATION:

It is used as anti-corrosion reactor in pharmaceutical, chemical and light industries.

FEATURES:

This reactor is lined by a special silicate glass through high temperature. It has both the strength of metal and anti-corrosion of glass. It can resist corrosion from various organic and inorganic acid, weak base liquid and organic solvents except for hydrofluoric acid.

TECHNICAL DATA:

| OVERALL DIMENSIONS(mm): | | |
|--|--------------------|-------|
| POWER(kW): | | 4.0 |
| Nominal diameter of stirring axle(mm): | | 80 |
| | in jac k et | 0.6 |
| Nominal pressure (MPa): | in container | 0.25 |
| Heating area of jacket(m2): | | 7.2 |
| Nominal volume(L): | | 2,000 |

OTHER SPECIFICATIONS:

K50, K100, K200, K300, K500, K1000, K1500 K3000, K5000 产品名称:搪玻璃开式反应罐

型号规格:K2000

- **主要用途:**用于医药、化工、轻工等行业的耐腐蚀反应容器。
- 工作原理及结构特点:是将特种硅酸盐玻璃经高温搪烧 于金属胎体表面上的压力容器。它具有金属的强度, 又具有玻璃的耐腐蚀性能。对于除氢氟酸以外的几 乎所有无机酸、有机酸、弱碱液、有机溶剂等,均有极 强的抗腐蚀作用。

主要技术参数:

| 公称容积(L): | 2000 |
|--------------|-------|
| 夹套换热面积: | 7.2 |
| 公称压力:容器内 | 0. 25 |
| 夹套内 | 0.6 |
| 搅拌轴公称直径(mm): | 80 |

配套电机(kW):

4.0

外形尺寸(mm):

其它规格:K50,K100,K200,K300,K500,K1000,K1500, K3000,K5000

COMMODITY: GLASS LINED REACTOR

MODEL: F 3000

APPLICATION:

It is used as anti-corrosion reactor in pharmaceutical, chemical and light industries.

FEATURES:

This reactor is lined by a special silicate glass through high temperature. It has both the strength of metal and anti-corrosion of glass. It can resist corrosion from various organic and inorganic acid, weak base liquid and organic solvent except for hydrofluoric acid.

TECHNICAL DATA:

| Nominal volume(L): | | 3,000 |
|---|--------------|-------|
| Heating area of jacket(m2): | | 9.74 |
| Nominal pressure (MPa): | in container | 0.25 |
| | in jacket | 0.6 |
| Nominal diameter of stirring axle (mm): | | 95 |
| | | |
| POWER(kW): | | 5.5 |
| OVERALL DIMENSIONS(mm). | | |

OTHER SPECIFICATIONS:

F 4000, F5000, F8000, F10000

产品名称:搪玻璃闭式反应罐

型号规格:F3000

主要用途:用于医药、化工、轻工等行业耐腐蚀反应容器。

工作原理及结构特点:是将特种硅酸盐玻璃经高温搪烧 于金属胎体表面上的压力容器,它具有金属的强度 和玻璃的耐腐蚀性能。对于除氢氟酸以外的几乎所 有无机酸、有机酸、弱碱液、有机溶剂等,均有极强的 抗腐蚀作用。

主要技术参数:

| 公称容积(L): | 3000 |
|---------------|------|
| 夹套换热面积(m²): | 9.74 |
| 公称压力(MPa):容器内 | 0.25 |
| 夹套内 | 0.6 |
| 搅拌轴公称直径(mm): | 95 |

配套电机(kW):

外形尺寸(mm):

其它规格:F4000,F5000,F8000,F10000

7

COMMODITY: HORIZONTAL CLOSED PEELER CENTRI-FUGE

MODEL: WFG 1050

APPLICATION:

It is used for separating liquid and solid in pharmaceutical, chemical, foodstuff and light industries, especially for high cleanness material or the material that can produce pernicious gas.

FEATURES:

The machine is whole closed type peeler discharge centrifuge. PLC controlled system is adopted to convert centrifugal frequency and speed. The whole process can be set by technological requirement. For the harmful or easily oxidized material, nitrogen can be filled in to protect the material.

TECHNICAL DATA:

| Diameter of bowl(mm): | 1,050 |
|-----------------------------------|-----------------------------------|
| Working volume(L): | 210 |
| Allowable loading(kg): | 230 |
| Speed of bowl(r/min): | 0~1,400 |
| Separation factor: | 0~1,145 |
| Filtering area (m ²): | 1.65 |
| POWER(kW); | 22 |
| OVERALL DIMENSIONS(mm): | $3,350 \times 1,700 \times 1,600$ |

OTHER SPECIFICATIONS:

产品名称:卧式封闭刮刀离心机

型号规格:WFG 1050

- **主要用途**:用于医药、化工、轻工、食品等行业的液、固分离。尤其适用分离洁净度要求高或有毒害气体逸出的物料。
- 工作原理及结构特点:该机采用全密闭型,刮刀卸料。PLC 自动控制,变频调速。进料、分离、洗涤、甩干、出料过程、转鼓转速、时间可根据工艺要求设定。对于有毒害易氧化物料可在机内充氮保护。并配有回收母液装置,便于有机溶剂回收。

主要技术参数:

| 转鼓直径(mm): | 1050 |
|--------------|--------|
| 工作容积(l): | 210 |
| 许用装料量(kg): | 230 |
| 转鼓转速(r/min): | 0~1400 |
| 分离因数: | 0~1145 |
| 过滤面积(m²): | 1.65 |
| 配套电机(kW): | 22 |

 $3350 \times 1700 \times 1600$

其它规格:

外形尺寸(mm):

COMMODITY: THREE-COLUMN MANUAL TOP DIS-CHARGE CENTRIFUGE

MODEL: SS 800

APPLICATION:

It is suitable for chemical, foodstuff, pharmaceutical, mining and metallurgy and environmental protection industries, etc.

FEATURES:

The product flows into the centrifugal force field in the bowl from the upper feeding pipe and is forced to be deliquated and separated on the bowl wall. The liquid phase is centrifuged into the housing space through solid phase product and gauze from perforate bowl wall and is dislodged out of the liquid outlet from machine bottom pan. Solid phase product would be left behind in the bowl. After machine stopping, it would be discharged manually from the top.

TECHNICAL DATA:

| Internal diameter of bowl(mm): | 800 |
|----------------------------------|-------------------|
| Filtering area(m ²): | 1 |
| Working volume of bowl(L): | 100 |
| Speed of bowl (r/min): | 1,200 |
| Max. loading volume(kg): | 135 |
| Max. separation factor of bowl: | 640 |
| POWER(kW): | 5. 5 |
| OVERALL DIMENSIONS(mm): | 1,670×1,270×1,000 |

OTHER SPECIFICATIONS:

SS450, SS600, SS600-NC, SS800-NC, SS1000

产品名称:三足式上部卸料离心机

型号规格:SS 800

主要用途:用于悬浮液的液固分离及纤维状物料脱水。

工作原理及结构特点:物料在离心力作用下,液相经滤布和转鼓壁上的开孔,被甩离转鼓,由排液口排出,滤饼留在转鼓内,分离完成后停机,人工上部卸出。

主要技术参数:

| 转鼓内径(mm): | 800 |
|------------|--------------------------------|
| 过滤面积(m²): | 1 |
| 工作容积(1): | 100 |
| 转速(r/min): | 1200 |
| 最大装料量(kg): | 135 |
| 分离因数: | 640 |
| 配套电机(kW): | 5.5 |
| 外形尺寸(mm): | $1670 \times 1270 \times 1000$ |
| | |

其它规格:

CHINA PHARMACEUTICAL EQUIPMENT 11

COMMODITY: TUBULAR SEPARATOR

MODEL: GQ 105A

APPLICATION:

GQ type is used for separation of various suspension which is difficult to separate, especially for the separation of solid-liquid suspension which has low concentration, small granularity and very small difference of solid-liquid gravity.

FEATURES:

Fine solid phase particles with great density are settled on the inner wall of the bowl gradually to form the residue layer to be discharged manually after the machine stops. The liquid phase flows out of the liquid outlet at the bowl's upper parts after clarification.

TECHNICAL DATA:

| Internal diameter of bowl(mm): | | 105 |
|--------------------------------|---|--------|
| Speed of bowl(r/min): | | 16,000 |
| Max. separation factor: | • | 15,025 |
| Settling volume(L): | | 5.5 |
| | | |
| POWER(kW): | | 2. 2 |

 $800 \times 1,000 \times 500$

OTHER SPECIFICATIONS:

OVERALL DIMENSIONS(mm):

GQL75, GQQ105, GQL125, GQL150

产品名称:管式分离机

型号规格:GQ 105A

主要用途:用于分离各种乳浊液。尤其适于三相比重差甚微的液一液分离及含有少量杂质的液一液一固分离。

工作原理及结构特点:物料由进液口进入高速转鼓,离心力使密度不同的两相液体分层。密度大的液相形成外环,小的形成内环,微量杂质沉积到鼓壁上。分层后两相液体同时向转鼓上部流动,由各自排液口排出。鼓壁上的沉积物停机后卸出。

主要技术参数:

转鼓内径(mm):

105

转速(r/min):

16000

分离因数:

15025

沉积容积(l):

5.5

配套电机(kW):

2.2

外形尺寸(mm):

 $800 \times 1000 \times 500$

其它规格:

CHINA PHARMACEUTICAL EQUIPMENT 13

COMMODITY: STARCH CENTRIFUGE

MODEL: LLC 350

APPLICATION:

This machine is suitable for the separation of suspension, solid phase gravity of which is bigger than its liquid phase gravity, especially for separation of maize starch and protein.

FEATURES:

Raw starch product flows into the screw from the charge-in inlet, and then into between the bowls by passing through the holes of screw wall, starch is settled on the inside wall of the bowl by the centrifugal force field and it is pushed to the top into the case from the residue opening of outside wall of the screw, to be diluted by the dilute water and then it flows out the starch outlet, settled raw starch liquid overflows through the overflow opening at the bottom and then drains out from the liquid outlet.

TECHNICAL DATA:

| Big end internal diameter of bowl(mm): | 350 |
|--|-----------------|
| Speed of bowl(r/min): | 3,000 |
| Separation factor of bowl: | 1,800 |
| POWER(kW): | 15 |
| OVERALL DIMENSIONS(mm); | 1,310×830×1,495 |

OTHER SPECIFICATIONS:

产品名称:淀粉离心机

型号规格:LLC 350

主要用途:用于固相比重大于液相比重的悬浮液分离,更适合用于玉米淀粉与蛋白质分离。

工作原理及结构特点:淀粉原浆由进料口流入螺旋器内腔料仓,在离心力场的作用下,淀粉沉向转鼓内壁, 靠螺旋推料器将淀粉推向转鼓小端,由排渣口排入外壳,再经稀释水稀释,由排出口流出。沉降后的原浆液由转鼓底上的溢流口流出机外。

主要技术参数:

沉降转鼓大端直径(mm):

350

沉降转鼓转速:

3000

分离因数:

1800

配套电机(kW):

15

外形尺寸(mm):

 $1310 \times 830 \times 1495$

其它规格:

CHINA PHARMACEUTICAL EQUIPMENT 15

COMMODITY: DRY DISC SEPARATOR

MODEL: DRY 400

APPLICATION:

It is used for separation of emulsion and suspension which contains 1-5% solid product, for example, separation of gum, nitroglycerine and extraction of tetracycline, bacterium, penicillin, and biological alkali.

FEATURES:

The product liquid flows into the bottom of bowl's wall through the bottom edge of disc seals and their holes. Solid product is settled into the inner surface wall of the bowl and forms into the residue by the centrifugal force field and then the residue will be discharged manually after the machine stops. Light liquid flows along up the outer tapered surface of discs to the main shaft center, passed through the light liquid centripetal pump and drains out from the light liquid outlet, Heavy liquid flows along the inner tapered surface of disc down to the bowl's wall and then passed through the heavy liquid centripetal pump and drains out from the heavy liquid outlet.

TECHNICAL DATA:

| Internal diameter of bowl(mm): | 420 |
|--------------------------------|---------------------------------|
| Speed of bowl(r/min): | 6,525 |
| Separation factor: | 10,000 |
| Production capacity (m³/h): | 5 |
| POWER(kW): | 11 |
| OVERALL DIMENSIONS(mm): | $1,500 \times 900 \times 1,500$ |

OTHER SPECIFICATIONS:

DRY 500

产品名称:碟片分离机

型号规格:DRY 400

主要用途:用于乳浊液及固形物含量为 1~5%的悬浮液的分离。

工作原理及结构特点:料液由上部中心进料管流至转鼓底部,经碟片座下面的分流孔趋向转鼓壁,在离心力场的作用下,比液体重的固相物沉向转鼓内壁,待停机后人工卸除。轻液向轴心流动至上部,经轻液向心泵排出,重液经重液向心泵排出。

主要技术参数:

转鼓内径(mm):

420

转鼓转速(r/min):

6525

分离因数:

10000

生产能力(m³/h):

5

配套电机(kW):

11

外形尺寸(mm):

 $1500 \times 900 \times 1500$

其它规格:

COMMODITY: THREE-PHASE HORIZONTAL SCREW DISCHARGE CENTRIFUGE

MODEL: LWS 320×1280

APPLICATION:

It is used for treating any concentration suspension in which solid particles are fine and creamy i.e. density fluctuation of suspension can not influence the separating result. It can also deal with the separation of the three-phase mixture (liquid-liquid-solid) of solid density being greater than liquid one and the separation of liquid-solid separation.

FEATURES:

The product liquid is brought into the bowl through a centripetal pump, solid product is settled into outer circular layer by the centrifugal force. Due to the different gravity of liquid phase in the bowl, two kinds of light and heavy liquid flows are formed and divided into two parts by the baffle at the big end of the bowl, light liquid flows out through the centripetal pump at the light liquid outlet and heavy liquid out at heavy liquid outlet.

TECHNICAL DATA:

| Diameter and length of bowl(mm): | $320 \times 1,280$ |
|-------------------------------------|---------------------------------|
| Speed of bowl(r/min): | 4,000 |
| Differential speed of screw(r/min): | 30 |
| Separation factor: | 2,870 |
| POWER(kW): | 15/18.5 |
| OVERALL DIMENSIONS(mm): | $2,400 \times 1,300 \times 800$ |

OTHER SPECIFICATIONS:

产品名称:三相卧式螺旋卸料离心机

型号规格:LWS 320×1280

主要用途:用于医药化工等行业液一液一固分离。

工作原理及结构特点:料浆经向心泵导管加入转鼓,在离心力作用下,固相物(滤渣)沉降到最外环,由于螺旋与转鼓的转差,沉渣被推至转鼓小端干燥区,由排渣口排出,转鼓内液相的浮浊液由于比重不同,形成轻、重两液环,被转鼓大端挡板将轻、重液分开,轻液经向心泵由轻液口流出,重液由重液口排出。

主要技术参数:

转鼓直径和长度(mm):320×1280转鼓转速(r/min):4000螺旋差动转速(r/min):30分离因数:2870

配套电机(kW):

15(18.5)

外形尺寸(mm):

 $2400 \times 1300 \times 800$

其它规格:

COMMODITY: QUICK-OPEN FILTER PRESS

MODEL: KY-500A

APPLICATION:

It is applicable to the pharmaceutical, chemical, foodstuff and light industries for separating solid from liquid in the suspensions.

FEATURES:

The pressurized suspension passes through the filter plate covered with the filter medium to separate the solid from the liquid. The equipment adopts the quick-open construction, horizontal filter plate and specially-designed self-sealing device.

TECHNICAL DATA:

| Filter area (m ²): | 1 |
|--------------------------------|------------------------------|
| Filter layer(pcs): | 7 |
| Working pressure(MPa): | ≤ 0. 4 |
| Working temperature (°C): | ≤ 200 |
| Cylinder dimension(mm): | \emptyset 510 \times 360 |

POWER(kW):

OVERALL DIMENSIONS(mm);

 \emptyset 518×814

OTHER SPECIFICATIONS:

GKY-600

产品名称:快开式压滤机

型号规格:KY-500A

- 主要用途:用于制药、化工、轻工、食品等行业中悬浮液的液固分离。
- 工作原理及结构特点:带压悬浮液通过敷有过滤介质的 滤板达到液固分离的目的。设备采用快开门结构,水 平式滤板、特殊设计的自密封。

主要技术参数:

过滤面积(m²):

1

过滤层数(层):

7

工作压力(MPa):

≤0.4

工作温度(℃):

≤200

简体规格(mm):

 \emptyset 510 \times 360

配套电机(kW):

外形尺寸(mm):

 \emptyset 518×814

其它规格:GKY-600

COMMODITY: PLATE-FRAME STERILE FILTER

MODEL: LY-2.5B

APPLICATION:

It is mainly used in sterile filtration, clarification, and purification processes in pharmaceutical, chemical, food and light industries.

FEATURES:

When pressurized fluid passes through sterile filtering media sandwiched between the plates and the frames, which are alternate arranged, the particles and impurities are removed, so that clear sterile filtration is obtained.

TECHNICAL DATA:

| Filtering area(m ²): | 2. 5 |
|----------------------------------|--------|
| Working pressure(MPa): | ≤0.42 |
| Filtering temperature (°C): | <100 |
| pH of filtered material: | 2.5~10 |

POWER(kW): 1.1 OVERALL DIMENSIONS(mm): $1,161 \times 511 \times 860$

OTHER SPECIFICATIONS:

产品名称:无菌过滤器

型号规格:LY-2.5B

主要用途:用于医药、化工、食品、轻工等行业对液体过滤,达到除菌、提纯、澄清的目的。

工作原理及结构特点:液体物料加压后进入交错排列的 滤板、滤框间,经夹在二者之间的除菌板时,微粒及 杂质被截留,从而使液体澄清、除菌。

主要技术参数:

过滤面积(m²):

2.5

工作压力(MPa):

≪0. 42

过滤温度(℃);

<100

过滤物料 pH 值:

 $2.5 \sim 10$

配套电机(kW):

1.1

外形尺寸(mm):

 $1161 \times 511 \times 860$

其它规格:LY-0.8

COMMODITY: PLATE-FRAME PRESSURE FILTER

MODEL: BAY 20/635-25

APPLICATION:

It is suitable for filtering suspension in pharmaceutical, chemical, light and metallurgical industries.

FEATURES:

Product liquid is pumped into compression chamber by the pump. Liquid phase is drained out from bottom through residue layer and filtration fabric under the pressure and solid is left behind on the filter frame and would be manually removed at regular intervals.

TECHNICAL DATA:

| Filter area (m ²): | 20 |
|--------------------------------|-----------------------------------|
| Quantity of filter plate(pcs): | 25 |
| Quantity of filter frame(pcs): | 26 |
| Total volume in frame(L): | 252 |
| Operating pressure (MPa): | 0.5 |
| POWER(kW): | 0.6 |
| OVERALL DIMENSIONS(mm): | $3,940 \times 1,220 \times 1,270$ |

OTHER SPECIFICATIONS:

BMY 84/1000-45A, BAJ60/80, BMJ80/1000, BAS350-2M²

产品名称:板框压滤器

型号规格:BAY 20/635~25

主要用途:用于制药、化工、轻工、冶金等行业,过滤悬浮液。

工作原理及结构特点:料液由料泵输入压缩室,在压力作用下,液相通过渣层及滤布,由下部排出,固相则留在滤框内,人工定期清除。

主要技术参数:

外形尺寸(mm):

计滤面积(m2).

| TT NO HI WY (III): | 20 |
|---------------------|-----|
| 滤板数(片): | 25 |
| 滤框数(片): | 26 |
| 框内总容积(1): | 252 |
| 工作压力(MPa): | 0.5 |
| | |
| 配套电机(kW): | 0.6 |

20

 $3940 \times 1220 \times 1270$

其它规格:BMY84/1000-45A,BAJ60/80,BMJ80/1000 BAS350-2m²

COMMODITY: PRESSURIZED REVERSE FLOW EXTRACTING MACHINE SET

MODEL: DT. V1 (m³)

APPLICATION:

It is used for extraction, evaporation and sterilization of herbal medicine.

FEATURES:

With its special new type sealing structure, this equipment is the third generation of tanks for extracting traditional Chinese medicine. The tank body is of vertical structure that makes drug discharging be extremely convenient and can be used independently or in combination with other tanks. It is even more effective for large quantity medicine extraction, thus increasing efficiency, saving energy, yielding high concentration liquid medicine and shortening the time for both extraction and concentration.

TECHNICAL DATA:

| Designed pressure (MPa): | Container | 0.15 |
|----------------------------|--------------------------|-------|
| | Jac k et | 0.3 |
| Designed temperature (°C): | Container | 127 |
| | Jacket | 143 |
| Volume(m³): | | 1 |
| POWER(kW): | | |
| OVERALL DIMENSIONS(mm): | $3,800\times3,100\times$ | 5,200 |

OTHER SPECIFICATIONS:

DT. V 1. 5, 2, 3, 5, 6 (m^3)

产品名称:加压逆流提取机组

型号规格:DT・V1(m³)

主要用途:用于中草药提取、蒸制、灭菌。

工作原理及结构特点:中草药装入罐内后,通入蒸汽加热、提取。配附件后,可完成回流提取、强制循环提取、提油、酒精提取和回收残渣中的酒精。

主要技术参数:

设计压力(MPa):

罐内 0.15、夹套 0.3

设计温度(℃):

罐内 127、夹套 143

容 积(m³):

1

配套电机(kW):

0.6

外形尺寸(mm):3800×3100×5200~4740×3580×8500

其它规格:DT·V1.5,2,3,5,6(m³)

COMMODITY: MINIATURE VERSATILE EXTRACTING AND CONCENTRATING MACHINE SET

MODEL: DT. V 100(L)

APPLICATION:

It is used for extraction, evaporation and sterilization of herbal medicine.

FEATURES:

This equipment is mainly used for middle and small-sized pharmaceutical factories, hospitals and scientific research institutes, developing new medicines or manufacturing small quantity but high-grade medicines.

TECHNICAL DATA:

Volume(L):

100

Working pressure (MPa):

0.15

Jacket pressure (MPa):

0.15

POWER(kW):

OVERALL DIMENSIONS(mm): $1.200 \times 1.600 \times 1.800$

OTHER SPECIFICATIONS:

DT. V (L) 200, 300

产品名称:小型立式提取、浓缩机组

型号规格:DT·V100(L)

主要用途:用于中草药提取、蒸制、灭菌。

工作原理及结构特点:本机主要适用于中、小型药厂、医院、科研机构、开发新药或小批量生产高等级的药品。

主要技术参数:

容 量(L):

100

工作压力(MPa):

0.15

线性夹套压力(MPa):

0.15

配套电机(kW):

外形尺寸(mm):

 $1200 \times 1600 \times 1800$

其它规格:

COMMODITY: PRODUCTION LINE FOR MAKING TRADITIONAL CHINESE MEDICINE BY DYNAMIC WATER LEACHING

MODEL: 50

APPLICATION:

This production line is used for dynamic water leaching of traditional Chinese medicine.

FEATURES:

This production line consists of dynamic water leaching unit, three-stage liquid-solid separation unit, low temperature vacuum concentration unit and spray drying unit. The whole line is operated under sealed conditions, to improve both yield and quality of the medicines.

TECHNICAL DATA:

| Raw material feed(kg/batch): | 50 |
|---------------------------------|-------------------|
| Max. feed water (L/batch): | 600 |
| Temperature for leaching (°C): | 95±5 |
| Max. concentration speed(kg/h): | 300 |
| Mode for drying: | centrifugal spray |
| Liquid feed (kg/h): | 20 |

POWER(kW):

OVERALL DIMENSIONS(mm):

OTHER SPECIFICATIONS:

100, 200, 300

产品名称:中药动态水提取生产线

型号规格:50

主要用途:用于中药材的水提取。

工作原理及结构特点:本生产线由动态水提取、三级液固 分离、低温真空浓缩、喷雾干燥组成。生产过程为全 密封状态,提高药品的收率和质量。

主要技术参数:

原料处理量(kg/批):

50

最大添加水量(L/批):

600

提取温度(℃):

 95 ± 5

最大浓缩速度(kg/h):

300

干燥方式:

离心喷雾

料液处理量(kg/h):

20

配套电机(kW):

外形尺寸(mm):

其它规格:100,200,300

COMMODITY: MINIATURE MULTI-PURPOSE EXTRACTING CONCENTRATING & RECOVERING DEVICE

MODEL: DT 300

APPLICATION:

For the extraction and concentration of effective substance from herbal medicine and juice, widely used in pharmaceutical, chemical and foodstuff industries.

FEATURES:

The complete system of extraction and concentration serves the functions as: warmly submerging, filtering, extracting fragrance oil, evaporating and recovering solvent, etc., and it can also keep normal pressure and can boost and reduce pressure, it can also operate in reverse flow. Quick evaporation and low-temperature concentration are the two major features of the system.

TECHNICAL DATA:

| Concentrator capacity(L): | 300 |
|---|-----------------------------------|
| Concentrator vacuum(MPa): | 0.05~0.08 |
| Extracting & concentrating temp(°C): | 45~ 100 |
| Heating power for extracting(kW): | 6 |
| Heating power for concentrating (kW): | 6 |
| POWER(kW): | 12 |
| OVERALL DIMENSIONS(mm): | $3,300 \times 2,400 \times 2,500$ |

OTHER SPECIFICATIONS:

DT 50, DT 100, DT 200, DT 500

产品名称:小型多能提取浓缩回收装置

型号规格:DT 300

主要用途:用于医药、化工、食品等行业中草药、果汁等有效物质的提取和浓缩。

工作原理及结构特点:提取、浓缩组成一个系统,具有温浸、芳香油提取、蒸制及溶媒回收等功能。可常压、减压、加压和回流煎煮。蒸发速度快,浓缩温度低。

主要技术参数:

提取罐容积(l): 300
 浓缩真空度(MPa): 0.05~0.08
 提取、浓缩温度(℃): 45~100
 提取电加热功率(kW): 6
 浓缩电加热功率(kW): 6

配套电机(kW):

12

外形尺寸(mm):

 $3300 \times 2400 \times 2500$

其它规格:DT50,DT200,DT300,DT500

COMMODITY: FAST SETTING ALCOHOL RECOVERY TOWER

MODEL: JH--400

APPLICATION:

This device is used in recovering the thin alcohol in pharmaceutical, light and chemical industries, for distillation of methanol and the other mediums.

FEATURES:

Diluted alcohol in distilling caldron is evaporated by heating through coil pipe, the materials will go up through distilling tower to the top and then into condenser. Qualified thick alcohol can be recovered after condensing.

TECHNICAL DATA:

| Recovering capacity(kg/h): | 120~200 |
|----------------------------------|---------------|
| Recovering density(%): | 90~95 |
| Disposing capacity(kg/h): | 200~300 |
| Volume of distilling caldron(L): | 1500 |
| Proportion of alcohol loss(%): | ≤ 2. 3 |
| Height of tower(m): | 8 |

POWER(kW):

OVERALL DIMENSIONS(mm): $3,560 \times 1,400 \times 9,500$

OTHER SPECIFICATIONS:

产品名称:快装酒精回收塔

型号规格:JH-400

主要用途:用于制药、化工、轻工等行业稀酒精回收,也用于甲醇等溶媒的蒸馏。

工作原理及结构特点:蒸馏釜内的稀酒精液,通过盘管加 热汽化,经蒸馏塔内填料上升至塔顶并导入冷凝器, 凝结成合格浓酒精后回收。

主要技术参数:

| 回收能力(kg/h): | 120~200 |
|-------------|---------|
| 回收浓度(%): | 90~95 |
| 处理量(kg/h): | 200~300 |
| 蒸馏釜容积(l): | 1500 |
| 酒精损耗率(%): | €2.3 |
| 塔高(m): | 8 |

配套电机(kW):

外形尺寸(mm): 3560×1400×9500

其它规格:

COMMODITY: TRI EFFECT ENERGY SAVING EVAPO-RATING & CONCENTRATING DEVICE

MODEL: SJZN-1000

APPLICATION:

It is widely used for material concentration in pharmacy, foodstuff, dairy industry, chemical industry and light industry.

FEATURES:

- 1. Steam can be utilized many times, over-heat self evaporating of condensed water and efficient heat preservation of evaporator, greatly reduce the hot steam consumption. It consumes less than 0. 4kg of hot steam as 1kg of water is evaporated.
- 2. It lowers the temperature of material evaporation and meets the requirment of thermal sensitive material concentrating by using vacuum evaporating operation.
- 3. The final material gravity can reach 1.35 \sim 1.4 by using its first degree to make further concentration and raise the ratio of concentration.

800~1,000

 $6.000 \times 1.500 \times 3.600$

TECHNICAL DATA:

Evaporating capacity (kg/h):

| e(MPa): 0. 05~0. 09(on the p | oiezometer) | |
|-------------------------------|--|--|
| vaporation temperature: | | |
| 16.7KPa (125mmHg) | 95° | |
| 61.3KPa (460mmHg) | 70° | |
| 85.3KPa (640mmHg) | 55° | |
| Steam consumption(kg/h): <400 | | |
| POWER (kW): | | |
| | vaporation temperature; 16.7KPa (125mmHg) 61.3KPa (460mmHg) 85.3KPa (640mmHg) | |

OTHER SPECIFICATIONS: SIZN-500

OVERALL DIMENSIONS (mm):

产品名称:三效节能蒸发浓缩器

型号规格:SJZN-1000

主要用途:广泛用于制药、食品、化工、轻工业的物料浓缩。

工作原理及结构特点:一次蒸汽可多次利用,冷凝水的过热自蒸发及蒸发器的有效保温,而使加热蒸汽消耗量大为降低,蒸发 1kg 的水消耗加热蒸汽不到0.4kg。降低了物料蒸发温度,满足热敏性物料的浓缩要求。可利用一效作进一步浓缩,提高浓缩比,使最终物料比重达到1.35~1.40。

主要技术参数:

蒸发能力(kg/h):

800~1000

加热蒸汽压力(MPa):

0.05~0,09(壳压)

真空度与蒸发温度:

I 效 16.7 KPa(125mmHg) 95℃

I效 61.3 KPa(460mmHg) 70℃

■效 85.3 KPa(640mmHg) 55°C

配套电机(kW):

外形尺寸(mm):

 $6000 \times 1500 \times 3600$

其它规格:SJZN--500

COMMODITY: VACUUM REDUCED PRESSURE CON-CENTRATING TANK

MODEL: 500L

APPLICATION:

With no pollution to the productive environment, the concentrator is characterized by easy operation and reliable use and suitable for concentrating and distilling the liquid materials, such as traditional Chinese, foodstuff, etc.

FEATURES:

This equipment consists mainly of the concentrator, condensing cooler, vapor-liquid separator, material accepting tank, etc. The concentrator features the structure of sleeve. In order to reduce the boiling point of materials, the concentrator is equipped with a vacuum-pumping port so that customers can adjust the vacuum pureness according to materials' sensitivity to heat by themselves.

TECHNICAL DATA:

| Capacity(L): | 500 |
|---------------------------|-------------|
| Evaporation capacity(kg): | 100 |
| Steam pressure (MPa): | ≤0.15 |
| Vacuum(KPa): | ≤ 60 |

POWER(kW):

OVERALL DIMENSIONS(mm): $2,100\times1,200\times3,400$

OTHER SPECIFICATIONS:

产品名称:真空减压浓缩罐

型号规格:500L

- 主要用途:本设备对生产环境无任何污染,用于浓缩及蒸馏液态物料,如中药、食品等,具有易于操作,使用安全可靠的特点。
- 工作原理及结构特点:本设备主要由浓缩罐、冷凝器、汽一液分离器、物料接收罐等组成。浓缩罐带有套筒,为了降低物料的沸点,浓缩罐配有真空泵装置,这样可按不同物料对热敏度的要求,自行调节罐内真空度。

主要技术参数:

容 量(L):

500

蒸发量(kg):

100

蒸汽压力(MPa):

≤0.15

真空度(kPa):

≤60

配套电机(kW):

外形尺寸(mm):

 $2100 \times 1200 \times 3400$

其它规格:

COMMODITY: SCRAPER FILM EVAPORATOR

MODEL: GZ-2m²

APPLICATION:

Widely used in pharmaceutical, petroleum, chemical and foodstuff industries, etc., especially for evaporating and concentrating the thermo-sensative, high-viscosity, easy foaming and easily filthy materials.

FEATURES:

The liquid material flows from the top to the centrifugal distributor and is evenly scattered on the internal wall of the cylinder with jacket. Affected by centrifugation of 240(300 r/min, the scraper slipper moves clinging to the cylinder wall to scrap the material into thin film, push it down to evaporate quickly so as to prevent charring.

TECHNICAL DATA:

| Evaporating area(m ²): | 1 | 2 |
|------------------------------------|---|--------|
| Vacuum (MPa): | | 0.086 |
| Max. discharge viscosity (cPa. s): | | 50,000 |
| Jacket pressure(MPa): | | 0.3 |

| POWER(kW); | 3 |
|--------------------------|-------------------------------|
| OVERALL DIMENSIONS (mm): | $950 \times 830 \times 3,200$ |

OTHER SPECIFICATIONS:

GZ-0. 8m²

产品名称:刮板薄膜蒸发器

型号规格:GZ-2m2

主要用途:适用于医药、石油、化工、食品等行业,高粘度、 热敏性、易起泡、易结垢物料的蒸发浓缩。

工作原理及结构特点:料液由顶部切向进入离心分配器,被均匀分散在带夹套的简体内壁上,刮板滑块在 300 转/分的离心作用下紧贴筒壁运动,将物料刮成薄膜,向下推送,迅速蒸发,避免结焦。

主要技术参数:

蒸发面积(m²):

2

直空度(MPa):

0.086

出料最大粘度(cPa·s):

50000

夹套压力(MPa):

0.3

配套电机(kW):

3

外形尺寸(mm):

 $950 \times 830 \times 3200$

其它规格:GZ-0.8m²

COMMODITY: MULTI-FU

MULTI-FUNCTIONAL EXTRACTING AND VACUUM CONCENTRATE-RE-

CLAIMING SET

MODEL: TN 50

APPLICATION:

This set is suitable for hospitals, pharmaceutical plants, foodstuff and chemical industries for reclaiming organic solvent as well as for extracting and concentrating Chinese herbal medicine, fruit juice, etc.

FEATURES:

The set has extracting and concentrating-two independent operating systems. It possessed the functions of soaking, percolating, steaming as well as extracting aromatic oil and reclaiming organic solvent. It can keep normal pressure, and can boost & reduce pressure, it can also operate in reverse flow. Concentration proceeded on the condition of low temperature and high evaporating speed.

TECHNICAL DATA:

| Extracting capacity(L): | 50 |
|--|---------------------------------|
| Working pressure in pot(MPa): | 0.15 |
| Working pressure in jacket(MPa): | 0.15 |
| Concentrating capacity(L): | 50 |
| Vacuum(MPa): | ≥0.03 |
| Working pressure in heating pipe(MPa): | 0. 15 |
| POWER(kW); | 1.85 |
| OVERALL DIMENSIONS(mm): | $970 \times 1,110 \times 1,900$ |

OTHER SPECIFICATIONS:

TN 100, TN 200

产品名称:多能提取•真空浓缩机组

型号规格:TN 50

主要用途:适用于医药、化工、食品等行业中草药、果汁等有效成份的提取和浓缩。

工作原理及结构特点:提取、浓缩组成一个系统,具有温浸、渗漉、芳香油提取、蒸制及溶媒回收等功能。可常压、减压、加压和回流煎煮,蒸发速度快,浓缩温度低。

主要技术参数:

| 提取罐容量(l): | 50 |
|-------------|--------------|
| 工作压力(MPa): | 0.15 |
| 夹套压力(MPa): | 0.15 |
| 浓缩罐容量(l): | 50 |
| 真空度(MPa): | ≥0.03 |
| 加热管压力(MPa): | 0.15 |
| 配套电机(kW): | 1.85 |
| 外形尺寸(mm): | 970×110×1900 |

其它规格:TN100、200

COMMODITY: HEAT CYCLE OVEN

MODEL: RXH-14-B

APPLICATION:

It is used in pharmaceutical, foodstuff and chemical industries, for drying wet materials.

FEATURES:

Using axial fan and heat exchanger, the heat exchange is carried out between paralleled layers of air flow and raw material, so that the moisture can be removed and the raw materials can be dried.

TECHNICAL DATA:

| Drying capacity(kg/batch): | 100 |
|---------------------------------|-------|
| Vapor consumption (kg/h): | 20 |
| Wind volume(m ³ /h): | 1,400 |
| Temperature difference (°C): | 2 |
| | |
| POWER(kW): | 1.1 |

 $2,430\times1,200\times2,433$

OTHER SPECIFICATIONS:

RXH-27-B, RXH-41-B, RXH-54-B RXH-5-C, RXH-25-A

OVERALL DIMENSIONS(mm):

产品名称:热风循环烘箱

型号规格:RXH-14-B

主要用途:用于制药、食品、化工等行业产品的干燥。

工作原理及结构特点:利用轴流风机、换热器等使空气呈平行层流与物料进行热交换,除去水份,使物料干燥。

主要技术参数:

 干燥能力(kg/次):
 100

 蒸汽耗量(kg/h):
 20

 风 量(m³/h):
 1400

 上下温差(℃):
 2

配套电机(kW):

1.1

外形尺寸(mm):

 $2430 \times 1200 \times 2433$

其它规格: RXH-27-B, RXH-41-B, RXH-54-B, RXH-5-C, RXH-25-A

COMMODITY: HIGH-SPEED CENTRIFUGAL SPRAY

DRYER FOR CHINESE HERBAL

MEDICINE

MODEL: ZLPG-40

APPLICATION:

This spray dryer is used in leaching concentration liquid of Chinese herbal medicine.

FEATURES:

The dry air is cleaned through primary-efficiency, mid-efficiency and high-efficiency filters, enters the atomizer and sprays out the material. The material will be dried in the drying room and be collected in the clean room after going through the spiral segregator. The temperature of incoming and outgoing air can be adjusted automatically. This equipment can handle the materials which contains sugar and which is of low soft-point.

TECHNICAL DATA:

| Chamber diameter(mm): | 4,000 |
|-------------------------------------|--------|
| Power of atomizer(kW): | 4 |
| Rotation speed of atomizer (r/min): | 15,800 |
| Diameter of disk(mm): | 150 |

POWER(kW):

29

OVERALL DIMENSIONS(mm):

OTHER SPECIFICATIONS:

ZLPG-45, ZLPG-50

产品名称:中药用高速离心喷雾干燥机

型号规格:ZLPG-40

主要用途:用于中药行业提取浓缩液的喷雾干燥。

工作原理及结构特点:干燥空气经初效、中效、高效过滤器净化,进入雾化器将物料雾化喷出,在干燥室内完成干燥,经旋风分离器在净化间收料。进出风温度自动控制。能较好地处理含有糖份及软化点较低的物料。

主要技术参数:

 干燥室直径(mm):
 4000

 雾化器功率(kW):
 4

 雾化器转速(r/min):
 15800

雾化盘直径(mm): 150

配套电机(kW):

29

外形尺寸(mm):

其它规格:ZLPG-45,ZLPG-50

COMMODITY: FLUID-BED CYCLONE CATCHING DRY-ER

MODEL: FGX-120

APPLICATION:

It is used in pharmaceutical, foodstuff and chemical industries for drying powder and granular materials.

FEATURES:

Based on the principle of fluid-bed dryer, cyclone separator is adopted in this machine to catch powder, where powder trapped by the first stage cyclone separator returns to the main unit through vortex gas-seal valve so as to dry simultaneously a batch of wet powder trapped in early stage of operation and to overcome the trouble with cleaning powder in the filtering bag of fluid-bed dryer during operation and shortcoming of fiber shedding.

TECHNICAL DATA:

| Production capacity(kg/time): | 80~160 |
|-------------------------------------|-----------------------------------|
| Air-in temperature (°C): | 20~120 |
| Moisture in finished product(%): | 0.5~3 |
| Obtainable ratio(%): | 99 |
| Pressure of gas source(MPa): | 0.6 |
| Compressed air consumption(m³/min): | 0.6 |
| POWER(kW): | 11 |
| OVERALL DIMENSIONS(mm): | $1,800 \times 1,800 \times 2,800$ |

OTHER SPECIFICATIONS:

FGX3, FGX10, FGX30, FGX60, FGX300, FGX1000

产品名称:流化床旋风捕集式干燥机

型号规格:FGX-120

主要用途:用于医药、食品、化工等行业粉状、颗粒状物料的干燥。

工作原理及结构特点:本机在流化床干燥机原理基础上, 采用旋风分离器捕集粉末,第一级旋风分离器捕集 的粉末,通过涡旋气封阀回到主机内,这样避免了作 业初始阶段捕集到的湿粉末无法一个批次同时干燥 的问题,克服了流化床干燥机过滤袋在作业中清理 粉末的麻烦和脱落纤维的缺点。

主要技术参数:

| 生产能力(kg/次): | 80~160 |
|-----------------|--------------------------------|
| 进风温度(℃): | 20~120 |
| 成品水份(%): | 0.5~3 |
| 收得率(%): | 99 |
| 气源压力(MPa): | 0.6 |
| 压缩空气耗量(m³/min); | 0.6 |
| 配套电机(kW): | 11 |
| 外形尺寸(mm): | $1800 \times 1800 \times 2800$ |

其它规格: FGX3, FGX10, FGX30, FGX60, FGX300, FGX1000

COMMODITY: VACUUM FREEZE DRYER

MODEL: TL-LYOEPIC 40

APPLICATIONS:

It is used to produce lyophilizing products, for instance, blood products, vaccine, antibiotics, biological products, biochemical medicines, gene engineering products, etc.

FEATURES:

This machine consists of a drying chamber containing shelves which may be either heated or cooled, a refrigerated vapor condenser with the associated refrigeration group vacuum pumping system, hydraulic rubber stopper pressing system, micro-computer auto-control system, and CIP, SIP. The design to eliminate surface welds (which obtains National Patent) together with high quality components selection improve performance and reliability.

TECHNICAL DATA:

| Useable shelf area (m²): | 40.5 |
|----------------------------------|-----------------------------------|
| The number of shelves: | 18±1 |
| Lowest temperature of shelf(°C): | -50 |
| Condenser capacity (kg): | 810 |
| POWER (kW): | 180 |
| OVERALL DIMENSIONS (mm): | $1,200 \times 2,600 \times 4,000$ |

OTHER SPECIFICATIONS:

产品名称:真空冷冻干燥机

型号规格:TL-LYOEPIC 40

主要用途:用于医药行业血液制品、疫苗、抗生素、生物制品、基因工程等产品的真空冷冻干燥。

工作原理及结构特点:本机由冻干箱、冷阱、蒸汽冷凝器、液压压塞、微机自控系统及 CIP、SIP 等组成。精良的设计消除了表面大焊缝现象(获得国家专利),加之选用高品质元件精心制作,大大改进了产品性能和可靠性。

主要技术参数:

板层面积(m²):

40.5

板层数目:

18 + 1

板层最低温度(℃):

-50

冷凝器凝冰能力(kg):

810

配套电机(kW):

180

外形尺寸(mm):

 $1200 \times 2600 \times 4000$

其它规格:

COMMODITY: FREEZE VACUUM DRYER

MODEL: GDA-3

APPLICATIONS:

It is used to produce lyophilizing product, for instance, biological product medicine and biochemical medicine, foods etc.

FEATURES:

The liquid products are frozen in the low temperature, then the frozen products are lyophilized, i. e, by heating the ice of products, is sublimed to vapor in the vacuum.

TECHNICAL DATA:

| Refrigeration speed(h): | (Normal to -50°) <1.5 |
|-------------------------------------|-------------------------------------|
| Min shelf temperature(℃): | -55 |
| Condenser refrigeration speed (min) | : (normal to -40° C) <30 |
| Min condenser temperature (°C): | -70 |
| Vacuumed speed(min): | (to 10 Pa) < 30 |
| Highest attained vacuum(Pa); | 0.01 |
| Vacuum leakness (Pa/l sec): | <5 |
| Max condensed water(kg): | 60 |
| • | |
| POWER (kW): | 22 |
| OVERALL DIMENSIONS (mm): | $2,800 \times 2,100 \times 2,500$ |

OTHER SPECIFICATIONS:

产品名称:冷冻真空干燥机

型号规格:GDA--3

主要用途:用于医药、生化行业生产冻干制品。

工作原理及结构特点:将物品中的水份在低温下冻结,真空,加热,使固态的水升华,达到干燥物品的目的。

主要技术参数:

| 搁板面积(m²): | 3. 24 |
|-------------------|--------------------------------|
| 搁板尺寸(mm): | 915×915 |
| 搁板数量(个): | 5 |
| 搁板间距(mm): | 100 |
| 搁板温度(℃): | $-45 \sim 70$ |
| 捕水器温度(℃): | -60 |
| 冷凝器每批最大捕水量(kg/批): | 50 |
| 极限真空(Pa): | <2.7 |
| 配套电机(kW): | 25 |
| 外形尺寸(mm): | $3000 \times 1600 \times 2800$ |
| | |

其它规格:

COMMODITY: VIBRATING SCREENING MACHINE

MODEL: ZS-350

APPLICATION:

It is used for screening the different size granule.

FEATURES:

It consists of hopper, vibrating chamber, connecting shaft, and electric motor. When the force of adjustable eccentric weight driven by the motor acts on the centre of main shaft, it may be out of balance and causes eccentricity forces, so the regular eddy of material formed in the vibrating chamber is changed.

TECHNICAL DATA:

| Production capacity(kg/h): | 60~500 |
|------------------------------------|--------|
| Screen meshes(mesh): | 12~200 |
| Main shaft rotating speed (r/min): | 1,380 |

POWER(kW): 0.55 **OVERALL DIMENSIONS(mm):** 540×540×1,060

OTHER SPECIFICATIONS:

产品名称:振荡筛

型号规格: ZS-350

主要用途:用于筛分颗粒大小比例不等的物料。

工作原理及结构特点:本筛由料斗、振荡室、联轴器、电机等组成,可调节的偏心重锤经电机驱动传递到主轴中心线,在不平衡状态下产生离心力,使物料改变在筛内形成的轨道旋涡。

主要技术参数:

生产能力(kg/h):

 $60 \sim 500$

过筛目数(目):

12~200

主轴转速(r/min):

1380

配套电机(kW):

0.55

外形尺寸(mm):

 $540 \times 540 \times 1060$

其它规格: ZS-515, ZS-650

COMMODITY: PULSE VACUUM STERILIZER

MODEL: XGI. P

APPLICATION:

It is used in medical and biological organizations to sterilize apparatus, utensils, fabric and liquids.

FEATURES:

The equipment applies vacuum pump eliminating chamber gas, filling steam to sterilize, and drawing off steam again. It can keep the sterilized subject dry, then filling the filter sterilized gas and eliminating vacuum. PLC controlled system is adopted.

TECHNICAL DATA:

| Working pressure(MPa): | 0.21 |
|------------------------------|---------------------------|
| Working temperature (°C): | 132 |
| Steam source pressure (MPa): | 0.29~0.59 |
| Water source pressure(MPa): | 0.15~0.3 |
| Chamber volume(mm): | $1,700\times610\times910$ |

POWER(kW):

5.5

OVERALL DIMENSIONS(mm):

 $2,600\times1,200\times1,900$ (double doors)

OTHER SPECIFICATIONS:

XGI. YS

产品名称:脉动真空灭菌器

型号规格:XGI.P

- **主要用途:**用于医院、生物制品单位对器械、器皿、织物或液体进行消毒灭菌。
- 工作原理及结构特点:本设备用真空泵脉动抽空气方式 排除灭菌室内空气,充入蒸汽灭菌,再用真空泵抽去 蒸汽,可保持被灭菌物品的干燥,然后充入过滤无菌 空气,消除真空。PLC 自动控制,预置程序。

主要技术参数:

工作压力(MPa):

0.21

工作温度(℃):

132

汽源压力(MPa):

 $0.29 \sim 0.59$

水源压力(MPa):

 $0.15 \sim 0.3$

灭菌室容积(mm):

 $1700 \times 1610 \times 910$

配套电机(kW):

5.5

外形尺寸(mm):

 $2600 \times 1200 \times 1900$

其它规格:XG1.YS

COMMODITY: ELECTRONIC STERILIZER FOR TRADI-TIONAL CHINESE MEDICINE

MODEL: GM20-1

APPLICATION:

It is used to sterilize medicinal herbs, semi-finished products or finished medicine, also for natural health foods and beverages.

FEATURES:

The machine makes use of the electromagnetic wave radiation osmosis and the polarization of medium for sterilization. It is sure that the effective component is not damaged, the sterilizing speed is fast, and it is operated easily and economically. All of the machines are in line with the GMP standard.

TECHNICAL DATA:

| Treated capacity(kg/h): | 100 |
|-------------------------------|-------|
| Treating duration(min/batch): | 15~25 |
| Output power(kW): | 13~20 |
| Power consumption(kVA): | 13~26 |

POWER(kW):

2.0

OVERALL DIMENSIONS (mm):

1,800(1,500(600 (main machine) 3,300(800(1,500(subsidiary machine)

OTHER SPECIFICATIONS:

GM10-1, GM5-1

产品名称:食品中药灭菌机

型号规格:GM 20-1

主要用途:用于中药材、半成品及成药的干燥、灭菌。

工作原理及结构特点:利用高频电子辐射浸透和介质极 化原理,对物料进行灭菌处理,具有灭菌快、操作方 便、费用低等特点。

主要技术参数:

处理量(kg/批):

100

处理时间(min/批):

15~25

输出功率(kW):

13~20

电源(kVA):

13~26

配套电机(kW):

2.0

外形尺寸(mm):

1800×1500×600(主机)

3300×800×1500(辅机)

其它规格:GM10-1,GM5-1

COMMODITY: MIXER

MODEL: CH200-2

APPLICATION:

It is used for evenly mixing wet and dry material in pharmaceutical, chemical and foodstuff industries.

FEATURES:

By means of mechanical driving, material is revolving back and forth to get mixing evenly, and it can also be discharged automatically.

TECHNICAL DATA:

| Working volume(L): | 200 |
|-----------------------------------|---------------------------------|
| Revolving angle(°): | 108 |
| Total volume(L): | 250 |
| Diameter of agitating blade (mm): | 525 |
| Agitating speed(r/min): | 20 |
| POWER(kW): | 5. 5 |
| OVERALL DIMENSIONS(mm): | $1,752 \times 660 \times 1,190$ |

OTHER SPECIFICATIONS:

CH50, CH100, CH150, WH150, WH200

产品名称:槽式混合机

型号规格:CH 200-2

主要用途:用于制药、化工、食品等行业,均匀混合粉状干湿物料。

工作原理及结构特点:通过机械传动,搅拌桨旋转,推动物料往复翻动,均匀混合,并能自动倒料。

主要技术参数:

| 工作容积(l): | 200 |
|--------------|------|
| 翻转角度: | 108° |
| 总容积(l): | 250 |
| 搅拌桨叶直径(mm): | 525 |
| 搅拌转速(r/min): | 20 |

配套电机(kW):

5.5

外形尺寸(mm):

 $1752 \times 660 \times 1190$

其它规格:CH50、CH100、CH150、WH150、WH200

COMMODITY: MULTI-DIRECTION MOVEMENT MIXER

MODEL: HD-400

APPLICATION:

It is used in pharmaceutical, chemical, light, foodstuff and agricultural processing industries, etc.

FEATURES:

The barrel makes self rotation, while the shaft rotating and at the same time swinging up and down, with the swinging speed equal to 4 times of the spinning speed, so as to make the particles material be mixed under the centrifugal action-free condition, and to decrease the gravity segregation.

TECHNICAL DATA:

| POWER(kW): | 4 |
|---------------------------|-----|
| Shaft revolution (r/min): | 20 |
| Max. loading weight(kg): | 200 |
| Max. loading volume(L): | 300 |
| Mixing barrel volume(L): | 400 |

OVERALL DIMENSIONS(mm): $1,500\times2,200\times1,500$

OTHER SPECIFICATIONS:

HD-5, HD-20, HD-100, HD-200, HD-600, HD-800, HD-1000

产品名称:多向运动混合机

型号规格:HD-400

主要用途:用于制药、化工、食品、轻工等行业混合粉状或颗粒状物料。

工作原理及结构特点:混合桶自转时也同时以 4 倍于自转速度进行上下摆动,物料在无离心力作用下混合,减少了比重的偏析。该机具有三维运动,增加了混合桶内的混合点,混合均匀度高,装载系数大。

主要技术参数:

| 混合桶容积(1): | 400 |
|--------------|-----|
| 装料容积(l): | 300 |
| 装料重量(kg): | 200 |
| 主轴转速(r/min): | 20 |

配套电机(kW):

4

外形尺寸(mm):

 $1500 \times 2200 \times 1500$

其它规格: HD—5, HD—20, HD—100, HD—200, HD—600, HD—800, HD—1000

COMMODITY: OSCILLATING GRANULATOR

MODEL: YK 160

APPLICATION:

It is used for granulating wet powdered material in pharmaceutical, chemical and foodstuff industries, etc.

FEATURES:

By means of reciprocal of the powder scraper, the materials are extruded to form granules through the screen.

TECHNICAL DATA:

Production capacity(kg/h):

300~400

Diameter of powder scraping roller(mm):

160

Oscillating frequency of powder scraping roller (times/min):

65

POWER(kW):

1.5

OVERALL DIMENSIONS (mm):

 $1,030 \times 450 \times 1,100$

OTHER SPECIFICATIONS:

YK 140-A, YK 160-A, YK 160-B

产品名称:摇摆式颗粒机

型号规格: YK 160

主要用途:用于医药、化工、食品等行业,将潮湿的粉状物料制成颗粒。

工作原理及结构特点:通过机械传动,使刮粉轴往复摆动,将物料从筛网中挤出,制成颗粒。

主要技术参数:

生产能力(kg/h):

300~400

刮粉轴直径(mm):

160

刮粉轴摆动次数(次/min):

65

配套电机(kW):

1.5

外形尺寸(mm):

 $1030 \times 450 \times 1100$

其它规格:YK140-A,YK160-A,YK160-B

COMMODITY: EXTRUSION GRANULATOR

MODEL: JZL 60

APPLICATION:

It is applicable to the pharmaceutical, chemical, foodstuff and forage industries for producing granules of different sizes wetly.

FEATURES:

Wet materials are directly extruded into granules through the screw and net plate. The granular products possess the advantages of the regular shape, tight construction, with little fine powder, anti-moisture and long storage period.

TECHNICAL DATA:

Production capacity (kg/h):

Finish product size (mm):

Center distance of 2 screws (mm):

Extrusion shaft speed (r/min):

15~16, stepless regulating

POWER (kW):

1.5

OVERALL DIMENSIONS (mm):

890×650×1,460

OTHER SPECIFICATIONS:

产品名称:挤压造粒机

型号规格:JZL 60

主要用途:用于医药、化工、食品、饲料等行业,湿法制作 各种规格的颗粒。

工作原理及结构特点:湿物料通过螺杆及网板硬性挤压 形成颗粒。所制颗粒形状规则,质地紧密,细粉少,不 易吸湿,产品保存期长。

主要技术参数:

生产能力(kg/h):

20~50

成品粒度(mm):

 \emptyset 0.5~3

两螺杆中心距(mm):

60

挤压轴转速(r/min):

15~16 无级变速

配套电机(kW):

1.5

外形尺寸(mm):

 $890 \times 650 \times 1460$

其它规格:

COMMODITY: HIGH SPEED MIXER/GRANULATOR

MODEL: SHK-220A

APPLICATIONS:

This equipment is used for mixing powder materials and making them into granule, such as tablets, capsules and medicines which can be taken after being mixed with boiling water, and also can be used for material mixing in food and chemical industries.

FEATURES:

Material is stirred and evenly mixed by using stirring paddle. Then atomized bonding agent is added. The even wet granule is formed by the both actions of stirring paddle and cutter. This equipment adopts explosion-proof structure.

TECHNICAL DATA:

| Effective volume (L): | 220 |
|--------------------------|-----|
| Productivity (kg/batch): | 100 |

Time for making granule (min/batch): dry and mix:

about 2 min

about 1~4min

making granule:

| | - | |
|--|-------|------|
| Mesh number for making granule (mesh): | 30 | ~100 |

Stirring speed (r/min): 80/270
Cutting speed (r/min): 1,500/3,000

Cutting speed (r/min): 1,500/3,000 **POWER (kW)**: 16.5

OVERALL DIMENSIONS (mm): $2,195\times970\times1,570$

OTHER SPECIFICATION:

SHK-220B; SHK-220FB; SHK-110; SHK-20A; SHK-20TS

产品名称:湿法混合颗粒机

型号规格:SHK-220A

主要用途:用于制药等行业,将粉状物料混合与制粒。

工作原理及结构特点:原料经搅拌桨搅拌,使其混合均匀。加入雾状粘合剂,在搅拌桨和切刀同时作用下形成均匀的湿颗粒。

主要技术参数:

有效容积(l):

220

生产能力(kg/批):

100

制粒时间(min/批):

干混约 2/制粒约 1~4

制粒目数(目):

30~100

搅拌速度(r/min):

180/270

切刀速度(r/min):

1500/3000

配套电机(kW):

16.5

外形尺寸(mm):

 $2195 \times 970 \times 1570$

其它规格: SHK-220B, SHK-220FB, SHK-110, SHK-20A, SHK-20TS

COMMODITY: FLUID-BED DRYING GRANULATOR

MODEL: FL-120A

APPLICATIONS:

It is used for one-step granulating and granules coating in pharmaceutical, foodstuffs and chemical industries.

FEATURES:

The machine utilizes high speed hot air flow to make the seeds powder on the fluid-bed be suspended and fluidized, then turns it into porous granules by spraying liquid materials onto the surface of the seeds. The machine combines multiple functions of mixing, spraying, drying, granulating and coating into one. The inside surface of the machine is very smooth and clean without dead corner and it meets 'GMP' requirement.

TECHNICAL DATA:

Production capacity (kg/h):

100~400

Diameter of container (mm).

1,200

Operation temperature (C):

Automatic regulation forms indoor temperature to 160°C

Moisture content:

< 0.2

POWER (kW):

18.5

OVERALL DIMENSIONS (mm): $2.300 \times 1.530 \times 3.230$

OTHER SPECIFICATION:

FL-5A FL-15A FL-30A FL-60A FL-200A FL-300A FL-500A

产品名称:沸腾干燥制粒机

型号规格:FL-120A

主要用途:用于医药、食品、化工等行业的制粒及干燥。

工作原理及结构特点:高速热气流使流化床内种子粉末 成流态化,液态物料喷成雾状凝聚在种子表面结成 多孔状颗粒。集混合、喷雾干燥、制粒、颗粒包衣功能 于一体,设备内表面光洁、无死角。

主要技术参数:

处理能力(kg/h):

 $100 \sim 400$

容器直径(mm):

1200

操作温度(℃):

室温~160(自动调节)

产品水份(%):

≤0.2

配套电机(kW):

18.5

外形尺寸(mm):

 $2300 \times 1530 \times 3230$

其它规格:FL-5A,FL-15A,FL-30A,FL-60A,FL-200A,FL-300A,FL-500A

COMMODITY: FLUID-BED GRANULATOR

MODEL: FL 60 B

APPLICATIONS:

It is applicable to the pharmaceutical, food, chemical and light industries for material granulating and drying.

FEATURES:

The powdery materials are fluidized and mixed evenly by means of the negative pressure in the hermetic vessel. By spraying the atomized bonding agent into the vessel, the powdery materials inside will be moist and coagulate, then be dried by the hot air to form into porous granules.

TECHNICAL DATA:

| Production capacity (kg/batch): | 60 |
|---------------------------------|----------------------|
| Finish product grain (mm): | Ø1.13∼Ø0.173 |
| Operation temperature (°C): | room temperature 100 |
| POWER (kW): | 8 |

OVERALL DIMENSIONS (mm): $2,200\times1,400\times3,350$

OTHER SPECIFICATION:

FL-3 FL-60C FL-120C FL-120B

产品名称:沸腾制粒机

型号规格:FL 60B

主要用途:用于制药、食品、轻工、化工等行业对物料制粒和干燥。

工作原理及结构特点:在密闭的容器中,通过负压使粉状物料呈流化状态,并达到均匀混合,然后喷入雾化粘合剂,使粉状的物料湿润,并产生凝聚,再经热空气干燥,形成多孔状颗粒。

主要技术参数:

生产能力(kg/批):

60

成品粒度(mm):

 \emptyset 1.13 \sim \emptyset 0.173

操作温度(℃):

室温~100

配套电机(kW):

8

外形尺寸(mm):

 $2200 \times 1400 \times 3350$

其它规格:FL3,FL60C,FL120C,FL120B

COMMODITY: FLUID-BED ROTARY GRANULATING/ COATING MACHINE

MODEL: FXLY-60

APPLICATIONS:

It is used in pharmaceutical, foodstuffs and chemical industries for making compact ball granules, pellets and fine & micro granules as well as coating of pellets.

FEATURES:

By the action of rotary disk, the fluidized powder formed in fluid-bed spirally rotates along periphery of fluid-bed, where, adhesive liquid is sprayed on the powder so that it is agglutinated. Then by the centrifugal action, it rolls constantly along smooth wall surface until it is formed into compact granules.

TECHNICAL DATA:

| Capacity(L): | 10 | 5 |
|-------------------------------------|----------------------------------|---|
| Steam consumption(kg/h): | 18 | 0 |
| Compressed air consumption(m³/min): | 0. | 9 |
| Number of nozzles(piece): | | 2 |
| | | |
| POWER (kW); | 11.5 7. | 5 |
| OVERALL DIMENSIONS (mm): | $1,800 \times 2,400 \times 3,20$ | 0 |

OTHER SPECIFICATION:

FXLY-5 FXLY-30 FXLY-100 FXLY -200

产品名称:流化床旋转制粒/包衣机

型号规格:FXLY-60

主要用途:用于医药、化工、食品等行业制备致密球形颗粒、丸粒、精细微粒颗粒及丸粒的包衣。

工作原理及结构特点:流化床内形成流态化的粉末在旋转圆盘作用下,沿流化床周边以螺旋运动的方式旋转,粘合液喷洒其上,使其胶结,再由于离心力作用使其不断沿光滑壁面滚动,使颗粒达到致密球形。

主要技术参数:

| 容量(l): | 105 |
|-----------------|-----|
| 蒸汽耗量(kg/h): | 180 |
| 压缩空气耗量(m³/min): | 0.9 |
| 喷嘴数量(个): | 2 |

配套电机(kW):

11.5,7.5

外形尺寸(mm):

 $2000 \times 2400 \times 5200$

其它规格:FXLY-5 FXLY-30 FXLY-100 FXLY

COMMODITY: GRINDING AND GRANULATING MACHINE

MODEL: FZ-150

APPLICATIONS:

It is used in pharmaceutical, foodstuffs and fine chemistry processing industries, not only for grinding viscose, gummy, heat-sensitive, moist and hard, but also for brittle materials.

FEATURES:

The material drops into the cone -shaped screening chamber, owing to the shearing between the high-speed rotating blade and the sieve surface, the material is ground into a smaller size and discharged through the sieve meshes. The ground size is decided by utilizing sieves with different mesh and adjusting the gap between the blades.

TECHNICAL DATA:

| OVERALL DIMENSIONS (mm): | $650 \times 480 \times 1,200$ |
|-------------------------------------|-------------------------------|
| POWER (kW): | 0.55 |
| Regulation range of speed(r/min): | 30~3650 |
| Chamber temperature increased (°C): | ≤ 40 |
| Grinding grain fineness (mm): | 0.3~3.35 |
| Production capacity(kg/h): | 12~150 |

OTHER SPECIFICATION.

FZ-300 FZ-450 FZ-700 FZ-1000

产品名称:粉碎整粒机

型号规格:FZ-150

主要用途:用于制药、食品、精细化工等行业的物料粉碎, 尤其适合粉碎热敏性药物、废药片及废胶囊,并可用 于湿料制粒、整粒。

工作原理及结构特点:物料由进料口进入筛分室,受到回转刀与筛网之间的剪切作用,被粉碎出网。颗粒大小可通过回转刀与网面之间的距离及筛网目数决定。该机发热量小,粉尘少,粉碎粒度均匀,便于清洗。

主要技术参数:

生产能力(kg/h):

 $12 \sim 150$

粉碎粒度(mm):

 $0.3 \sim 3.35$

整粒室温升(C):

≤40

调速范围(r/min):

30~3650

配套电机(kW):

0.55

外形尺寸(mm):

 $650 \times 480 \times 1200$

其它规格:FZ-300 FZ-450 FZ-700 FZ-1000

COMMODITY: ROTARY TABLET PRESS

MODEL: ZPW21A

APPLICATIONS:

It is mainly used in pharmaceutical industry and also in the chemical, foodstuffs, electronic, pottery, plastic and powdery metallurgy industries. The machine can press not only round tablets but also different geometrical shape tablets, double-layered and annular tablets. The tablets may have the impressed letters on both sides.

FEATURES:

Through mechanical transmission, it makes upper and lower punches fixed in turret along track go up and down, and finishes feeding, filling, pressing tablet out with continuous action. Pressure is controlled by hydraulic system. ZPW21A is a new designed type to meet GMP standard. The parts which contact the material are made of stainless steel and their surface is specially processed to avoid pollution.

TECHNICAL DATA:

| Punches: | 21 |
|-------------------------------|-------------------------------|
| Production capacity: | 15,000~45,000pcs/h |
| Max. pressure (KN): | 80 |
| Max. pre-pressure (KN): | 20 |
| Max. depth of filling (mm): | 18 |
| Max thickness of tablet (mm): | 9 |
| Turret speed (r/min): | 12~36 |
| POWER (kW): | 3 |
| OVERALL DIMENSIONS (mm): | $900 \times 840 \times 1,640$ |

OTHER SPECIFICATION.

ZP100 (ZPY136, ZPY134, ZPY129, ZPY124, ZPY150 ZP19, ZP21, ZP25, ZP33, ZP33B, ZP35, ZP35A

产品名称:旋转式压片机

型号规格:ZPW 21A

- 主要用涂:用于制药、化工、食品、电子、陶瓷、塑料、粉末 冶金等行业将颗粒状物料压制成圆形片、异形片、双 层片、环形片及双面刻字片。
- 工作原理及结构特点:通过机械传动、使转盘上的上、下 冲杆沿轨道作升降运动,完成加料、充填、压片、出片 的连续动作。本机压片力由液压系统设定。ZPW21A 是按药品 GMP 要求设计的新机型,除保留原机的基 本压制功能外,在压制室上作了改进,压制室采用四 门有机玻璃封闭结构,与外界环境隔离,室内与物料 接触的各构件均采用不锈钢制成或者表面经过特殊 防腐处理,避免了物料受污染。

主要技术参数:

生产能力(万片/h):

1.5 \sim 4.5

冲模数(副):21 最大充填深度(mm): 18

最大压片力(kN):80 最大片剂厚度:

最大压片直径(mm),20 转台转速(r/min): 12~36

配套电机(kW):

3

9

外形尺寸(mm):

 $900 \times 840 \times 1640$

其它规格: ZP100 系列(ZPY136、ZPY134、ZPY129、 ZPY124, ZPY150), ZP19, ZP21, ZP25, ZP33, ZP33B、ZP35、ZP35A

COMMODITY: ROTARY TABLET PRESS

MODEL: ZP33

APPLICATIONS:

It is mainly used in pharmaceutical industry and also in the field of chemistry, foodstuffs, electronics for pressing various granular raw materials into tablets or medicines.

FEATURES:

The upper and lower punches move up and down along the tracks when the turret is in rotation, so as to achieve the goal of filling, pressing tablets. The tablets can be discharged continuously. The mechanism of the machine is double-press type.

TECHNICAL DATA:

| Production capacity(pcs/h): | 43,000~110,000 | |
|-------------------------------|-------------------------------|--|
| Max. dia. of tablet(mm): | 12 | |
| Max. depth of filling (mm): | 15 | |
| Max. thickness of tablet(mm): | 6 | |
| POWER (kW): | 2. 2 | |
| OVERALL DIMENSIONS (mm). | $930 \times 900 \times 1.600$ | |

OTHER SPECIFICATION:

ZP15 ZPT-15 ZP19 ZP19A ZP33A ZP33B ZP33 ZP35 ZP35A ZP37 YH20 产品名称:旋转式压片机

型号规格:ZP 33

主要用途:用于制药、化工、电子、食品等行业将颗粒状原料压制成片剂或块剂。

工作原理及结构特点:通过机械传动,使转盘上的上、下冲杆沿轨道作升降运动,完成充填、压片、出片的连续动作。本机为双压式结构。

主要技术参数:

| 生产能力(万片/h): | 4.3~11 |
|-------------|--------|
| 压片力(kN): | 40 |
| 最大压片直径(mm): | 12 |
| 充填深度(mm): | 15 |
| 片剂厚度(mm): | 6 |
| | |

配套电机(kW):

2.2

外形尺寸(mm):

 $930 \times 900 \times 1600$

其它规格:ZP15,ZPT—15,ZP19,ZP19A,ZP33A,ZP33B, ZP33C,ZP35,ZP35A,ZP37,YH20

COMMODITY: HIGH SPEED ROTARY TABLET PRESS

MODEL: GZPK37A

APPLICATIONS:

For pressing granular raw materials into tablets in pharmaceutical, chemical, food and electrical industry.

FEATURES:

The machine is automatically controlled by computer. The computer control system monitors the tablet production on both sides. With different menu function, it completes the collection of trial production parameters, control of production data, automation adjustment, rejection of unqualified tablets, regular printing / sampling, trouble alarming.

TECHNICAL DATA:

| Production capacity(pcs/h): | 10,000~340,000 |
|-------------------------------|-------------------------------|
| Max. pressure (kN): | 100 |
| Max. pre-pressure(kN): | 20 |
| Max. dia. of tablet(mm): | 25 |
| Max. thickness of tablet(mm): | 11 |
| | |
| POWER (kW): | 11 |
| OVERALL DIMENSIONS (mm): | $900 \times 840 \times 1,640$ |

OTHER SPECIFICATIONS:

GZPK100 (138, 132, 126)

产品名称:高速旋转式压片机

型号规格:GZPK 37A

主要用途:用于制药、化工、食品、电子等行业将颗粒状物料压制成片剂。

工作原理及结构特点:本机采用计算机控制,每项菜单分别进行试生产参数采集和生产参数数据管理,具有自动调节、剔废、定时打印、采样及故障报警等功能。

主要技术参数:

| 生产能力(万片/h): | 1~34 |
|-------------|------|
| 最大主压力(kN): | 100 |
| 最大预压力(kN): | 20 |
| 最大压片直径(mm): | 25 |
| 最大片厚(mm): | 11 |
| | |

配套电机(kW):

11

外形尺寸(mm):

 $1600\times1520\times2200$

其它规格:GZPK100 高速系列旋转式压片机(138 型、132 型、126 型)

COMMODITY: FULLY AUTOMACTIC HIGH-SPEED TABLET PRESS

MODEL: GZP-28

APPLICATIONS:

It is used for pressing tablet with granular material.

FEATURES:

The rotation disc is driven by a stepless speed regulator to make the upper and lower heads move relatively in die, and then filling, pre-pressing, pressing, producing tablet, etc are completed. Through transmitting to press the signal by the computer, it can automatically control the weight of tablet, reject the waste, and print all kinds of statistical data.

TECHNICAL DATA:

| 28 |
|----------------|
| 50,000~181,500 |
| 5.5~12 |
| 80 |
| 16 |
| |

POWER (kW):

OVERALL DIMENSIONS (mm): $690 \times 980 \times 1,826$

OTHER SPECIFICATION:

GZP-28B

产品名称:全自动高速压片机

型号规格:GZP-28

主要用途:用于将颗粒状原料压制成片剂。

工作原理及结构特点:本机由无级变速器带动转盘转动,使上、下冲头在冲模中产生相对运动,完成充填、预压、加压、出片等动作,通过对压力信号的计算机处理,实现对片重的自动控制及剔废,可打印各种统计数据。

主要技术参数:

冲头数(副):28生产能力(万片/h):5~18.15压片直径(mm):5.5~12最大主压力(kN):80

最大充填深度(mm): 16

配套电机(kW):

外形尺寸(mm):

 $690 \times 980 \times 1826$

其它规格:GZP-28B

COMMODITY: STEPLESS VARYING SPEED STAINLESS SUGAR-COATING MACHINE

MODEL: TNW 800

APPLICATIONS:

It is used in pharmaceutical plant for producing large number of sugar-coated tablets, pills or ball-shaped foodstuffs, as well it can be used for polishing sugar-coated tablets.

FEATURES:

The stepless speed regulator is driven by motor to make the sugar coating pot finish coating. The rotation speed can be adjusted, due to which various tablets can be coated.

TECHNICAL DATA:

| Diameter of pot (mm): | Ø800 |
|--------------------------------|---------------------------------|
| Production capacity (kg): | 50 |
| Rotation speed of pot (r/min): | 0~50 |
| Power for heating (kW): | 6 |
| POWER (kW): | 1.47 |
| OVERALL DIMENSIONS (mm); | $800 \times 1,300 \times 1,650$ |

OTHER SPECIFICATION:

TNW200 TNW300 TNW400 TNW500 TNW1000 TNW1000I

产品名称:无级调速不锈钢糖衣机

型号规格:TNW800

- **主要用途**:用于制药厂大批量生产糖衣片、丸,滚制球形食品和化工产品,并可用于糖衣片打光。
- 工作原理及结构特点:电机驱动无级变速器使糖衣锅旋转完成包衣过程。转速可调,能满足各种药片在不同转速下的包衣要求。

主要技术参数:

锅体直径(mm):Ø800生产能力(kg/次):50锅体转速(r/min):0~50加热功率(kW):6

配套电机(kW):

1.47

外形尺寸(mm):

 $800 \times 1300 \times 1650$

其它规格:TNW200 TNW300 TNW400 TNW500 TNW1000 TNW1000II

COMMODITY: HIGH EFFICIENCY COATING MACHINE

MODEL: BGB-150

APPLICATIONS:

This equipment is used for coating sugar, organic film or aqueous film of Chinese herbal medicines and Western tablets in pharmaceutical industry.

FEATURES:

The pre-loaded tablets run in the coating pan and coating medium are sprayed on the tablet surface. Hot air is supplied for drying, so the even and sturdy protective film can be formed on the tablet surface. The coating process is fully hermetic and automatic. There are three spraying systems in order to meet the different requirements. Atomization air and air for blowing the atomization liquid into flat shape for the spraying gun can be separately adjusted. Specially designed structure makes the gun avoid blockage.

TECHNICAL DATA:

| Productivity (kg/batch): | 150 |
|--------------------------------------|-------------------------------------|
| Drum speed (r/min): | 3~18 |
| Hot air temperature (℃): | from indoor temperature to 80 |
| Steam consumption (kg/h): | 100 |
| Air consumption (m ³ /h); | 90 |
| POWER (kW): | 10 |
| OVERALL DIMENSIONS (mn | 1,680 \times 1,260 \times 2,000 |

OTHER SPECIFICATION:

BGB-10,BGB-75

产品名称:高效包衣机

型号规格:BGB-150

主要用途:用于医药行业为中、西药片包敷糖衣、有机薄膜或水溶薄膜。

工作原理及结构特点:片芯在包衣机旋转滚筒内作复杂的轨迹运动,包敷介质喷洒在片芯上,并供热风进行干燥,使片芯表面快速形成坚固均匀的保护膜。包衣过程全密闭,自动化。配有三种喷雾系统,可满足包制糖衣、肠溶衣或胃溶衣的需要。整机采用防爆结构。

主要技术参数:

生产能力(kg/批):

150

滚筒转速(r/min):

 $3 \sim 18$

热风温度(C):

室温~80

蒸汽耗量(kg/h):

100

空气耗量(m³/h):

90

配套电机(kW):

10

外形尺寸(mm):

 $1680 \times 1260 \times 2000$

其它规格:BGB-10,BGB-75

CHINA PHARMACEUTICAL EQUIPMENT 89

COMMODITY: HIGH-EFFICIENCY SUGAR COATING MACHINE

MODEL: BGB-75B

APPLICATIONS:

It is used for coating sugar, organic film and aqueous film of Chinese herbal medicines and Western tablets in pharmaceutical, foodstuff industry.

FEATURES:

The pre-loaded tablets run in the coating pan and coating medium is sprayed on the tablet surface. Hot air is supplied for drying, so the even and sturdy protective film can be formed on the tablet surface. It is controlled by PLC and has multi-option for choose. It has auto input and discharge unit.

TECHNICAL DATA:

| Production capacity (kg/times): | 75 |
|------------------------------------|-----------------------------------|
| Drum speed (rpm): | 4~20 |
| Hot air temperature ('C): | 0~80 |
| Filtering accuracy of hot air((m): | 0.5 |
| | |
| POWER (kW): | 1.1 |
| OVERALL DIMENSIONS (mm): | $1,320 \times 1,013 \times 1,630$ |

OTHER SPECIFICATION:

BGB-10 BGB-75A BGB-150A BGB-150B

产品名称:高效包衣机

型号规格:BGB-75B

主要用途:用于制药、食品行业片剂、丸剂、糖果包制糖衣、薄膜。

工作原理及结构特点:片芯在滚筒中作复杂的轨迹运动, 在负压状态下热风穿过片芯层,从底部排出。喷洒在 片芯上的包衣介质快速干燥,形成坚固光滑的表面 薄膜。该机工作过程由微机控制,备有多种程序可供 选择,并配有自动进出料装置。

主要技术参数:

生产能力(kg/次): 75 滚筒调速范围(r/min): 4~20 调温范围(℃): 0~80 空气过滤精度(μm): 0.5

配套电机(kW):

1.1

外形尺寸(mm):

 $1320 \times 1013 \times 1630$

其它规格:BGB—10 BGB—75A BGB—150A BGB—150B

CHINA PHARMACEUTICAL EQUIPMENT 91

COMMODITY: TABLET DEDUSTING UNIT

MODEL: CFQ-I

APPLICATIONS:

It is used for removing dust from tablet, piece and granule in pharmaceutical, chemical and food industries, etc.

FEATURES:

It is the necessary accessories of the high-speed tablet press, and can be used with the vacuum cleaner.

TECHNICAL DATA:

Output efficiency (tablets/h):

 $550,000(\emptyset 6 \text{ tablet})$

POWER (kW):

OVERALL DIMENSIONS (mm):

 \emptyset 320 \times 950

OTHER SPECIFICATION:

产品名称:除粉器

型号规格:CFQ-1

- 主要用途:用于医药、化工、食品等行业片、块、粒状产品 除尘。
- 工作原理及结构特点:本机是高速压片机配套设备,可与 气动吸尘装置配套使用。

主要技术参数:

输出效率(万片/h): 55(Ø6片)速度连续可调

配套电机(kW):

外形尺寸(mm):

 \emptyset 320 \times 950

CHINA PHARMACEUTICAL EQUIPMENT 93

COMMODITY: AMPOULE CLEANING-DRYING AND FILLING-SEALING MACHINE

MODEL: BXSZ1/20

APPLICATIONS:

It is applicable for producing liquid injection ampoules in the pharmaceutical industry.

FEATURES:

This set can complete such procedures as ultrasonic cleaning, multi-nozzle water and steam washing alternately, hot air LF sterilizing, LF purification, multi-nozzle filling, wire drawing and sealing etc.

TECHNICAL DATA:

Applicable size (ml):

 $1 \sim 20$

Production capacity (pieces/h):

1,2ml 16,000

POWER (kW):

44

OVERALL DIMENSIONS (mm): $7,500\times2,100\times2,600$

产品名称:安瓿洗、烘、灌、封联动机组

型号规格:BXSZ 1/20

主要用途:用于医药行业生产安瓿水针剂。

工作原理及结构特点:本机组可完成超声波清洗、多针水气交替冲洗、热层流消毒、层流净化、多针灌装和拉丝封口等工序。

主要技术参数:

适用规格(ml):

 $1 \sim 20$

生产能力(支/h):

1,2ml 16000

配套电机(kW):

44

外形尺寸(mm):

 $7500 \times 2100 \times 2600$

CHINA PHARMACEUTICAL EQUIPMENT 95

COMMODITY: AMPOULE WASHING, DRYING & SEALING PRODUCTION LINE

MODEL: ACSD 1-20

APPLICATIONS:

It is applicable for producing liquid injection of $1\sim20$ ml ampoules in the pharmaceutical industry.

FEATURES:

This line consists of three separate machines, such as QCA series ultrasonic ampoule washer, SZA 400/18 far-infrared or SZA 400/32 hot air circular sterilizer, and DGA 6//1-20 ampoule filling and sealing machine. It can carry out the processes of water rinsing, cleaning with ultrasonic and water vapor alternatively, drying, sterilizing, cooling, nitrogen charging, liquid filling and sealing. It adopts the laminar flow protection and frequency conversion and PLC as control techniques..

TECHNICAL DATA:

| Production capacity (pcs. /h): | 1~2 ml 15 ⋅ | ,000~16,000 |
|--------------------------------|-------------|-------------|
| | 5ml | 12,000 |
| | 10ml | 8,000 |
| | 20ml | 4,500 |

POWER (kW):

OVERALL DIMENSIONS (mm): $3,500 \times 2,300 \times 2,400$

OTHER SPECIFICATION:

产品名称:安瓿洗、烘、灌、封联动机

型号规格:ACSD 1-20

主要用途:用于医药行业生产1~20ml 安瓿水针剂。

工作原理及结构特点:本机组由 QCA 系列安瓿超声波清洗机、SZA400/18 型远红外或 SZA400/32 型热风循环杀菌干燥机、DGA6/1—20 型安瓿灌装封口机三台单机组成。可完成喷淋灌水、超声波及水气交替冲洗、干燥灭菌、冷却、充氮、灌液与封口等工序。设有百级层流保护,采用变频、PLC 等控制技术。

主要技术参数:

| 生产能力(支/h): | 1,2ml: | 15000~16000 |
|------------|--------|-------------|
| | 5ml: | 12000 |
| | 10ml: | 8000 |
| | 20ml: | 4500 |

配套电机(kW):

外形尺寸(mm):

 $8500 \times 2300 \times 2400$

CHINA PHARMACEUTICAL EQUIPMENT 97

COMMODITY: COMPUTERIZED AUGER FILLING MA-CHINE

MODEL: DLF 5 PC

APPLICATIONS:

It is used for filling crystalline powder, spray-drying powder, freeze-drying powder and powder with poor fluidity.

FEATURES:

The filling machine, which is fully computerized and applies Man-Machine talk and screen display, has eight menus for easy selection. It is easy to set and amend the parameters of the step motor and to choose the opposite rotation of the motor as necessary. If the augar touches the surface of the funnel, the computer will show A or B and alarm. There is no filling while no vial on the filling station.

TECHNICAL DATA:

| Capacity: | 0.06~2.0g/vial | 120vial/min |
|------------|--------------------------|---------------|
| | 2. $0\sim3$. 0 g/vial | 70~80vial/min |
| | 5. 0g/vial | 30~40vial/min |
| Tolerance: | 0.5∼2.0g | ±2~4% |
| | 150 g | ±3~6% |

POWER (kW):

OVERALL DIMENSIONS (mm): $1,450 \times 900 \times 1,700$

OTHER SPECIFICATION:

产品名称:电脑自控螺杆分装机

型号规格:DLF 5PC

- **主要用途:**用于结晶粉、喷雾干燥粉、冻干粉及流动性差的药粉的分装。
- 工作原理及结构特点:本机全部采用电脑控制,实现人机对话,屏幕显示设有8个菜单可任意选用、任意设定与修改步进电机运行参数和电机是否加反转,任何一个螺杆"碰壳"A/B,电脑立即会显示和报警。加有工位无瓶则相对应的分装头不下粉的功能。

主要技术参数:

| 生产能力: | 0.06~2.0 克装 | 120 瓶/分 |
|-------|-------------|-----------|
| | 2.0~3.0 克装 | 70~80 瓶/分 |
| | 5.0 克装 | 30~40 瓶/分 |
| 装量误差: | 0.5~2.0 克装 | ±2~4% |
| | 150mg | ±3~6% |

配套电机(kW):

2kW,380VAC,50Hz

 $1450 \times 900 \times 1700$

外形尺寸(mm):

CHINA PHARMACEUTICAL EQUIPMENT 99

COMMODITY: COMPUTER-CONTROLLED SCREW

TYPE POWDER BOTTLING MACHINE

MODEL: FZL2-120

APPLICATIONS:

It is used in pharmaceutical factories for automatic filling powder into bottles of diameters $\varnothing 22 \sim 33$ and sealing.

FEATURES:

Advanced micro computer program-controlled is adopted in this unit. It is equipment with such functions as , stopping filling while there is no bottle, and detecting when touching with bottle wall. The speed of main motor is regulated by frequency inverter. Some parts which contact powder are made of stainless steel. It is easy to operate and repair, and is safe and reliable.

TECHNICAL DATA:

Production capacity(bottles/min): 120

Measurement precision: $\pm 1\% \sim 3\%$

Sealing types: Automatic rolling sealing

Sealing ratio: 99.5%

POWER (kW): 0.8

OVERALL DIMENSIONS (mm): $1,800 \times 1,200 \times 1,550$

OTHER SPECIFICATION:

产品名称:微机数控螺杆粉剂灌装机

型号规格:FZL 2-120

- 主要用途:用于制药厂生产直径 22~32mm 西林瓶的粉 剂自动灌装及封口。
- 工作原理及结构特点:本机采用先进的微机程控,参数记忆断电不丢失,无瓶止灌,碰壳检测等功能。采用步进电机调速,与药粉接触的部分全部采用不锈钢件,易于操作及维护,使用安全可靠。

主要技术参数:

生产能力(瓶/分):

120

计量精度(%):

 $\pm 1\sim 3$

落塞方式:

机械自动滚塞

带塞率(%):

99.5

配套电机(kW):

0.8

外形尺寸(mm):

 $1800 \times 1200 \times 1550$

COMMODITY: ANTIBIOTIC POWDER PRODUCTION LINE

MODEL: HFZ-300

APPLICATIONS:

It can be used for filling and packaging antibiotic powder.

FEATURES:

This production line consists of an ultrasonic vial cleaning machine, sterilizing tunnels, drying box, powder filling machine, capping machine, labeling machine, semi-product inspection belts and four damper trays. The cleaning, drying, sterilizing, powder filling, rubber stopping, capping, inspection, labeling of vials and finished products can all be accomplished automatically.

TECHNICAL DATA:

| Production capacity (vials/min): | 280~3 00 |
|--|-----------------|
| Charge of powder (mg): | 250~500 |
| Filling tolerance(%): | ±5 |
| Sterilizing temperature (°C): | 350±7 |
| Broken rate of vials (%): | 0. 2 |
| Missing probability of rubber plug(%): | 0.3 |
| Defective for labeling(%): | 1 |

POWER (kW).

OVERALL DIMENSIONS (mm): $23,040 \times 3,000 \times 3,100$

OTHER SPECIFICATION:

产品名称:抗生素粉针生产联动线

型号规格:HFZ-300

主要用途:用于抗生素粉剂的充填、包装。

工作原理及结构特点:全线由超声波洗瓶机、隧道消毒烘干箱、分装机、轧盖机、贴签机及半成品检查台和四个缓冲转盘组成。自动完成洗瓶、烘干消毒、装粉、上胶塞、加铝盖、检查、贴签出成品全部过程。

主要技术参数:

| 生产能力(瓶/min): | 280~300 |
|--------------|---------|
| 装药量(mg): | 250~500 |
| 装量误差(%): | ±5 |
| 灭菌温度(℃): | 350±7 |
| 玻璃瓶破损率(%): | 0.2 |
| 盖胶塞失误率(%): | 0.3 |
| 贴标签不合格率(%): | 1 |
| | |

 $23040 \times 3000 \times 3100$

其它规格:

配套电机(kW):

外形尺寸(mm):

COMMODITY: RUBBER STOPPER WASHING STERIL-IZATION-DRYING MACHINE

MODEL: JS-140

APPLICATIONS:

It is used for washing, sterilizing, drying, and meets the GMP standard.

FEATURES:

The rubber stoppers are sucked into the vessel by vacuum, then they are washed by clean water and compressed air and be dried through the high temperature sterilizer. All rubber stoppers are sterilized as the whole procedure is finished in a hermetic vessel.

TECHNICAL DATA:

Vessel diameter (mm):

1,400

Treatment quantities (piece/batch):

180,000(stopper for antibiotics)

POWER (kW):

11

OVERALL DIMENSIONS (mm): $3,700\times1,650\times3,200$

OTHER SPECIFICATION:

JS-90 JS-100

产品名称:胶塞洗涤--灭菌--干燥联合机组

型号规格:JS-140

- 主要用途:用于医药、化工、食品等行业胶塞的洗涤、灭菌和干燥。
- 工作原理及结构特点:通过真空将胶塞抽入转筒,经压缩空气和水自动清洗,蒸汽高温消毒灭菌,热空气干燥,胶塞得到净化。全过程在一个封闭的容器内完成。

主要技术参数:

容器直径(mm):

1400

处理量(只/批):

180000(抗生素瓶胶塞)

配套电机(kW):

11

外形尺寸(mm):

 $3700 \times 1650 \times 3200$

其它规格:JS-90 JS-100

COMMODITY: PROGRAM CONTROL RUBBER STOPPER WASHER

MODEL: MJK-80

APPLICATIONS:

It is suitable for pharmaceutical, chemical, foodstuffs industries etc. to wash, sterilize and dry the rubber stopper.

FEATURES:

This machine is equipped with PLC and can automatically finish all process of washing, sterilizing and drying. Rubber stopper are loaded in separate small vessel and can be washed circularly. It has such functions: as displaying working state and trouble alarming, Chinese prompting and report printing etc.

TECHNICAL DATA:

Capacity(pcs/batch):

30,000

Stopper to be used:

A2 rubber stopper

POWER (kW):

25, 51

OVERALL DIMENSIONS (mm):

 $2,000 \times 2,500 \times 2,500$

产品名称:程控胶塞清洗机

型号规格:MIK-80

- **主要用途:**用于医药、化工、食品等行业胶塞的洗涤、灭菌和干燥。
- 工作原理及结构特点:本机采用可编程控制器,自动完成各种胶塞的清洗、硅化、消毒、干燥全过程、胶塞分仓装入、循环清洗。该机并具有工况显示、报警、中文揭示及报表打印等功能。

主要技术参数:

生产效率(只/次):

采用的胶塞:

30000(A2胶塞)

A。橡胶塞

配套电机(kW):

25.51

外形尺寸(mm):

 $2000 \times 2500 \times 2500$

COMMODITY: MULTI FUNCTIONAL FULLY-AUTO-MATIC ULTRASONIC RUBBER-CORK WASHER

MODEL: DQX-1B

APPLICATIONS:

As a specially-used equipment in medical line, the washer is suitable to wash, bake dry and sterilize medical infusion and oral-liquid rubber corks and Al-lid.

FEATURES:

- 1. The washer combines the multi-procedures originally used to process rubber corks as an one operation of washing, siliconizing, baking dry and sterilizing, uses a new-style programmable controller to control and completes all procedures once inside of it with a fully-sealed operation.
- 2. Washing with ultrasonic, further spraying by high pressure water.
- 3. Baking dry and sterilizing; to bake dry and sterilize under grade 100 of cleanness.

TECHNICAL DATA:

Capacity (pcs/box):

10,000~80,000

Time for washing (min/box): 10~20 (natural rubber cork)
10 (butyl rubber cork)

Time for bake dry and sterilizing (h/box): 2.5

Temperature for sterilizing(°C): 125(5 with circular wind

POWER (kW): 28

OVERALL DIMENSIONS (mm): $2,300 \times 1,620 \times 2,100$

OTHER SPECIFICATION:

产品名称:多功能全自动超声波胶塞清洗机

型号规格:DQX-1B

- 主要用途:本机是医药行业胶塞及铝盖清洗的专用设备, 主要用于清洗、烘干、消毒、药用胶塞、输液胶塞、口 服液胶塞及铝盖。
- 工作原理及结构特点:该机将原来的多道工序:洗胶塞、硅化、预烘、干燥及灭菌合为一体。采用新型可编程控制器,全封闭式操作,在机内一次完成全部工序。 先用超声波清洗,再用高压水进一步喷淋清洗。烘干及灭菌:在100级洁净度下进行烘干、灭菌。

主要技术参数:

生产能力(支/箱):

10000~80000

洗涤时间(分钟/箱):

10~20(天然胶塞)

10(丁基胶塞)

烘干、灭菌时间(小时/箱):

2. 5

灭菌温度(℃):

循环风 125±5

配套电机(kW):

28

外形尺寸(mm):

 $2300 \times 1620 \times 2100$

COMMODITY: WATER-BATH TYPE STERILIZER

MODEL: SR 2540/90

APPLICATIONS:

It is used for glass bottles, plastic bottles (bag) of infusion and injection in pharmaceutical industry.

FEATURES:

The facility consists of cabinet, hot water circulation pump, heat exchanger and computer-control cabinet. By using industrial computer, the circulation water in sterilizing cabinet pour out from top to down after being exchanged in heat exchanger. Sterilizing process of temperature going up, stabilizing and cooling down are completed automatically. The whole procedure is automatically controlled, color screen displays the real-time values of temperature, pressure, sterilizing time and default detection and the main parameters of sterilizing process are printed automatically.

TECHNICAL DATA:

| Production capacity (bottles/cabinet): | 2,540 (500ml) |
|--|-------------------|
| Nominal pressure (MPa): | 0.3 |
| Working temperature (°C): | 115 |
| Steam pressure (MPa): | 0.4~0.7 |
| Volume (m³): | 5.48 |
| POWER (kW): | 12 |
| OVERALL DIMENSIONS (mm): | 4,600×4,300×2,580 |

OTHER SPECIFICATION:

SR-600/90 SR-1200/90 SR-3000/90 SR-4000/90 SR-5000/90 SR-6000/90

产品名称:水浴式灭菌柜

型号规格:SR 2540/90

主要用途:用于医药行业玻璃瓶输液剂和水针剂的灭菌。

工作原理及结构特点:本设备由柜体、热水循环泵、热交换器和微机控制柜组成。微机控制灭菌柜的循环水通过热交换器加热、恒温、冷却、自上向下喷淋,完成灭菌,冷却出瓶过程。全过程自动控制,具有温度、压力数字显示、超限自动报警、F。值自动打印等功能。

主要技术参数:

| 生产能力(瓶/90min): | 2540(500ml) |
|----------------|-------------|
| 设计压力(MPa): | 0.3 |
| 工作温度(℃): | 115 |
| 蒸汽压力(MPa): | 0.4~0.7 |
| 设备容积(m³): | 5, 48 |

配套电机(kW):

12

外形尺寸(mm):

 $4600 \times 4300 \times 2580$

其它规格:SR-600/90 SR-1209/90 SR-6000/90 SR-3000/90 SR-4000/90 SR-5000/90

COMMODITY: QUICK COOLING STERILIZER

MODEL: KLM-B

APPLICATIONS:

It is used for sterilizing of infusion bottles and glass bottles in pharmaceutical industry.

FEATURES:

By using gradually cooled water to cool infusion bottles which have been sterilized, it makes the temperature of infusion bottles decline gradually and also avoids them being broken due to temperature suddenly declined. It is auto-controlled by PLC and its program can be preset.

TECHNICAL DATA:

| Pressure (MPa): | 0. 2 |
|------------------------------|---------------------------------|
| Temperature (°C): | 110~121 |
| Steam source pressure (MPa): | 0.25~0.69 |
| Water source pressure (MPa): | 0.15~0.3 |
| Volume (mm): | $1,500 \times 760 \times 1,100$ |
| Displacement time (min): | 5 |
| Sterilizing time (s): | 0~9,999 |
| | |

| POWER (kW); | 0.5, 0.75 |
|--------------------------|-------------------|
| OVERALL DIMENSIONS (mm): | 2,400×1,700×2,200 |

OTHER SPECIFICATION:

XG1 K KLM-C

产品名称:快速冷却灭菌器

型号规格:KLM-B

主要用途:用于医药行业对玻璃瓶输液剂的灭菌处理。

工作原理及结构特点:利用逐渐降温的冷却水对灭菌后的大输液瓶喷淋冷却,使之逐渐降温,避免了大输液瓶因温度骤降而造成的爆瓶现象。PLC自动控制,程序予置。

主要技术参数:

工作压力(MPa):

0.2

工作温度(℃):

110~121

汽源压力(MPa):

 $0.25 \sim 0.69$

水源压力(MPa):

 $0.15 \sim 0.3$

灭菌容积(mm):

 $1500 \times 760 \times 1100$

灭菌时间(s):

0~9999

置换时间(min):

5

配套电机(kW):

0.5,0.75

外形尺寸(mm):

 $2400 \times 1200 \times 2200$

其它规格:XG1 K KLM-C

COMMODITY: BOTTLE INFUSION PRODUCTION LINE

MODEL: PSY 100/500B

APPLICATIONS:

It is used in production for 100ml, 250ml and 500ml infusion bottles in pharmaceutical industry.

FEATURES:

Except for bottle-transferring machine, bottle-filling machine (linear inject type), plug-turning machine, and cap pressing machine are B-type. In this production line the other necessary machines are as same as A-type line. This line raise the speed of bottle transferring machine and be equipped with filter unit, nitrogen-filled unit, synchronizing filling device, automatic alarm and protection devices, etc.

TECHNICAL DATA:

| Production capacity (bottles/h): | 1,200~3,000 |
|-------------------------------------|--------------------|
| Filling error (%): | <2 |
| Number of nitrogen-filling feeders: | 8 |
| Number of bottle-filling feeders: | 8 |
| | |
| POWER (kW); | 13 |
| OVERALL DIMENSIONS (mm): | 28,000×3,400×2,000 |

OTHER SPECIFICATION:

PSY-100/500C PSY 100/500A PSY-100/500D

产品名称:玻璃瓶输液生产联动线

型号规格:PSY 100/500B

主要用途:用于医药行业生产 250ml 和 500ml 玻璃瓶输 液。

工作原理及结构特点:生产线输瓶机、灌装机(直线式)翻 盖机、轧盖机为B型机外,其余配套单机与A型联动 线相同,该线提高了输瓶机的速度,并设有终端过 滤、充氮、自动报警、停机保护等装置。

主要技术参数:

生产能力(瓶/小时): $1200 \sim 3000$ 灌装误差(%): ≤ 2 充氮头数(个): 8 灌装头(个): 8

配套电机(kW):

13

外形尺寸(mm): 28000×3400×2000(不含贴签机)

其它规格: PSY100/500C PSY 100/500A PSY 100/ 500D

COMMODITY: SEMI-AUTOMATIC CAPSULE FILLING MACHINE

MODEL: JTJ-II

APPLICATIONS:

It is suitable for pharmaceutical industry to fill hard capsules with powder and granule.

FEATURES:

The feeding, separation, filling and sleeving of capsules and the discharging of finished products can be accomplished automatically by mechanical driving.

TECHNICAL DATA:

| Capacity(capsules/h): | 6,000-14,000 |
|--------------------------------|--------------------|
| Capsules size: | 00",0",1",2",3",4" |
| Error of weight dosage(%): | less than ± 5 |
| Rate of finished products(%):: | more than 97 |

POWER (kW): 3.15 OVERALL DIMENSIONS (mm): $1,400 \times 950 \times 1,480$ 产品名称:半自动胶囊充填机

型号规格:JTJ-I

主要用途:适用于制药行业生产硬胶囊剂。

工作原理及结构特点:采用电器、气动联合控制,配备电子自动计数装置,自动完成播囊、分离、充填药料、套合、并输出成品。

主要技术参数:

生产能力(粒/h):

6000~14000

合格率(%):

97

适用胶囊:

 $00^{#}, 0^{#}, 1^{#}, 2^{#}, 3^{#}$

药量充填误差(%):

 $<\pm 5$

配套电机(kW):

3.15

外形尺寸(mm):

 $1400 \times 950 \times 1480$

COMMODITY: SEMI-AUTOMATIC CAPSULE FILLING MACHINE

MODEL: NCB-200B

APPLICATION:

Used for filling hard capsules in the pharmaceutical industry.

FEATURES:

The capsule- arranging, setting, separating, filling, sleeving and other processes of capsule-filling are carried out by mechanical drive.

TECHNICAL DATA:

| Production capacity(pcs/h): | 12,000 |
|-----------------------------|-----------------|
| Form of drug filled: | powder, granule |
| Specification of capsule: | 0",1",2",3" |

POWER(kW): 1.1 **OVERALL DIMENSIONS(mm):** 1,380×950×1,410 产品名称:半自动胶囊充填机

型号规格:NCB-200B

主要用途:用于制药行业生产硬胶囊剂。

工作原理及结构特点:通过机械传动,完成胶囊的排囊、 就位、分离、充填、锁紧等过程。

主要技术参数:

生产能力(粒/h):

12000

充填剂型:

粉剂、小颗粒

胶囊:

0 " 1 " 2 " 3 "

配套电机(kW):

1.1

外形尺寸(mm):

 $1380 \times 950 \times 1410$

COMMODITY: SEMI AUTOMATIC CAPSULE FILLING MACHINE

MODEL: JTJ-III

APPLICATIONS:

This machine meets the requirement of medicine and hygiene, which can fill capsules home-made or imported from abroad. It is the ideal equipment for pharmaceutical manufacturers to fill capsules till now.

FEATURES:

This machine combining electrical and pneumatic control, fitted with AC frequency conversion stepless speeder and equipped with an electronic automatic counter and touch-type button, is reliable to use, sensitive to control, and easy to operate and automatically performs such actions as: capsule positioning, separating, filling and locking respectively, so the finished products rate of 97% or more is obtained.

TECHNICAL DATA:

| Production capacity(granule/h): | 6,000~13,000 |
|---------------------------------|---------------------------------|
| Capsule size: | 0",1",2",3" |
| Filling types & forms: | powder, pellet |
| POWER (kW): | 2. 61 |
| OVERALL DIMENSIONS (mm): | $1.240 \times 860 \times 1.480$ |

OTHER SPECIFICATIONS:

产品名称:半自动胶囊充填机

型号规格:JTJ-I

主要用途:适用于国产和进口胶囊的充填。

工作原理及结构特点:本机采用电动、气动联合动作和交流变频无级变速,并配备电子计数装置,触摸式按键,动作灵敏、可靠,操作方便,能分别自动完成胶囊的就位、分离、充填、锁紧等动作。充填剂量准确,充填合格率高于97%,符合制药卫生要求。

主要技术参数:

生产能力(粒/h):

 $6000 \sim 13000$

胶囊:

 $0^{\#}, 1^{\#}, 2^{\#}, 3^{\#}$

充填剂型:

粉剂、小颗粒剂

配套电机(kW):

2.61

外形尺寸(mm):

 $1240 \times 860 \times 1480$

COMMODITY: HARD CAPSULE FILLING MACHINE

MODEL: YJF 15

APPLICATIONS:

It can be used to fill the Chinese herbal medicines and western powdery medicines or granules into capsules.

FEATURES:

The feeding, separation, filling and sleeving of capsules and discharging of finished products can be accomplished automatically by mechanical driving.

TECHNICAL DATA:

Production capacity (capsules/min):

1.500

Rate of feed (%):

98

Size of capsule:

0"1"2"3"4"5"00"

POWER (kW);

6.7

OVERALL DIMENSIONS (mm): $1,200\times1,050\times2,100$

OTHER SPECIFICATION:

产品名称:硬胶囊分装机

型号规格:YJF 15

主要用途:用于医药行业将粉、粒状物料充填入胶囊。

工作原理及结构特点:通过机械传动,自动完成排囊、分离、充填药料、套合、输出成品等过程。

主要技术参数:

生产能力(粒/min):

1500

胶囊上机率(%):

98

胶囊规格:

 $0^{#}, 1^{#}, 2^{#}, 3^{#}, 4^{#}, 5^{#}$

配套电机(kW):

6.7

外形尺寸(mm):

 $1200 \times 1050 \times 2100$

COMMODITY: FULLY AUTOMATIC CAPSULE FILLING MACHINE

MODEL: ZTJ-400

APPLICATIONS:

It is used for filling the capsules with powder or drug granule in pharmaceutical industry.

FEATURES:

The main shaft is driven by a stepless transmission to make the whole devices operate intermittently, and automatically complete vacuum separating the heads and hoods of the capsules, pharmaceutical material filling and product conveying, etc. It can control filling dosage, accumulative total output counting, accumulative total time counting, and reject the unqualified capsules in time.

TECHNICAL DATA:

| Production capacity (capsules/min): | 400 |
|-------------------------------------|----------------|
| Capsule acceptability (%): | >98 |
| Capsule specification: | 0" 1" 2" 3" 4" |
| Noise (dB): | <85 |

POWER (kW): 2.5 **OVERALL DIMENSIONS (mm):** 790×960×1,800

OTHER SPECIFICATION:

ZTJ-700 ZTJ-800A

产品名称:全自动胶囊充填机

型号规格:ZTJ-400

主要用途:用于医药行业将粉、粒状物料充填入胶囊。

工作原理及结构特点:通过无级变速器驱动主轴,使各机构作间歇运动,自动完成真空分离胶囊体、帽、药物充填、套合、输出成品等。可自动控制充填剂量、产量累积计数、累积计时及排除废胶囊。

主要技术参数:

生产能力(粒/min):

400

胶囊上机率(%):

> 98

胶囊规格:

0",1",2",3",4"

噪音(dB):

<85

配套电机(kW):

2.5

外形尺寸(mm):

 $790 \times 960 \times 1800$

其它规格:ZTJ-700 ZTJ-800A

COMMODITY: AUTOMATIC HARD CAPSULE FILLING MACHINE

MODEL: CFM 600

APPLICATION:

It is used for filling powder/pellet into the empty hard gelatin capsules.

FEATURES:

The filling machine is running intermittently. All of operation such as separation, powder/pellet dosing, capsule joining and capsule rejection are worked automatically. The double leak-free seal is used to prevent main shaft from wear and being damaged. The frequency control and slow loading are used to improve the life of the machine.

TECHNICAL DATA:

| Capacity (pcs/min): | 600 |
|------------------------|--------|
| Filling percentage(%): | > 98 |
| Capsule size: | 00"~5" |
| Noise(dB): | < 85 |

| POWER(kW): | 2.5 |
|-------------------------|-------------------------------|
| OVERALL DIMENSIONS(mm): | $700 \times 900 \times 1.800$ |

OTHER SPECIFICATIONS:

CFM 400, CFM 800

产品名称:全自动胶囊充填机

型号规格:CFM 600

主要用途:用于医药行业将粉、粒状物料充填入空胶囊。

工作原理及结构特点:本机采用蜗轮间歇分度机构,所有的生产过程,如体帽分离、粉、粒灌装充填、体帽锁合及成品推出均实现全自动化。充填机构具有双层防渗漏密封圈,有效地防止了主轴磨损过快、抱死现象的出现,整机运转采用无级变频调速及慢加负载控制,提高了机器使用寿命。

主要技术参数:

生产能力(Pcs/min):

600

充填精度(%):

>98

可充填胶囊规格:

 $00^{\#}\sim5^{\#}$

噪音(dB):

<85

配套电机(kW):

2.5

外形尺寸(mm):

 $700 \times 900 \times 1800$

其它规格:CFM 400 CFM 800

COMMODITY: CAPSULE POLISHING MACHINE

MODEL: JMJ-1

APPLICATIONS:

It is used for cleaning dust which adheres to the capsules.

FEATURES:

Capsules enter the cylinder of machine and pass through the polish shaft and brush, which have a helix motion along the guide spring. Dust will be brushed down and sucked by the cleaner. DC serve motor is fitted with the machine. The velocity of rotation can be adjusted.

TECHNICAL DATA:

Capacity(capsules/min):

3,000.

Capsule size:

00" 1" 2" 3" 4"

POWER (kW):

0.185

OVERALL DIMENSIONS (mm):

 $1,310\times405\times1,200$

OTHER SPECIFICATION:

产品名称:胶囊磨光机

型号规格:JMJ-1

主要用途:用于清除附着在胶囊外壳上的粉尘。

工作原理及结构特点:胶囊经磨光轴毛刷的带动进入简体,沿着导向弹簧作螺旋运动,粉尘被刷掉后由吸尘机吸走,从而达到磨光的目的。

主要技术参数:

生产能力(粒/min):

3000

适用范围:

00",1",2",3",4"

配套电机(kW):

0.185

外形尺寸(mm):

 $1310 \times 405 \times 1200$

COMMODITY: SOFT GELATIN ENCAPSULATOR

MODEL: RJNJ-2

APPLICATIONS:

It is suitable for producing medicine, cosmetics and food, etc.

FEATURES:

It is an automatic production line for soft capsule manufacturing, which comprises a main machine, a capsule conveyer, tumbler dryers, an electric control panel, a gelatin tank and material tank, etc. By adopting a frequency inverter, the die rollers of the main machine can be regulated form 0 to 5 rpm steplessly. This machine has the advantage of high level of automation and high productivity, and is accurate in dosage filling.

TECHNICAL DATA:

Rotating speed of die rollers(rpm):

0~5 (stepless speed regulation)

Dimension of die rollers (mm):

Ø103×152

Dosage filling tolerance(%):

士2

POWER (kW);

7

OVERALL DIMENSIONS (mm):

 $880 \times 640 \times 1.900$

OTHER SPECIFICATION:

产品名称:软胶囊机

型号规格:RJNJ-2

主要用途:用于医药、化妆品、食品等行业生产软胶囊。

工作原理及结构特点:本机由主机、输送机、干燥机、电控箱等组成。可以将油类、液体、混悬液或糊状物定量喷注于两条连续生成的明胶带之间,经两个圆筒状成型模具滚压、切断而成球形、圆柱形、橄榄形等形状和大小各异的密封软胶囊。

主要技术参数:

滚模转速(r/min):

0~5(变频无级调速)

滚模尺寸(mm):

 \emptyset 103×152

装量差异(%):

+2

配套电机(kW):

7

外形尺寸(mm):

 $880 \times 640 \times 1900$

COMMODITY: AUTOMATIC SEAMLESS ENCAPSULATOR

MODEL: DWJ-1

APPLICATIONS:

It is widely used for producing the seamless soft capsules in pharmaceutical, foodstuffs and cosmetics industries.

FEATURES:

The encapsulator adopts the advanced pulsation constriction technology, eliminating shortcomings of the three rod pump encapsulator. When the material solution supplied by a high precision pump and the self-concentric liquid column at the dropper, intermittent pulsation is applied to make the liquid column to break up into evenly dosed droplets, which automatically form spherical gelatin capsules under the action of surface of the liquid. After cooling, congelation and separation, the capsules are automatically ejected form the machine. The encapsulator has a feature of accurate filling dosage, high production rate and automation.

TECHNICAL DATA:

| Dosage filling of capsule(mg): | 20~500 |
|--------------------------------|---------------------------------|
| Production capacity (pills/h): | 7,200~43,200 |
| POWER (kW): | 1.5/5.0 |
| OVERALL DIMENSIONS (mm): | $1.080 \times 800 \times 2.040$ |

OTHER SPECIFICATION:

产品名称:全自动无缝滴丸机

型号规格:DWJ--I

主要用途:用于医药、食品、化妆品等行业生产无缝软胶丸。

工作原理及结构特点:采用脉冲切割技术,当高精度泵供给的料液与带自控的明胶液在滴头处形成同心组合流柱时,脉动液体"刀"均匀将其切割。在液体表面张力作用下自然形成球状软胶丸。经冷却、凝固、分离被自动排出机外。

主要技术参数:

胶丸装量(mg):

20~500

生产能力(丸/h):

7200~43200

平均/最大耗电量(kW):

1.5/5.0

配套电机(kW):

外形尺寸(mm):

 $1080\times800\times2042$

COMMODITY: GRANULATING AND COATING MACHINE

MODEL: BZJ-1000F

APPLICATIONS:

It can be used for manufacturing of pills and granulating of acetic medicines and slow release medicines, and can also be used in the fields of food, pesticide, detergent and other chemical products.

FEATURES:

This machine has three basic functions: seed beading, granulating and coating. It can produce various globular granules with better globosity, more uniform size and higher qualification rate of percentage.

TECHNICAL DATA:

| Rotor diameter (mm): | 1,000 |
|---------------------------------------|-------------------|
| Min. input of seed beads(kg): | 5 |
| Max. output of globular granules(kg): | 50 |
| Granulation diameter (mm): | 0.25~22 |
| Max. enlargement: | K=2 |
| POWER (kW): | 20 |
| OVERALL DIMENSIONS (mm); | 1,800×1,400×1,400 |

OTHER SPECIFICATION:

产品名称:包衣造粒机

型号规格:BZJ-1000F

主要用途:用于制药行业。完成起母、造粒、包衣及初步干燥等操作。可制造缓释型微胶囊、不同药物分层的组合药粒,也可填充硬胶囊剂。该设备还可用于食品、农药、化妆品、颜料、化工和粉末冶金等领域。

工作原理及结构特点:在旋转的主机转子内,输入一定量的母粒,由于离心力和摩擦力的作用,在定子与转子间的过渡曲面上形成涡旋回转运动的粒子流,按一定规律喷射雾化的浆液和粉状物料,则母粒越滚越大,越滚越圆,成为真球度很高的颗粒。

主要技术参数:

转子直径(mm):

1000

生产能力(kg/批):

50

造粒直径(mm):

 $0.25 \sim 22$

母粒输入量(kg):

5

配套电机(kW):

20

外形尺寸(mm):

 $1800 \times 1400 \times 1400$

COMMODITY: AUTOMATIC PILL MAKING MACHINE

MODEL: ZW-120

APPLICATIONS:

It is suitable for producing honey pills, water pills & condensed pills in pharmaceutical industry.

FEATURES:

The even mixed drug material is put into the taper hopper, after being extruded by the screw pusher, the fine drug stripes with same diameter are formed. The self-controlled guide pulley carries the drug strip to the cutting unit and finishes the cutting procedure continuously.

TECHNICAL DATA:

| Production capacity (kg/h): | 80~120 |
|-----------------------------|--------|
| Diameter of pills (mm): | 4~12 |
| Weight of pill (g): | 3,6,9 |
| | |

POWER (kW): 4.5 **OVERALL DIMENSIONS (mm):** 1,200×850×2,000

OTHER SPECIFICATION:

ZW-10 ZW-20 ZW-80

产品名称:中药自动制丸机

型号规格:ZW-120

主要用途:用于医药行业蜜丸、水丸和浓缩丸的生产。

工作原理及结构特点:将混合均匀的药料投入到锥形料 斗内,在螺旋推进器的挤压下,推出多条相同直径的 药条,在自控导轮的控制下,同步进入制丸刀,连续 制成大小均匀的药丸。

主要技术参数:

生产能力(丸/h):

80~120

制丸直径(mm):

4~12

制丸重量(g):

3,6,9

配套电机(kW):

4.5

外形尺寸(mm):

 $1200 \times 850 \times 2000$

其它规格: ZW-10, ZW-20, ZW-80

COMMODITY: THREE ROLLER HONEY BOLUS MA-CHINE

MODEL: ZTM80-3

APPLICATIONS:

It is used for processing bolus which is made of powdery Chinese herbal medicine and honey in Chinese herbal medicine factory. The weight of bolus is from 3g to 10g.

FEATURES:

After putting lump made of powdered Chinese medicine and honey into screw press, it is squeezed into strip which can be cut in short sticks and they are pushed to the rollers automatically to form the round shape.

TECHNICAL DATA:

| Production capacity(grains/h): | $(3\sim 9g)$ 28,00 $\sim 19,000$ |
|-------------------------------------|----------------------------------|
| Weight accuracy of bolus (%): | ±2 |
| Speed of cutting strip (times/min): | 15 |
| Length of strip (mm/one): | 576 |
| | |
| POWER (kW): | 4 |
| OVERALL DIMENSIONS (mm): | $1,700 \times 700 \times 1,300$ |

产品名称:三辊蜜丸机

型号规格:ZTM 80-3

主要用途:用于中药厂制做 3、6、9 克蜜丸。

工作原理及结构特点:中药粉剂经搅拌合坨后由挤出器 挤出连续药条,并经自动切条、推条、进入模辊制成 圆丸。

主要技术参数:

生产能力(丸/h): (3~9g)2800~19000 丸重允差(%): ±2 切条次数(次/min): 15 切条长度(mm/条): 576

配套电机(kW):

4

外形尺寸(mm):

 $1700 \times 700 \times 1300$

2

COMMODITY: OINTMENT FILLING MACHINE

MODEL: RHG1-20B

APPLICATIONS:

It is used in the field of pharmacy, food and cosmetics for filling such material as paste, ointment and glutinous flow liquid into metallic tube.

FEATURES:

It is equipped with heating tank, pneumatic unit and can complete folding, sealing and lot number printing. It fills as tubes are available and stops filling without tubes. The whole procedure is detected by means of electron.

TECHNICAL DATA:

Temperature keeping power(kW):

| Compressed air (be provided)(L/min): | 200 |
|--------------------------------------|---------------------------------|
| Pressure(kg/cm ²): | 6~8 |
| POWER (kW): | 1.1 |
| OVERALL DIMENSIONS (mm): | $1,200 \times 900 \times 1,800$ |

产品名称:金属软管灌封机

型号规格:RHG.1-20B

- **主要用途**:用于制药、食品、化妆品等行业糊状、膏体及粘稠流体的金属软管灌封。
- 工作原理及结构特点:本机包括加热罐、气动装置并能完成折叠、封口及批号印刷等步骤,有管置于灌装头时可自动灌装,无管时则停止灌装,整个工作程序由电子监控。

主要技术参数:

保持温度用电力(kW):

2

所需供的压缩空气(L/min):

200

压力(kg/cm²):

 $6 \sim 8$

配套电机(kW):

1.1

外形尺寸(mm):

 $1200 \times 900 \times 1800$

COMMODITY: AUTOMATIC ULTRASONIC TUBE FILLER

MODEL: ZGC 40

APPLICATIONS:

It is used in pharmaceutical and chemical industries. It can fill the unguent in a certain scope of viscosity into plastic tube with high quality and high efficiency.

FEATURES:

This machine adopts intermittent rotating working table structure, and it has functions as follows: automatic feeding, identifying and aiming marks, filling, sealing, trimming and discharging the tubes. It also has assisting functions: no tube-no filling, sealing protection, lower pressure protection and safety protection.

TECHNICAL DATA:

| Tube diameter (mm): | Ø15~Ø30 |
|------------------------------------|-----------------|
| Tube length(mm): | 60~140 |
| Dosage filling volume (ml): | 10~200 |
| Max production capacity(tube/min): | 30~40 |
| POWER (kW): | 1. 2 |
| OVERALL DIMENSIONS (mm): | 1,250×980×1,800 |

OTHER SPECIFICATION:

产品名称:全自动超声波软膏灌装机

型号规格:ZGC 40

主要用途:本机主要应用于医药、日化等行业。可实现将 具有一定粘稠度的膏体高效高质地贮于热塑性塑管 的包装形式。

工作原理及结构特点:本机采用间歇回转工作台的结构, 具有自动上管、自动识标对标、自动灌装、自动封尾、 自动修齐以及自动出管等系列功能。另外还具有识 管灌装,保护封尾,欠压保护及安全防护等辅助功 能。

主要技术参数:

| 软管直 径(mm): | Ø15~Ø30 | |
|-------------------|---------|--|
| 软管高度(mm): | 60~140 | |
| 灌装调量(ml): | 10~200 | |
| 最大生产速度(支/分): | 30~40 | |

配套电机(kW):

1.2

外形尺寸(mm):

 $1250 \times 980 \times 1800$

COMMODITY: SUPPOSITORY DRUG FILLING MA-CHINE

MODEL: SJ-II

APPLICATIONS:

It is widely used in pharmaceutical industry for producing suppository drug.

FEATURES:

The well compound drug powder and carrier are melted in a filling barrel which is at constant temperature, then being automatically and quantitatively filled to the packing strip package. Through the condenser, the drug powder and carrier are solidified to take shape. Through the sealing, counting, cutting units, the final suppository drug is finished.

TECHNICAL DATA:

| Production capacity (pcs/h): | 3,600 |
|-----------------------------------|-------------------------------|
| Volume of the barrel (L): | 20 |
| Volume per piece (ml/pcs): | . 0~3.5 |
| Filling tolerance (%): | ± 0.1 |
| Range of cutting pcs (pcs/plank): | 1~7 |
| POWER (kW): | 2 |
| OVERALL DIMENSIONS (mm): | $2,500\times2,600\times1,700$ |

OTHER SPECIFICATION:

产品名称:栓剂灌装机组

型号规格:SJ-I

主要用途:用于医药行业生产栓剂。

工作原理及结构特点:将配好的药粉和基质倒入恒温灌 装机内融合,并自动计量注入包装带,在冷凝室固化 成形,经封口、齐边、计数切断,制成栓剂。

主要技术参数:

生产能力(粒/h):

3600

装料桶容量(1):

20

每粒装量(ml/粒):

 $0 \sim 3.5$

装量误差(%):

 ± 0.1

剪切粒数(粒/板):

 $1\sim7$

配套电机(kW):

2

外形尺寸(mm):

 $2500 \times 2600 \times 1700$

COMMODITY: AUTOMATIC LIQUID-FILLING LINE

MODEL: YZ25/500

APPLICATIONS:

This machine is used for filling liquid materials, as tincture, syrup, electuary, oil, water, wine and general emnulsion liquid and floating liquid in pharmaceutical, foodstuffs and chemical industries.

FEATURES:

It can automatically arrange bottles, convey bottles, meter filling, arrange caps, turn caps, label bottles, print serial numbers etc. All single machines independently drive and the speed is regulated steplessly.

TECHNICAL DATA:

| Capacity(bottles/min): | 20~80 |
|------------------------|----------------|
| Size(ml): | 25~500 |
| Metering error(%): | $\leq \pm 0.1$ |

POWER (kW): 4.77 **OVERALL DIMENSIONS (mm):** 9,900×1,800×2,000

OTHER SPECIFICATION:

产品名称:液体灌装自动线

型号规格:YZ 25/500

主要用途:用于制药、化工、食品等行业灌装酊水、糖浆、 膏汁、油类、水类及一般乳状悬游液。

工作原理及结构特点:生产线由 GCB4A 型液体灌装机、 SN8 型塞内塞机、XG12 型旋盖机及 ZT20/1000 型 转鼓贴标机组成。可自动完成理瓶、输瓶、计量灌装、 理塞、塞塞、理盖、旋盖、贴标、印批号等工序。各单机 独立传动,无级调速。

主要技术参数:

生产能力(瓶/min):

 $20 \sim 80$

装量(ml):

25~500

计量误差(%):

 ± 0.1

配套电机(kW):

4.77

外形尺寸(mm):

 $9900 \times 1800 \times 2000$

FULLY-AUTOMATIC COMMODITY:

PRODUCTION

LINE FOR FILLING ORAL LIQUID

SEALING WITH ROLL-ON CAPS

MODEL: YLG-3

APPLICATION:

It is suited to combine package-producing oral liquid for glass bottle-washing, bottle-sterilizing, bottle-filling and capsealing.

FEATURES:

This production line consists of the XP-3 ultrasonic cleaning machine, the HX-3 hot air sterilizing tunnel, the GS-3 fully automatic volumetric filling machine with sealing mechanism for roll-on caps. All of the processes (washing, sterilizing, filling, sealing) are finished automatically. Advanced photoelectricity control, PLC technique and automatic safe protection are adopted.

TECHNICAL DATA:

Production rate(bottles/h):

8,000

Nominal filling volume(ml):

5~30 (adjusted by stepless regulator)

Volume tolerance (%):

+3

Broken bottle rate(%): ≤ 2 (standard bottles and caps)

POWER(kW):

Max. /Min. 52/32

OVERALL DIMENSIONS (mm):

 $9.500 \times 4.700 \times 1.900$

OTHER SPECIFICATIONS:

产品名称:易拉盖口服液全自动灌装生产线

型号规格:YLG-3

主要用途:用于生产玻璃瓶口服液剂。

工作原理及结构特点:生产线由 XP-3 超声波洗瓶机、HX-3 隧道式烘干消毒箱、GS-3 全自动灌装锁口机三台单机组成,自动完成洗瓶、烘干消毒、灌液、盖盖、锁口、贴签、打印等工序。可配 TQ-2、TQ-4 自动贴标签机完成贴标,隧道及灌液锁口区达到 100级洁净度。

主要技术参数:

生产能力(瓶/h):

8000

灌装量(ml):

 $5 \sim 30$

灌装精度(%):

+3

破损率(%):

≤2

最大功耗(kW):

52/32

配套电机(kW):

外形尺寸(mm):

 $9500 \times 4700 \times 2900$

COMMODITY: LYING ULTRASONIC BOTTLE CLEAN-ING MACHINE

MODEL: CXP 5-10

APPLICATION:

The machine is mainly used to wash the inside and outside of 5-10ml antibiotics bottles.

FEATURES:

The bottles are delivered softly, then infeeded by rotary drum. The bottles are cleaned by ultrasonic firstly then by water for injection (WFI), ions water and purify compressed air exchangely. The machine can be used individually or combined with a production line.

TECHNICAL DATA:

Production capacity (bottles/min): 60~100

Bottle size(ml): 5~10 (antibiotics bottles)

Ultrasonic power consumption(W): 500

Pressure of purified air(MPa): 0.1~0.15

POWER(kW): 0.6 **OVERALL DIMENSIONS(mm):** 1,200×960×1,100

OTHER SPECIFICATIONS:

产品名称:卧式超声波洗瓶机

型号规格:CXP5~10型

主要用途:用于 5ml 至 10ml 抗生素玻璃瓶的内外清洗。

工作原理及结构特点:本机采用轻拨送瓶,转鼓间歇进瓶,用超声波粗洗,后经注射用水,蒸馏水,净化压缩空气交换冲洗完成洗瓶。本机可单机或联线使用。

主要技术参数:

生产能力(瓶/min):

60~100

适用瓶规格:

5~10ml 抗生素瓶

超声波消耗功能(W):

500

空气压力(MPa):

 $0.1 \sim 0.15$

配套电机(kW):

0.6

外形尺寸(mm):

 $1200 \times 960 \times 1100$

COMMODITY: UNIVERSAL COMMINUTOR

MODEL: 30B II

APPLICATION:

It is applicable to the pharmaceutical, chemical, foodstuff and agriculture industries for breaking the dried fragile material.

FEATURES:

Because of the high-speed relative movement between movable toothed disc and fixed toothed disc, materials are broken through comprehensive actions such as teeth shocking, shearing, friction and material bumping.

TECHNICAL DATA:

| Production capacity (kg/h): | 100~200 |
|------------------------------------|-------------------------------|
| Spindle speed(r/min): | 3,800 |
| Finish product grain fineness(mm): | 0.17~0.12 |
| Movable toothed disc diameter(mm): | 318 |
| POWER(kW): | 5. 5 |
| OVERALL DIMENSIONS(mm): | $720 \times 590 \times 1,400$ |

OTHER SPECIFICATIONS:

30B III

产品名称:万能粉碎机

型号规格:30B I

- **主要用途**:适用于医药、化工、食品、农业及粮食等行业, 粉碎干燥的脆性物料。
- 工作原理及结构特点:本机利用活动齿盘和固定齿盘间的高速相对运动,使被粉碎物经齿冲击、剪切、摩擦及物料彼此间撞击等综合作用,达到粉碎目的。

主要技术参数:

生产能力(kg/h):

100~200

主轴转速(r/min):

3800

成品粒度(mm):

 $0.17 \sim 0.12$

活动齿盘直径(mm):

318

配套电机(kW):

5.5

外形尺寸(mm):

 $720 \times 590 \times 1400$

其它规格:30B I

7.5

COMMODITY: HIGH EFFICIENT MILL

MODEL: GFSJ-16A

APPLICATION:

It is widely used for crushing material in pharmaceutical, foodstuff and chemical industry, etc.

FEATURES:

The material is crushed by impacting and shearing of high speed knives which have a sharp chin in the one side and an impacting fly-bar on the other side. The fineness of powder may be achieved by changing the different screen mesh.

TECHNICAL DATA:

POWER(kW):

| Production capacity(kg/h): | 50~300 |
|-----------------------------------|-------------------|
| Finished product grain fineness: | 12~120 |
| Adjustable rotating speed(r/min): | 1,000/2,400/4,400 |
| | |

OVERALL DIMENSIONS (mm): $1,100 \times 860 \times 1,520$

OTHER SPECIFICATIONS:

GFSJ-8A; GFSJ-18, GFSJ-32A

产品名称:高效粉碎机

型号规格:GFSJ-16A

主要用途:用于制药、食品、化工等行业粉碎物料。

工作原理及结构特点:采用一面为快口刀刃,另一面为冲击式刀片,使被粉碎物料经高速旋转刀片冲击、剪切而粉碎,粒度大小通过筛网获得。

主要技术参数:

生产能力(kg/h):

 $50 \sim 300$

成品粒度(目):

12~120

可调转速(r/min):

1000/2400/4000

配套电机(kW):

7.5

外形尺寸(mm):

 $1100 \times 860 \times 1520$

其它规格:GFSJ-8A,GFSJ-18,GFSJ-32A

COMMODITY: ULTRAFINE POWDER PULVERIZER

MODEL: CWF-600

APPLICATION:

This machine is suitable to pulverize the herbal medicine, dyes, food and ore, for instance, calcium carbonate, coconut shell, cassava, ginger, chili, mixed powder, rice, sand and so on in pharmaceutical, foodstuff and chemical industries, etc.

FEATURES:

A waveform abrasive wheel is used for pulverization. Pneumatic concentration is adopted for powder graining. Without screen and mesh, the grain size of powder may be adjusted according to the operator's instruction. The production process is carried out in a sealing system. It cancels any contamination from dust. In the pulverizer, circulating water cooling is adopted.

TECHNICAL DATA:

| Production capacity(kg/h): | 30~250 |
|---------------------------------------|--------|
| Finish product grain fineness (mesh): | 60~250 |
| Diameter of grinding chamber (mm): | Ø600 |

| POWER(kW): | 23. 8 |
|--------------------------|---------------------------------|
| OVERALL DIMENSIONS (mm): | $1,980 \times 865 \times 2,215$ |

OTHER SPECIFICATIONS:

CWF-360; CWF-900

产品名称:超微粒粉碎机组

型号规格:CWF--600

主要用途:用于制药、食品、化工等行业,粉碎中药、粮食、 矿砂等物料。

工作原理及结构特点:粉碎方式是采用浪纹式磨轮碾磨,物料磨碎后风送旋风离器分离,并经出料机构进入包装袋内。当粉粒吸送到旋风离器时,一部分小微粒沉积下来,落到集粉漏斗中,另一部分通过回风管返回粉碎机的磨室内。整个生产过程在封闭状态下完成,避免了污染。粉碎细度可任意调节。主要采用循环水冷却。

主要技术参数:

生产能力(kg/h):

 $30 \sim 250$

成品粒度(目):

 $60 \sim 250$

磨室直径(mm):

Ø600

总功率(kW):

23.8

配套电机(kW):

外形尺寸(mm):

 $1980 \times 865 \times 2215$

其它规格:CWF-360,CWF-900

60.55

COMMODITY: SELECTING PARTICLE MILL ASSEMBLY

MODEL: TF-700

APPLICATION:

It is used for grinding material in pharmaceutical industry, foodstuff industry, chemical industry, etc., and it is especially suitable for grinding material containing fiber.

FEATURES:

Materials are fed to mill by feeder, and evenly spread out around grinding room by fan, and sheared and impacted by revolving body at high speed and teeth circle, and ground under the condition of violent turbulence caused by revolving body. Then all of the ground products are transmitted to the winnower, qualified product is transmitted to collector by wind, meanwhile unqualified material is transmitted again to grinding room by wind.

TECHNICAL DATA:

POWER(kW):

| Production capacity(kg/h): | 80~1,000 |
|--|----------|
| Finish product grain fineness (Mesh): | 40~350 |
| Qualified ratio of product ($\%$, based on 100 mesh) | . 95∼100 |
| Spindle speed(r/min): | 3500 |
| | |

OVERALL DIMENSIONS(mm): $8,000 \times 5,000 \times 4,800$

OTHER SPECIFICATIONS:

F-400 mill; XFG series circulate boiling drier; D74-10T cutter.

产品名称:分粒式粉碎机组

型号规格:TF-700

主要用涂:用于中药、食品、化工、农药、各种饲料的粉碎 加工,尤其对含纤维的原料加工效果更佳。

工作原理及结构特点:物料通过自动加料器输送到粉碎 机中,被风板均匀地散布到粉碎室周围,在打板与牙 板间被剪切和冲击,在机内形成激烈涡流将物料粉 碎,并在气流作用下送到风选器内,经风板的作用, 将大颗粒择出再送回粉碎室,合格物料被吹入集料 装置,完成收粉。

主要技术参数:

生产能力(kg/h):

800~1000

成品粒度(目):

 $40 \sim 350$

一次出粉合格率(%): 95~100(以 100 目为准)

主轴转速(r/min):

3500

配套电机(kW):

60.55

外形尺寸(mm):

 $8000 \times 5000 \times 4800$

其它规格:F-400,XFG 系列,D74-107

COMMODITY: ULTRA-CRYOGENIC MILL

MODEL: DFJ-2

APPLICATION:

It is used for grinding materials in pharmaceutical, foodstuff and chemical industries, etc., especially, for grinding thermosensitive, oil-containing, sticky materials which are tenacious in ambient temperature.

FEATURES:

Materials to be ground are cooled below their brittle points by cryogen and ground to demanded fineness without destroying original ingredient. It consists of pre-cooling barrel, the spiral conveyer, low-temperature grinder, low-temperature air-blower, whirl separator, powder filter, electric control section etc.

TECHNICAL DATA:

| Production capacity(kg/h): | 10~180 |
|---------------------------------------|----------------|
| Finish product grain fineness (Mesh): | 80~500(option) |
| Feeding granularity(mm): | <5~10 |
| Temperature for grinding (°C): | 0~-196 |

POWER(kW): 6.75 **OVERALL DIMENSIONS(mm):** 2,160×1,860×2,300

OTHER SPECIFICATIONS:

WFS-15; WFS-30; WFJ-250

产品名称:超低温粉碎机

型号规格:DFJ-2

主要用途:用于制药、食品、化工等行业物料的粉碎,尤其适用于粉碎热敏性、油性、粘性、常温下软韧性物质。

工作原理及结构特点:用冷媒将被粉碎物冷却到脆化点以下,经粉碎机粉碎到所需的细度,原有成份不被破坏,本机由贮料预冷筒,螺旋推进器、低温粉碎机、低温鼓风机、旋风分离器、粉尘过滤器、电控柜等组成。

主要技术参数:

生产能力(kg/h):

10~180

成品粒度(目):

80~500(任选)

进料粒度(mm):

 $<5\sim10$

粉碎温度(℃):

 $0 \sim -196$

配套电机(kW):

6.75

外形尺寸(mm):

 $2160 \times 1860 \times 2300$

其它规格:WFS-15,WFS-30,WFS-250

COMMODITY: SPRAY AND RINSE MACHINE WITH HIGH-PRESSURE WATER

MODEL: BGX 500A

APPLICATION:

The machine is an advanced equipment of rinsing. It is suitable for washing the medicinal materials and foodstuffs etc.

FEATURES:

Through machine drive, the medicinal materials are continuously moved in by stepless speed regulating net belt. The upper four and the under four swing-rinsing water pipes sway back and forth to rinse every side with high-pressure water, the direction of which is crisscross to inlet direction of the medicinal materials.

TECHNICAL DATA:

| Production capacity(t/h): | 0.25~2.5 |
|--------------------------------------|----------|
| Working time(s): | 15~150 |
| Water pressure(MPa): | 0.2~0.3 |
| Water consumption(t/h): | 2~4 |
| Swing angle of spraying water pipes: | 0°~90° |

POWER(kW):

0.8

OVERALL DIMENSIONS(mm):

 $3,300 \times 850 \times 1,200$

产品名称:高压水喷射式洗净机

型号规格:BGX 500A

主要用途:用于洗净中药材、蔬菜、水果表面的泥土、污物。

工作原理及结构特点:用交错布置的 8 个高压水喷头,往 复摆动连续冲洗在网状传送带上的药材。无喷射死 角,洗净度高,破损率低。

主要技术参数:

生产能力(t/h):

 $0.25 \sim 2.5$

冲洗时间(s):

 $15 \sim 150$

水压(MPa):

 $0.2 \sim 0.3$

耗水量(t/h):

2~4

喷头摆动角度:

0°~90°

配套电机(kW):

0.8

外形尺寸(mm):

 $3300 \times 850 \times 1200$

COMMODITY: RECIPROCATING SLICING MACHINE

MODEL: WQY 240-2

APPLICATION:

It is suited to cut roots, stems, leaves, grass and barks of herbal medicines.

FEATURES:

The herbal medicines are transferred to the infeeder by the chains, then the reciprocating cutter cuts the herbal medicines to its required thickness. The thickness of the finished products is controlled by the stepless speed regulator on the infeeder. The machine can also be used for cutting similar material in other fields.

TECHNICAL DATA:

Production capacity (kg/h):

 $120 \sim 600$

Thickness of the slice(mm).

 $1\sim20$ (adjusted by stepless speed regulator)

Size of cutter gate(mm):

 240×50

POWER(kW):

3

OVERALL DIMENSIONS(mm): $1,500\times1,000\times1,200$

产品名称:往复式切药机

型号规格:WQY 240-2

主要用途:用于根、茎、叶、草、皮类中药饮片的切制。

工作原理及结构特点:药材用链条输送至进料口,由往复式 切刀切制成所需要厚度的饮片。无级调节进给机构可切制不同厚度饮片。

主要技术参数:

生产能力(kg/h):

120~600

切片厚度(mm):

1~20(无级调节)

刀片尺寸(mm):

 240×50

配套电机(kW):

3

外形尺寸(mm):

 $1500 \times 1000 \times 1200$

COMMODITY: ROTATING TYPE SLICING MACHINE

MODEL: QY 120-4

APPLICATION:

It is suited to cut roots, stems, leaves, grass and barks of herbal medicines.

FEATURES:

The herbal medicines are transferred to the infeeder by the chains, then the rotating cutter cuts the herbal medicines to its required thickness, The machine can also be used for cutting similar material in other fields.

TECHNICAL DATA:

| Production capacity(kg/h): | 70~700 |
|---------------------------------|---------------------|
| Thickness of the slice(mm): | 1, 1.5, 2, 3, 6, 12 |
| The number of cutters (pieces): | 3 |
| Size of cutter gate(mm): | 120×50 |
| | |
| POWER(kW): | 3 |

 $1,500 \times 750 \times 1,100$

OTHER SPECIFICATIONS:

QY 120-3, QY140-4, QY160-4

OVERALL DIMENSIONS(mm):

产品名称:旋转式切药机

型号规格:QY 120-4

主要用途:用于根、茎、叶、草、皮类中药饮片的切制。

工作原理及结构特点:上、下链条将物料输送至切药口, 旋转的切刀将物料切成规定厚度的饮片。

主要技术参数:

生产能力(kg/h):

70~700

切片厚度(mm):

1,1.5,2,3,6,12

切片数(片):

3

刀片尺寸(mm):

 120×50

配套电机(kW):

3

外形尺寸(mm):

 $1500 \times 750 \times 1100$

其它规格:QY120-3,QY140-4,QY160-4

COMMODITY: MULTI-EFFECT TUBULAR WATER DISTILLER

MODEL: LD 2000

APPLICATION:

It is used for making distilled water (WFI) or pure steam in pharmaceutical, chemical, foodstuff and electronic industries, etc.

FEATURES:

It is based on the principle of multi-effect evaporation of thin falling film on the pipe, where latent heat in steam can be utilized. Steam and liquid is separated totally by using S type baffle. Liquid distributor has special structure to ensure film forming even at low pressure. PLC controlled system is adopted. Water quality is conformed to the standard set by the pharmacopoeia of China for use in medical injection.

TECHNICAL DATA:

| Output of distilled water (L/h) : | 2,000 |
|---|-------------------------------|
| Pressure of steam(MPa): | 0. 3 |
| Steam consumption(kg/h): | 600 |
| Consumption of raw water(kg/h): | 2,150 |
| Consumption of cooling water (kg/h) : | 800 |
| POWER(kW): | 2. 57 |
| OVERALL DIMENSIONS(mm): | $3,600\times1,200\times4,000$ |

OTHER SPECIFICATIONS:

LD 500, LD 1000

产品名称:列管式多效蒸馏水机

型号规格:LD 2000

主要用途:用于医药、化工、食品、电子等行业制取蒸馏水 或纯蒸汽。

工作原理及结构特点:采用列管降膜与多效蒸发原理,蒸汽中的潜热多次被利用。汽、液分离用 S 型挡板,分离彻底。液体分布器采用特殊结构,保证在低压力下液体也能很好成膜。整机运行 PLC 全自动控制。水质符合《中国药典》1995 年版"注射用水"规定指标。

主要技术参数:

| 蒸馏水产量(kg/h): | 2000 |
|--------------|------|
| 蒸汽压力(MPa): | 0.3 |
| 蒸汽耗量(kg/h): | 600 |
| 原料水用量(kg/h): | 2150 |
| 冷却水用量(kg/h): | 800 |

配套电机(kW):

2.57

外形尺寸(mm):

 $3600 \times 1200 \times 4000$

其它规格:LD500,LD1000

COMMODITY: THERMOCOMPRESSION DISTILLER

MODEL: RY 1000

APPLICATION:

It is used for making distilled water in pharmaceutical (WFI), foodstuff and electronic industries.

FEATURES:

This distiller is based on the principle of heat pump evaporating, where condensation and evaporation use one heat exchanger without cooling water, so more than 90% energy can be saved. PLC controlled system is adopted. The distilled water is conformed to the standard set by pharmacopoeia of China for use in medical injection.

TECHNICAL DATA:

| Production capacity(L/h): | 1,000 |
|--------------------------------------|-------|
| Steam pressure (MPa): | 0.1 |
| Steam consumption(kg/h): | 20 |
| Temperature of distilled water (°C); | 30~90 |
| Consumption for raw water(L/h): | 1,020 |
| POWER(kW): | 22 |

 $3,550 \times 2,200 \times 3,010$

OTHER SPECIFICATIONS:

OVERALL DIMENSIONS(mm):

RY 500, RY 2000

产品名称:热压式蒸馏水机

型号规格:RY 1000

主要用途:用于医药、食品、电子等行业制取蒸馏水。

工作原理及结构特点:采用热泵蒸发原理,蒸发与冷凝共用一个换热器,不用冷却水,节能 90%以上。整机采用 PLC 实现全自动控制。蒸馏水质符合《中国药典》 1995 年版"注射用水"规定指标。

主要技术参数:

生产能力(kg/h): 1000

蒸汽压力(MPa): 0.1

蒸汽耗量(kg/h): 20

蒸馏水温度(℃): 30~90(可调)

原料水耗量(kg/h): 1020

配套电机(kW):

22

外形尺寸(mm):

 $3550 \times 2200 \times 3010$

其它规格:RY500,RY2000

COMMODITY: ELECTRODIALYSIS-MIXED BED PURE WATER MAKING EQUIPMENT

MODEL: 10T

APPLICATION:

It is suited to make pure water for pharmaceutical production (powder injection, liquid drug injection).

FEATURES:

High cost performance is due to the application of desalting electrodialysis and mixed bed. Recovery of concentrate water helps save water and money invested. Automatic suction and simultaneously regeneration decrease second pollution. All the electric apparatus are of central control. Automatic stop acts as the water reaches a certain level. It is possible to apply remote control.

TECHNICAL DATA:

| Pure water output(t/h): | 10 |
|--|----------|
| Filter precision (µm): | 1 |
| Pure water resistivity (M $\Omega \cdot cm$): | ≥2 (25℃) |
| Number of bacteria: | < 10 |
| | |
| POWER(kW): | 30 |
| OVERALL DIMENSIONS(m ²): | 120 |

产品名称:电渗析一混床纯水制备生产线

型号规格:10t

- 主要用途:用于医药工业生产过程用纯水的制备(粉针剂、水针剂等)。
- 工作原理及结构特点:由于采用脱盐电渗析和混床使之具有离心性能。浓水回收有利于节水和节约投资。自动吸收并同时再生减少了二次污染。全部电气设备集中控制。具有定水位自动停板功能,可以进行遥控。

主要技术参数:

纯水产量(t/h):

10(其余规格:0.5~30

纯水电阻率($M\Omega \cdot cm$):

 $\geq 2(25^{\circ}\text{C})$

过滤精度(μm):

1

细菌(个):

< 10

总功率(kW):

30

占地面积(m²):

120

配套电机(kW):

外形尺寸(mm):

COMMODITY:

ELECTROMAGNETIC INDUCTION SEALING CAP MACHINE

MODEL: DSPF-I

APPLICATION:

It is used on the bottle's mouth of plastic bottle or glass bottle to be sealed with aluminum foil in pharmaceutical, foodstuff and chemical industries, etc. It is effective to prevent from overflow, moisture, mildew, fake and damage to the goods in the bottle.

FEATURES:

The main functions are to cover and turn the bottle cap with preset aluminum foil, make it pass through electromagnetic induction head following conveyer belting without contact, and cause aluminum foil to be sealed on mouth. This machine has the function of self-detecting and rejecting bottles, whose cap have not a aluminum foil. The electromagnetic wave to be used will not leak out, the strength of electromagnetic radiation in near field is <10V/m.

TECHNICAL DATA:

Production capacity(bottle/min):

 $10\sim100$ (stepless adjustable)

Bottle sizes (mm): Max. \emptyset 70×200 (plastic or glass bottle) Bottle cap sizes (mm): Max. \emptyset 70, Min \emptyset 15

POWER(kW): $2.6\sim3$

OVERALL DIMENSIONS(mm): $1,500 \times 500 \times 1,200$

(conveyer and main machine)

OTHER SPECIFICATIONS:

DSPF-II is a portable type.

产品名称:电磁感应式铝箔封口机

型号规格:DSPF-1

- **主要用途**:用于制药、食品等行业将铝箔密封在塑料瓶或玻璃瓶瓶口上。可有效地防止溢出、变潮、发霉、假冒和损坏。
- 工作原理及结构特点:将预置铝箔的瓶盖旋合在瓶口上,使其随输送带通过电磁感应加热器,铝箔被无接触热封在瓶口上。本机具有瓶盖内无铝箔自检和剔废功能,所用电磁波不外泄,近场电磁辐射强度<10V/m。

主要技术参数:

生产能力(瓶/分):

 $10 \sim 100$

(根据瓶口规格无级可调)

瓶子规格(mm):

最大Ø50×200(塑、玻瓶)

瓶盖规格(mm):

最大Ø70、最小Ø15

配套电机(kW):

 $2.06 \sim 3$

外形尺寸(mm):

输送主机:1500×500×1200

COMMODITY: CAPSULE PRINTER, TABLET PRINTER

MODEL: AP-III

APPLICATION:

Ink printing on tablet or capsule surface.

FEATURES:

Multi-directional or lateral printing on capsule/tablet surface is based on the technique of automatic numeral information display and on the principle of an automatic mechanical divider.

TECHNICAL DATA:

Production capacity(tablets or capsules/h): $50,000 \sim 80,000$ Diameter of applicable tablet(mm): $5 \sim 8$ Model of applicable capsule: 00° , 0° , 1° , 2° , 3° , 4° , 5°

POWER(kW):

0.25

OVERALL DIMENSIONS(mm):

 $900 \times 610 \times 640$

OTHER SPECIFICATIONS:

Model DP, model RP

产品名称:胶囊、药片印字机

型号规格:AP-Ⅱ

主要用途:用于胶囊、糖衣片表面印字。

工作原理及结构特点:采用全自动数显技术,利用机械等分原理对胶囊进行轴向、不定向(帽体随机)转移方式印字。

主要技术参数:

生产能力(万粒/h):

 $5 \sim 8$

药片直径(mm):

 $5 \sim 8$

胶囊规格:

 $00^{+}, 0^{+}, 1^{+}, 2^{+}, 3^{+}, 4^{+}, 5^{+}$

配套电机(kW):

0.25

外形尺寸(mm):

 $900 \times 610 \times 640$

其它规格:DP 系列、RP 系列

COMMODITY: DRUM LABEL-STICKING MACHINE

MODEL: ZT 20/ 1000

APPLICATION:

It is suited to stick label on columnar containers of diameter 20-100mm.

FEATURES:

By mechanical drive, the bottles on the conveying truck are carried to the big-vacuum drum, the labels are drawn from label-case by small-vacuum drum, after printed and pasted, they are carried to big-vacuum drum. The labels absorbed on the big-vacuum drum are sticked on the bottles, and fastened through rolling and pressing.

TECHNICAL DATA:

| Production capacity(bottle/min): | 20~80 |
|----------------------------------|------------------|
| Volume of bottles(ml): | 20~1,000 |
| Size of the label(mm): | 128×120 |
| Vacuum(Pa): | 7,993 |

2, 05 POWER(kW): OVERALL DIMENSIONS(mm): $2,800 \times 775 \times 1,100$

产品名称:转鼓贴标机

型号规格:ZT 20/1000

主要用途:用于医药、食品、电子等行业圆柱形容器贴标签。

工作原理及结构特点:通过机械传动,瓶子由输瓶轨道经变距螺杆送至大真空转鼓,小真空转鼓从签盒中吸出标签,经打批号及涂浆送至大真空转鼓,吸附在大真空转鼓上的标签贴在瓶身上,经搓、滚粘牢。

主要技术参数:

生产能力(瓶/min):

20~80(无级调速)

容器规格(ml):

20~1000(圆柱形瓶)

标签规格(mm):

 128×110

工作真空度(Pa):

7993

配套电机(kW):

2, 05

外形尺寸(mm):

 $2800 \times 775 \times 1100$

COMMODITY: FULLY-AUTOMATIC LABELER FOR AN-TIBIOTIC BOTTLES

MODEL: TQ-3

APPLICATION:

It is suited to label and print production date and invalid date for antibiotic bottles automatically.

FEATURES:

- 1. Rectilinear bottles transfer and rectilinear label sucker are adopted, so machine works reliably and quickly.
- 2. Labeling and printing stop for absent bottle.
- 3. All of the processes (transferring bottles, transferring label, gluing on label, labeling on bottle, rolling depression, printing) are finished automatically in a definite sequence.

TECHNICAL DATA:

| Labeling rate(bottles/min): | 70~300 |
|-----------------------------|------------------------------------|
| Labeling tolerance(mm):: | ± 0.5 |
| Broken bottle rate(%): | <0.01 |
| Bottle size applied(mm): | \emptyset 22 \sim 25 \pm 0.5 |
| | |
| POWER(kW): | 2.5 |
| OVERALL DIMENSIONS(mm): | $3,750 \times 1,630 \times 1,140$ |

产品名称:抗生素瓶全自动标签粘贴机

型号规格:TQ-3

主要用途:用于抗生素玻璃瓶贴标签、打印批号。

工作原理及结构特点:输瓶为直线式运转平稳。缺瓶自动停止贴签和打印。按预定程序自动完成送瓶、送签、涂胶、印字、贴签、滚压牢实等全过程。

主要技术参数:

贴签速度(瓶/min):

 $70 \sim 300$

贴签误差(mm):

 ± 0.5

破瓶率(%):

< 0.01

瓶子直径(mm):

 $22 \sim 25 \pm 0.5$

功耗(kW):

2.5

配套电机(kW):

外形尺寸(mm):

 $4750 \times 1630 \times 1140$

COMMODITY: GRAIN AUTOMATIC PACKAGING MA-CHINE

MODEL: SB-80

APPLICATION:

It is suited to pack grain materials with small package in pharmaceutical and foodstuff industries, etc.

FEATURES:

The machine makes flat bags with three sides sealed, it can automatically finish bag making, measuring, filling, sealing, printing batch number, notch-tearing, cutting and counting processes.

TECHNICAL DATA:

| Production capacity(bag/min): | 55~90 |
|-------------------------------|-----------|
| Range of measurement(ml): | 10~50 |
| Bag size(mm): | L: 50~120 |
| | W: 60~85 |

Packing materials: paper/polyethylene, cellophane/polyethylene, Boop/ polyethylene, aluminum/ polyethylene, polyester/plated aluminum/polyethylene, foil, etc.

POWER(kW): 0.79 **OVERALL DIMENSIONS(mm):** $625 \times 751 \times 1,558$

OTHER SPECIFICATIONS:

SB-60P

产品名称:颗粒自动包装机

型号规格:SB-80

主要用途:用于医药、食品等行业颗粒物料的小袋包装。

工作原理及结构特点:该机是三边封合的扁平袋型包装机可自动完成制袋、计量、充填、封合、打印批号、加易撕口、切断及计数等全过程。

主要技术参数:

生产能力(袋/min):

55~90

计量范围(ml):

10~50

制袋尺寸(mm):

长 50~120、宽 60~85

包装材料:

纸/聚乙烯、玻璃纸/聚乙烯

BOOP/聚乙烯、铝箔/聚乙烯

涤纶/镀铝/聚乙烯等。

配套电机(kW):

0.79

外形尺寸(mm):

 $625 \times 751 \times 1558$

其它规格:SB-60P

183

COMMODITY: BLISTER PACKING MACHINE

MODEL: DPA 250

APPLICATION:

It is suited to package tablets and capsules.

FEATURES:

The PVC film heated is formed by vacuum pumped system in forming roller. After filling products, it is sealed with aluminum foil in hot-pressing roller. Then it is printed with batch number and cut into plates.

TECHNICAL DATA:

Production capacity (tablets/min):

 $800 \sim 1,104$

(capsules/min):

 $800 \sim 920$

Packing material: PVC: width: 250mm, thickness: 0.25mm PTP AL foil:

width: 250mm, thickness: 0.02mm

POWER(kW):

7.5

OVERALL DIMENSIONS(mm):

 $2,460\times1,030\times1,485$

OTHER SPECIFICATIONS:

DPB 210, DPT 130, DPT 220

产品名称:泡罩包装机

型号规格:DPA 250

主要用途:用于医药行业包装片剂、胶囊。

工作原理及结构特点:PVC 加热后,经模辊真空吸塑成 型,填入药品与铝箔辊热压封合,打印批号冲切成板 块。

主要技术参数:

生产能力(粒/min): 800~1104(糖衣片、素片)

800~920(胶囊)

包装材料尺寸(mm):

PVC 宽 250、厚 0.25

铝箔 宽 245~250、厚 0.02

配套电机(kW):

7.5

外形尺寸(mm):

 $2460 \times 1030 \times 1485$

其它规格:DPB210,DPT130,DPT220

COMMODITY:

AUTOMATIC

ALUMINUM-PLASTICS

PACKER

MODEL: DPT-140

APPLICATION:

Used for packing tablets, capsules, suppositories and medical utensils in pharmaceutical industry.

FEATURES:

PVC after having been heated will form bubble cap under negative pressure, and after the medicine material is fed in, the bubble cap will be heat sealed through knurled roll and aluminum foil, and then it will be punched into finished plates.

TECHNICAL DATA:

Production capacity (pcs/h):

 $50,000 \sim 120,000$

Packaging materials:

Medicinal Aluminum foil (PTP)

Polyvinyl Chloride (PVC)

Dialysis Paper

POWER(kW):

0.55

OVERALL DIMENSIONS(mm):

 $1,900 \times 1,000 \times 1,450$

产品名称:自动铝塑包装机

型号规格:DPT-140

- **主要用途:**用于医药行业包装片剂、胶囊、栓剂、医疗器 具。
- 工作原理及结构特点:经加热的 PVC 在负压下泡罩成型,加料后经网纹辊与铝箔热封,冲切成板块成品。

主要技术参数:

生产能力(万粒/h):

 $5 \sim 12$

包装材料:

药用铝箔(PTP)

聚氯乙烯(PVC)

透折纸

配套电机(kW):

0.55

外形尺寸(mm):

 $1900 \times 1000 \times 1450$

COMMODITY: PLATE AL-PLASTIC & AL-AL PACK-AGING MACHINE

MODEL: DPP 170-S

APPLICATION:

Used for packaging various kinds of tablets, capsules, gummy pills, suppositories and special tablets in pharmaceutical industry and foodstuff industry.

FEATURES:

The machine is a multi-functional blister packaging machine. It uses upright-pressed forming and heat sealing. The cold press forming for aluminum blister, the photoelectricity detection, optical fiber, and light guide technique are adopted. It can distinguish the tablet missed plate, automatically. The work process is controlled by a Programmable Logical Controller. The change of Al-Al packing to Al-Plastic packing is carried out only by changing forming tools.

TECHNICAL DATA:

| Max. capacity (times/min): | 40 |
|----------------------------|---------------------------------|
| Max. plate size (mm): | 154×104 |
| Max. forming depth(mm): | 10 |
| Power for heating(kw): | 6 |
| POWER(kW): | 1.5 |
| OVERALL DIMENSIONS(mm); | $2,445 \times 980 \times 1,450$ |

OTHER SPECIFICATIONS:

DPP 170, DPT 250

产品名称:平板铝塑、铝一铝包装机

型号规格:DPP170-S

主要用途:用于医药行业包装片剂、胶囊、胶丸、栓剂及异形片。

工作原理及结构特点:本机是一台多功能泡罩包装机,增加了铝泡罩冷压成型和光电检测功能。采用光导纤维和光导技术可自动识别缺药板型。整机由一台可编程控制器控制,铝一铝和铝塑包装的更换,通过改变成型模具实现。

主要技术参数:

冲切次数(次/min):

40

最大板型尺寸(mm):

 154×104

成型深度(mm):

10

加热功率(kW):

6

配套电机(kW):

1.5

外形尺寸(mm):

 $2445 \times 980 \times 1450$

其它规格:DPP170 DPT250

COMMODITY: MULTI-FUNCTION AND PLATE TYPE AUTO AL-PLASTIC BUBBLE CAP PACK-ING MACHINE

MODEL: DPP-200K

APPLICATION:

For packing tablets and capsules etc. in pharmaceutical industry.

FEATURES:

It adopts the stepmotor, plate-type mould, pressure formation and frequency control. With this machine, two-color chromatograghy printing on front and back of aluminum foils at same time, direction and color separation and arrangement of capsule can be made. It can fill various medicine (for taking in one course of treatment) in the same sheet at same time. It is also suitable for packing of aluminum/aluminum.

TECHNICAL DATA:

| Punch frequency (times/min): | 30~80 |
|------------------------------|---------|
| Max. shaping area (mm): | 190×130 |
| Max. shaping depth(mm): | 18 |
| Travel(mm): | 40~140 |

POWER(kW):

OVERALL DIMENSIONS(mm): $2,070 \times 960 \times 1,750$

OTHER SPECIFICATIONS:

DPP-250B

产品名称:多功能辊板式自动铝塑包装机

型号规格:DPP-200K

主要用途:用于医药、食品等行业包装片剂、胶囊等。

工作原理及结构特点:采用板式模具、正压成型、无边冲裁、电子计数,配"封闭式"凹版印刷及吸尘装置。无级变速,并具有双色套印及压易折线功能。本机可同时完成铝箔正反面的双色套印、胶囊、粗细和颜色的区分排列,胶囊可充入各种药物,还可用于铝一铝包装。

主要技术参数:

冲裁频率(次/min):

 $30 \sim 80$

最大成型面积(mm):

 190×130

最大成型深度(mm):

18

行程(mm):

 $40 \sim 140$

配套电机(kW):

6. 1

外形尺寸(mm):

 $2370 \times 950 \times 1600$

其它规格:DPP-250B

COMMODITY: AUTOMATIC HIGH-SPEED DOUBLE A-LUMINUM FOIL PACKER

MODEL: ZGL-160

APPLICATION:

Applied to double aluminum packing of tablet and capsule, etc, in pharmaceutical industry.

FEATURES:

The packer features adopting frequency control of motor speed, punching the array feeding material into required plates after it has been heat sealed by the moulding wheel, and automatic temperature regulation, automatic counting of finished product and other functions.

TECHNICAL DATA:

| Production capacity(pcs/h): | 22,000~38,000 |
|-----------------------------|---------------------------------|
| Max. packing width(mm): | 160 |
| Punching size(rows): | 1~5 |
| Air pressure(MPa): | 0.4 |
| | |
| POWER(kW); | 0.4 |
| OVERALL DIMENSIONS(mm): | $1,200 \times 850 \times 1,650$ |

产品名称:自动高速双铝箔包装机

型号规格: ZGL-16

主要用途:用于医药行业片剂、胶囊等的双铝箔包装。

工作原理及结构特点:采用变频调速,整列送料经模轮热 封后冲成所需板块,具有温度自动调节、成品自动计 数等功能。

主要技术参数:

生产能力(万粒/min):2.2~3.8最大包装宽度(mm):160冲切规格(排):1~5空气压力(MPa):0.4

配套电机(kW):

0.4

外形尺寸(mm):

 $1220 \times 850 \times 1650$

COMMODITY: PLASTIC BOTTLE FILLING, COUNTING and SEALING ASSEMBLY LINE

MODEL: SBX-2

APPLICATION:

It is widely used in the pharmaceutical industry to bottle-pack the tablets, pills, capsules and the similar forms of medicines.

FEATURES:

The whole line consists of strengthening machine; computerized counting-filling machine; sealing machine; automatic caps screwing machine, and self-adhesive labeling machine. The compact line can complete all the procedure automatically for bottling the tablets, pills and capsules.

TECHNICAL DATA:

Production capacity (bottle/min): 60
Diameter of the bottle(mm): 25~70

POWER(kW):

3

OVERALL DIMENSIONS(mm):

 $7,000 \times 2,000 \times 1,800$

产品名称:塑料瓶包装生产线

型号规格:SBX-2

主要用途:用于医药行业片剂、丸剂、胶囊类药物的瓶装。

工作原理及结构特点:全线采用塑料瓶理瓶机、条板式电脑计数自动装瓶机、塑料瓶封口机、自动旋盖机、不干胶贴标机等单机组成。可自动完成理瓶、计数、装瓶、复合膜封口、旋盖、贴签等工序。

主要技术参数:

生产能力(瓶/min):

60

瓶子直径(mm):

 $25 \sim 70$

配套电机(kW):

3

外形尺寸(mm):

 $7000 \times 2000 \times 1800$

COMMODITY: PLASTIC-BOTTLE ASSEMBLY LINE

MODEL: BXSP 15/100

APPLICATION:

For bottling tablets, pills and capsules in the pharmaceutical industry.

FEATURES:

The line can complete procedures such as arranging bottles, conveying bottles, counting or quantitative filling, stuffing paper, arranging corks, stuffing inside-corks, sealing-wax, arranging caps, screwing caps, labeling, printing serial numbers, etc.

TECHNICAL DATA:

Production capacity(bottles/min): 40~80
Size(ml): 15~100(big mouth medicinal bottles)
Error of counting & metering for tablets(%): ≤0.2
Metering error for small pills: ≤3(weight meter)

POWER(kW): 7. 3 **OVERALL DIMENSIONS(mm):** 15,000×1,500×2,100

产品名称:塑料瓶瓶装自动线

型号规格:BXSP 15/100

主要用途:用于医药行业片剂、丸剂、胶囊类药物的瓶装。

工作原理及结构特点:由 UZPS1000 型理瓶机、LJ20/1000 型履带计数机(或 WP8 型小丸瓶装机)、SZ2 型塞纸机、FYS8 型压塞机、FNL 型绞龙式封蜡机、FXG8 型旋盖机、JNJ30/80 型不干胶贴签机(或ZT20/1000 型转鼓贴签机)组成生产线。全线采用单机组合,独立机械传动,无级调速,各单机装有电控装置以平衡全线速度。自动完成理瓶、输瓶、计数装瓶、塞纸、塞内塞、封蜡、理盖、旋盖、印批号、贴标签等工序。

主要技术参数:

生产能力(瓶/min):

40~80

瓶子规格(ml):

15~100

(大口药瓶)

计数计量误差(%):

≤0.2(片)、≤3(小丸)

配套电机(kW):

7.3

外形尺寸(mm):

 $15000 \times 1500 \times 2100$

COMMODITY: AUTOMATIC STRAPPER

MODEL: SK-1A

APPLICATION:

It is used for strapping the corrugated paper, wooden and plastic cases.

FEATURES:

The plastic strap is used as strapping material. The complete strapping procedure can be automatically finished except for carrying the parcels. The junction is tight due to heat sealing. It can auto-stop as something wrong happens to convey belt.

TECHNICAL DATA:

| Max. strapping size(mm): | 600×600 |
|---------------------------------|---------------------------------|
| Strapping force(N): | 700 |
| Strapping time per strap(s): | 3.5 |
| Height of working platform(mm): | 845 |
| Specification of PP strap(mm): | 13. 5 |
| POWER(kW): | 0.55 |
| OVERALL DIMENSIONS(mm): | $1,260 \times 636 \times 1,530$ |

OTHER SPECIFICATIONS:

SK-1D, SK88

产品名称:自动捆扎机

型号规格:SK-1A

主要用途:用于瓦楞纸箱、钙塑箱、木箱等类包装件的自动捆扎。

工作原理及结构特点:用塑料打包带做为捆扎材料,除送包外,所有捆扎动作均自动完成。塑料带接头热合牢固,如送带出现故障,能自动停机。

主要技术参数:

| 最大捆扎尺寸(mm): | 600×600 |
|-----------------|------------------|
| 捆扎力(N): | 700 |
| 每道捆扎时间(s): | 3.5 |
| 工作台高度(mm): | 845 |
| PP 塑料打包带规格(mm): | 13.5 |
| | |

配套电机(kW):

0.55

外形尺寸(mm):

 $1260 \times 636 \times 1530$

其它规格:SK-1D,SK88

COMMODITY: GLASS TUBE DRAWING PRODUCTION LINE

MODEL: BLGX10

APPLICATION:

This production line can be used to produce the neutral glass tubes for pharmaceutical industry.

FEATURES:

This production line consists of feeding channel, muffle furnace, danner machine, air-blow device, forming chamber, gas mixing station, cooling bridge, laser tester, tube drawing machine, pre-cutter, separator, cooling transfer, fine-cutting and smoothing machine, breaker and electrical control cabinet, etc.

TECHNICAL DATA:

| Production capacity(t/d): | 10 |
|-------------------------------|-----------------------------|
| Diameter of tube(mm): | 10~37.5 |
| Wall thickness of tube(mm): | 0.5~1 |
| Length(mm): | $700\sim1,400$ for cutting |
| | 1,400~2,000 for non-cutting |
| Drawing speed(m/min): | 7~280 |
| Length of production line(m): | 54~66 |
| POWER(kW): | |
| OVERALL DIMENSIONS(mm) |): |

产品名称:玻璃拉管生产线

型号规格:BLGX 10

主要用途:用于拉制医药用中性玻璃管。

工作原理及结构特点:该生产线由供料道、马弗炉、丹纳机、吹气装置、成型室、燃气混合站、冷却跑道、激光检测装置、牵引机、初切装置、分选机、冷却传送装置、精切圆口机、破碎机、电控箱等部分组成。

主要技术参数:

平均产量(t/d):

10

拉管直径(mm):

 $10 \sim 37.5$

拉管壁厚(mm):

0.5 \sim 1

以及(IIIII);

长度(mm): 中切 700~1400、不中切 1400~2000

牵引速度(m/min):

 $7 \sim 280$

生产线长度(m):

 $54 \sim 66$

配套电机(kW):

外形尺寸(mm):

COMMODITY: VERTICAL AMPOULE PRODUCTION LINE

MODEL: ZA 32

APPLICATION:

This production line can be used to produce neck breakable ampoules of $1\sim20$ ml.

FEATURES:

Glass tube loading, preheating, drawing, necking, bottom cutting, bottom sealing, flame cutting, scoring, applying color ring, annealing, detecting and output of passed ampoules, etc can be accomplished automatically.

TECHNICAL DATA:

Production capacity (ampoules/h):

 $4500 \sim 5700$ (for ampoules of 1, 2 and 5ml)

3800~4400 (for ampoules of 10 and 20 ml)

Percent of pass(%):

> 85

Heating value of gas(MJ):

> 16.74

POWER(kW):

5.98

OVERALL DIMENSIONS(mm):

1800×1800×1700 (for principal machine) The whole production line is 12,996mm in length.

产品名称:立式安瓿生产线

型号规格:ZA 32

主要用途:用于制造 1~20ml 曲颈易折安瓿。

工作原理及结构特点:该机可自动完成上管、预热、拉丝、制颈、切底、封底、切丝、划痕、涂色环、退火、检测、出合格安瓿等过程。

主要技术参数:

生产能力(支/h):

 $4500 \sim 5700(1,2,5ml)$

 $3800 \sim 4400 (10, 20 \text{ml})$

合格率(%):

> 85

煤气热值(MJ):

>16.74

配套电机(kW):

5.98

外形尺寸(mm):

主机 1800×1800×1700

生产线总长 12996

COMMODITY: GLASS TUBE VIAL MAKING LINE

MODEL: GZPX 40

APPLICATION:

It is used to make glass tube vials for the pharmaceutical industry.

FEATURES:

This line is composed of mechanical system, hydraulic system, compressed air supply system and gas system. Charging, heating, forming, cutting, and outputting of finished products can be finished automatically.

TECHNICAL DATA:

| Production capacity(vials/h): | 1,800~2,400 |
|-------------------------------------|--------------|
| Diameter of vial(mm): | 14~34 |
| Height of vial(mm): | ≤ 110 |
| Percent of pass(%): | 80 |
| Gas consumption(m ³ /h): | 3~5 |
| Heating value of gas(MJ/m³); | 12.5~16.7 |
| POWER(kW): | 43. 5 |

 $11,700 \times 2,400 \times 4,260$

OTHER SPECIFICATIONS:

OVERALL DIMENSIONS(mm):

产品名称:管制瓶生产线

型号规格:GZPX 40

主要用途:用于医药行业制造管制瓶。

工作原理及结构特点:生产线由机械系统、液压系统、空气动力系统、燃气系统组成,可自动完成上料、加热、成型、切片、输出成品全过程。

主要技术参数:

| 生产能力(只/h): | 1800~2400 |
|--------------|---------------------------------|
| 瓶子直径(mm): | 14~34 |
| 瓶子高度(mm): | €110 |
| 合格率(%): | 80 |
| 煤气消耗(m³/h): | 3~5 |
| 煤气热值(MJ/m³): | 12.5~16.7 |
| 配套电机(kW): | 43.5 |
| 外形尺寸(mm): | $11700 \times 2400 \times 4260$ |

COMMODITY: PRODUCTION LINE FOR HARD CAP-SULE

MODEL: JSZ-III

APPLICATION:

It is used to produce hollow hard capsule in pharmaceutical, chemical industries, etc.

FEATURES:

During the operation, the production line can automatically complete the process of oiling, glue-dipping, drying, demoulding, cutting-out, sleeving and product-conveying sequentially. The safety inter-locks, over-loading protection device and photo-sound warning system are equipped at 20 sections of the production line.

TECHNICAL DATA:

Capsule specifications:

0", 1", 2", 3", 4"

Production capacity(capsules/h):

42,000~52,000

Number of die set:

900 sets for the body (30 cones per set)

900 sets for the hood (30 cones per set)

POWER(kW):

4

OVERALL DIMENSIONS(mm):

 $12,000 \times 1,750 \times 2,500$

产品名称:硬胶囊生产自动线

型号规格:JSZ-I

主要用途:用于医药、化工等行业生产空心硬胶囊。

工作原理及结构特点:在机器运转过程中,依次自动完成进油、蘸胶、烘干脱模、切断、套合、输出成品。设有 20 个安全互锁、过载保护及声光报警。

主要技术参数:

胶囊规格:

 $0^{+}, 1^{+}, 2^{+}, 3^{+}, 4^{+}$

生产能力(万粒/h):

 $4.2 \sim 5.2$

模具数量:

体 900 条(每条 30 锥)

帽 900 条(每条 30 锥)

配套电机(kW):

4

外形尺寸(mm):

 $12000 \times 1750 \times 2500$

COMMODITY: TABLET HARDNESS TESTER

MODEL: YD-II

APPLICATION:

It is used to test the tablet hardness value for pharmaceutical detecting departments.

FEATURES:

The hardness tester is controlled by monolithic processor. The tablet is placed between the end of crown bar and pressure transducer. The crown bar feeds at constant speed to produce extrusion pressure. After changing-over the electrical signal, the greatest pressure value can be gotten and be printed out immediately.

TECHNICAL DATA:

| Full-scale range(N): | 9.806~196 |
|-----------------------|---------------------------|
| Testing accuracy(%): | ±1 |
| Tablets diameter(mm): | $1\sim35$ (can be preset) |
| Test time(s/tablet): | ≤3 0 |

POWER(kW): 0.065 **OVERALL DIMENSIONS(mm)**: 555×294×230

产品名称:片剂硬度计

型号规格:YD-1

主要用途:用于医药检测部门测定药片的硬度值。

工作原理及结构特点:该硬度计由单片机控制。将药片置 ,于顶杆头与压力换能器之间,顶杆恒速进给产生挤 压力,此电信号经转换后,即得到最大压力值、压力 值等,并由打印机打出测试结果。

主要技术参数:

量程(N):

9.806 \sim 196

测量精度(%):

+1

可测药片直径(mm):

1~35(可预置)

测量时间(s/片):

≪30

配套电机(kW):

0.065

外形尺寸(mm):

 $555 \times 294 \times 230$

COMMODITY: INTELLECTUAL PHARMACY DISSOLUTION TESTER

MODEL: RC-610 (100)

APPLICATION:

It is used to test medicine dissolution for pharmaceutical detecting departments.

FEATURES:

The rotation speed, temperature and timing of the tester are controlled by monolithic processor. The dissolution testing parameter on kinds of medicine or sampling tablet that have met China Pharmacopoeia (95 edition) can be written into EPROM chip and be called, displayed and run randomly. The simple sampling set is also provided with tester. The tester meets the standard of China Pharmacopoeia (95) and ZBC 95001-89 (Dissolution Tester).

TECHNICAL DATA:

Speed range(r/min): $20\sim200\pm2$ (1 rot/step, when preset) Temperature range($^{\circ}$):

 $20\sim40\pm0.1$ (0. 10170/step, when preset)

Timing (min): 0~999 (up to mine timing stages continuous)

POWER(kW):

1.3

OVERALL DIMENSIONS(mm);

 $1,100 \times 320 \times 450$

OTHER SPECIFICATIONS:

RC-610 (300, 600, 700 & 800)

产品名称:智能型药物溶出度测定仪

型号规格:RC-610(100)

主要用途:用于医药检验部门测定药物溶出度。

工作原理及结构特点:本测定仪由单片机控制转速、温度和定时。95 版药典规定的各种药品及标准片溶出试验参数写入 EPROM 芯片,可随机调取、显示、执行。随机附简易取样装置。本仪器符合《中国药典》1995年版及 ZBC95001—89《溶出试验仪》的有关规定。

主要技术参数:

调速范围(r/min): 20~200±2(预置时每步 1 转) 控温范围(℃): 20~40±0.1(预置时每步 0.1℃) 定时(min): 0~999(最多九个定时段连续执行) 总功率(kW): 1.3

配套电机(kW):

1.3

外形尺寸(mm):

 $1100 \times 320 \times 450$

其它规格:RC-610(300、600、700、800)

COMMODITY: INTELLECTUAL DISINTEGRATION TESTER

MODEL: BJ I

APPLICATION:

It is used to test disintegration rate within medicine disintegration time-limit.

FEATURES:

A shaft driven by two alternating motors controlled by monolithic processor drives two groups of basket to move up and down. The tester can determine the different tablet disintegration rate within the prescribed time limit, and can govern the water bath constant temperature or timing alarm as well.

TECHNICAL DATA:

Water temperature(℃):

 37 ± 0.3

Basket moving (up and down) speed(times/min):

30~32

POWER(kW):

0.06

OVERALL DIMENSIONS(mm):

 $365 \times 231 \times 382$

OTHER SPECIFICATIONS:

BJ(II), BJ(III)

产品名称:智能崩解仪

型号规格:BJ I

主要用途:用于测试药物在崩解时限中的崩解度。

工作原理及结构特点:通过两台交流电机驱动导杆带动两组吊篮上、下运动。该电机由单片机控制,可测定不同药片在规定时限中的崩解度,并可控制恒温水浴及定时报警。

主要技术参数:

水浴槽温度(C):

 37 ± 0.3

吊篮升降次数(次/min):

 $30 \sim 32$

配套电机(kW):

外形尺寸(mm):

 $365 \times 231 \times 382$

其它规格:BJI、BJI

COMMODITY: INTELLIGENT DISINTEGRATION TESTER

MODEL: ZBS-6D

APPLICATION:

This is a test instrument of medicines, which is specially used in testing the disintegration time limit of solid preparations such as tablets, capsules and pills.

FEATURES:

The instrument is an electron-mechanical product. It consists of a water bath, test baskets, beakers, a heating device and a microcomputer system for measurement and control. It is provided with automation functions of constant temperature, updown motion and timing controls. It accords with the regulations of the pharmacopoeia (USP XXIII), and it is easy to operate.

TECHNICAL DATA:

Raising-lowering frequency rate(runs/min): $30\sim32$ Moving distance(mm): 55 ± 2 When the basket is at the lowest position, the distance be-

When the basket is at the lowest position, the distance between the sieve and the bottom of the beaker is 25 ± 2 mm.

POWER(kW): 0.850 (Max) (50/60Hz, 220/110V) OVERALL DIMENSIONS(mm): 500×430×540

产品名称:智能崩解试验仪

型号规格:ZBS-6D

- 主要用途:用于医药行业检测片剂、胶囊、滴丸等固体制剂的崩解时限。
- 工作原理及结构特点:该仪器是一种机电一体化产品,由水浴、吊篮、烧杯、加热组件及微电脑测控装置构成,具有恒温、升降及定量等自动化功能,符合药典(US-PXX Ⅱ)的规定,并方便于操作。

主要技术参数:

往返频率(次/分):

 $30 \sim 32$

上下移动距离(mm):

 55 ± 2

吊篮降到最低位置时筛网距杯底 25±2

配套电机(kW):

50/60Hz,220/110V,0.85kWmax

外形尺寸(mm):

 $500 \times 430 \times 540$

COMMODITY: GEL STRENGTH TESTER

MODEL: JS-I

APPLICATION:

It is used to test gelatin strength on food additive and gelatin used for making capsule.

FEATURES:

Determining gel strength value by means of BLOOM is JS-I tester's measuring method. This is that the gelatin surface is pressed down by plunger (dia: 12.7mm) to produce the force (unit: g) which can make it cave-in 4mm. The gelatin concentration is 6.67% (gravity) and needs to be placed in jar, at 10 (C for 16-18 hours. Digital circuit and step-motor driving are adopted to ensure accuracy of 4mm.

TECHNICAL DATA:

| Gel strength value range(g): | 5~500 |
|-------------------------------|------------------|
| Gel strength accuracy(%): | $\pm 0.5 \pm 1g$ |
| Concave(4mm) error (mm); | ±0.1 |
| Water-basin temperature (°C): | 10 ± 0.1 |

POWER(kW):

1.2

OVERALL DIMENSIONS(mm):

265×230×530 (Strength testing table)

产品名称: 冻力测试仪

型号规格:JS-I

主要用途:供胶囊用明胶或食品添加剂作冻力测试用。

工作原理及结构特点:国际上均以 BLOOM 确定的冻力测试方法测定冻力,并定义冻力的量值 BLOOMg。其方法是用直径为 12.7mm 的柱塞压入冻胶表面,使其产生 4mm 深凹陷所需要的重量(以克表示)。冻胶浓度为 6.67%(重量计),并须在广口瓶中 10℃下置16~18 小时。本仪器结构上采用数字电路和步进电机驱动,以保证 4mm 的精确度。

主要技术参数:

冻力值测试范围(g):

5~500

冻力值精度(%):

 $\pm 0.5 \pm 1g$

凹陷 4mm 时误差(mm):

 ± 0.1

水浴箱水温(℃):

 10 ± 0.1

总功率(kW):

1. 2

配套电机(kW):

55BF003 步进电机

外形尺寸(mm):

265×230×530(冻力测试台)

COMMODITY: INTELLIGENT DRUG DISSOLUTION TESTER

MODEL: D-800

APPLICATION:

It is used in pharmaceutical, university and government laboratories for the dissolution determining of dosage form as tablet or capsule.

FEATURES:

The D-800 dissolution tester features 8 active spindles and an integrated external heater in a space-saving layout. The D-800 is designed to operate continuously according to USP XXIII specifications. The power lift head electronically seeks the correct height for paddles or baskets. The microprocessor controls the bath temperature within (0.1 (C and the spindle speed within (1% during the test.

TECHNICAL DATA:

| Spindle/vessel(pcs): | 8 |
|-------------------------|-------------------|
| Speed range (rpm): | $25\sim250\pm1\%$ |
| Temperature range (°C): | $25\sim50\pm0.1$ |
| Spindle diameter (mm): | 9. 525 |

| POWER(kW); | 0.8 (50/60 Hz, 110/220V) |
|-------------------------|-----------------------------|
| OVERALL DIMENSIONS(mm); | $630 \times 560 \times 880$ |

OTHER SPECIFICATIONS:

ZRS-6, ZRS-8, ZRS-8A

产品名称:智能药物溶出仪

型号规格:D-800

- **主要用途:**用于医药、院校和政府有关检验、试验部门确定片剂或胶囊剂的溶出度。
- 工作原理及结构特点: D—800 溶出度测定仪有 8 根转杆和一个紧凑的外置加热器, D—800 符合 USPXX II 的技术要求,适于连续工作。机头可电动寻找浆法或篮法的正确高度。微型计算机在测试过程中控制水箱温度精度±0.1℃,转杆转速精度±1%。

主要技术参数:

转杆/杯数:

8

转速范围:

 $20 \sim 250 \text{ RPM}: \pm 1\%$

温度范围(℃):

20~50°C:±0.1°C

转杆尺寸(mm):

9.525

配套电机(kW):

110/220V,60/50Hz,0.8kW

外形尺寸(mm):

 $630 \times 560 \times 880$

其它规格:ZRS-6,ZRS-8,RCZ-8A

COMMODITY: INTELLIGENT PARTICLE DETECTOR

MODEL: WI-2

APPLICATION:

It is used for detecting of the size and the number of undissolved particles in the injection.

FEATURES:

The particle detector, according to Coulter Technique, adopts Siphon control system composed of a vacuum pump, a negative pressure sensor and electromagnetic valves instead of mercury. It can take data of three channels at the same time. It is fashionable and easy to operate.

TECHNICAL DATA:

Detecting range((m):

 $5 \sim 40$

Set channels (μm) :

 $\geq 5, \geq 10, \geq 25$

Sample amount:

1ml

Counting time:

<30sec

Error in repeating count: $\leq 2\%$ (Particle count $\leq 5000/\text{ml}$)

POWER(kW):

0.07

OVERALL DIMENSIONS(mm):

 $355 \times 412 \times 308$

产品名称:智能微粒检测仪

型号规格:WJ-2

主要用途:用于测定注射液中不溶性微粒的大小和数量。

工作原理及结构特点:以库尔特计数技术为基础,采用了 真空泵、负压传感器、电磁阀组成的主动虹吸控制系 统,可同时测定三种不同粒径范围粒子的数量,操作 简单,方便实用。

主要技术参数:

测定范围(µm):

 $5 \sim 40$

设定通道(μm):

 $\geq 5, \geq 10, \geq 25$

采液量(ml):

1

计数时间(s):

€30

重复精度:

≤2%(粒子浓度≤5000 粒/ml)

配套电机(kW):

0.07

外形尺寸(mm):

 $355 \times 412 \times 308$

COMMODITY: BACTERIAL ENDOTOXINS TESTER

MODEL: BET-16

APPLICATION:

An automated test and analysis system is suitable for measuring the Bacterial Endotoxins by using Turbidimetric Kindtic Assay.

FEATURES:

The test and analysis system consists of responder and computer with software for controlling and analyzing. It measures the Bacterial Endotoxins by using Turbidimetric Kinetic Assay. It has an automated temperature controlling responder which detects and records turbidity changes for data processing later. Available range of Bacterial Endotoxins concentration is 0.001~100EU/ml.

TECHNICAL DATA:

| Temperature (°C): | 37.0 ± 0.3 |
|---------------------------|-------------------------|
| Test tubes: | 16 |
| Sample interval(s): | 10 |
| Total detective time (s): | 8,400 |
| | 3 |
| POWER(kW): | 0. 1 |
| OVERALL DIMENSIONS(mm): | $495\times290\times150$ |

产品名称:细胞内毒素测定仪

型号规格:Bet-16

主要用途:应用动态浊度法定量测定和分析细菌内毒素。

工作原理及结构特点:该系统由反应器、计算机以及运行于计算机上的控制分析软件组成利用动态浊度法的原理,能够自动恒温并检测浊度变化,能够测定 0. 001—100EU/ml 浓度范围的细菌内毒素。

主要技术参数:

温度(℃):

 37.0 ± 0.3

试管数(个):

16

采样间隔(秒):

10

最长定时(秒):

8400

配套电机(kW):

0.1

外形尺寸(mm):

 $495 \times 290 \times 150$


LIAOYANG PHARMACEUTICAL MACHINERY CO., LTD.

LEGAL REPRESENTATIVE: WU LIJIU

Liaoyang Pharmaceutical Machinery Co., Ltd. was set up in 1947 and it has the history of production of over 50 years. It is the largest pharmaceutical machinery manufacture in China.

Its main products are centrifugal separation machinery, glass—lined equipment, medicinal packing machinery, cycloidal reducer and design and manufacture of non-standard pressure vessels. Its products have reached the advanced level in China and some to the international level. The products are sold well more than 20 provinces, cities and autonomous regions. They are also exported into Pakistan, India, Thailand, Indonesia, Korea, South Korea, Spain and Argentina etc.



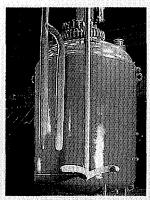




Address: No.2, Sec.1, Shengli Rd, Liaoyang City, Liaoning Province

Postal Code: 111000 Telephone: 2126405


Fax: 0419-2158106

Model QJT100/500A Cylinder-Tape Infusion-Fluid Bottle Cleaner

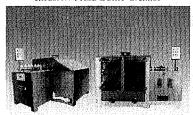


Glass Lined Reactor

Glass Lined Equipment

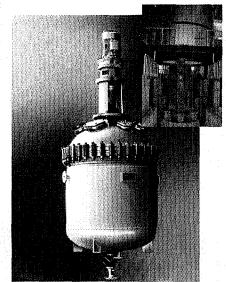

Model GC100/500B Glass Infusion Bottle Line Inject Type Filling Machine

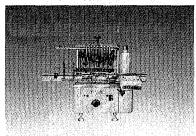


CHONGQING PHARMACEUTICAL
MACHINERY FACTORY

Address:No.35,Shapingba Main Street, Chongqing, P.R,China Tel:0086-23-65312352 Cable:1434 FAX:0086-23-65315646 Post Code:400030


Model FZJ100/500B Glass Infusion Bottle Rubber-plug Turning Machine




Model FGL100/1000B Glass Infusion Bottle Cap Pressing Machine

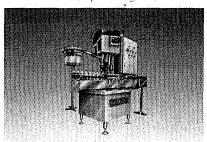

Model SR2540(3000° 1200 600)/90 Water-Bath Type Sterilizer Product Series.

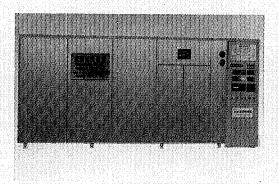

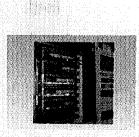

CHONGQING PHARMACEUTICAL MACHINERY FACTORY


Address:No.35,Shapingba Main Street, Chongqing, P.R.China


Tel:0086-23-65312352

Cable:1434


FAX:0086-23-65315646

Post Code:400030

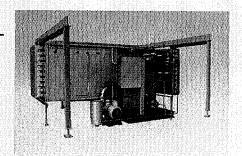

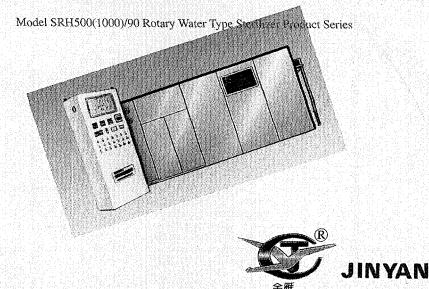

Ry Thermocompression Distiller

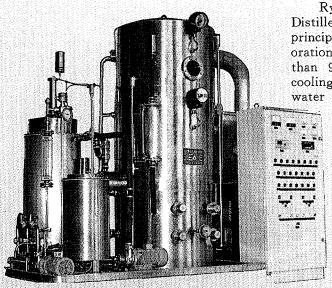

Ry Thermocompression Distiller is based on the principle of heat pump evaporation. It's saved more than 90% energy without cooling water. The distiller water is high quality, low

> cost and stable water quality. It conforms to the standard set by the People's Republic of China for use in medical injection. It's also suitable for use in chemical. petroleum, food, desalination of sea water, medical treatment and light industry.

Main Products

- 1) Air Lift Ferment Reactor
- 2) Thermocopression Distiller
- 3)Horizontal Closed Peeler Centrifuge 4)Doubble—cone Rotary Vacuum Dry-
- 5)Oscillating Granulator
- 6)Multi-effect Water Still
- 7)Glass lined Closed Reactor
- 8)Mixer

- 9)Powder Forming Machine
- 10) Sugar Coating Machine
- 11) Honey Refining Unit
- 12) Non standard Chemical Equip-
- 13) Motion Crystallizing Machine
- 14) Spraying Drying Tower
- 15) Auto Ion Exchange Water Treatment Machine

Brief Introduction

Baoji Pharmaceutical Machinery Works is one of the key manufactures of the State Medicine Administrative Bureau, specializing in producing pharmaceutical and chemical equipment. It has 30 years kindly service experience for state and abroad user in pharmaceutical industry. It also gets good name from user.

The works has chief production shops of riveting and welding glass lining . forging machine and heat treatment. It has a backbone of rich experience skilled business and excellent technology. It has perfect system of quality administration and advanged test means. The works research institute is in contact with the enterprises of the pharmaceutical chemical and college of the whole country. It devotes itself to the development of pharmary chemical machines and equipments new products and technique. It has put a lot of high technology products on the market.

Director:Li Guitang Post Code:721006


Add: 21 Chuan Shaan Road. Baoji. Shaanxi. Tel:0917-3315195 Fax. Cable: 2750

Shandong Xinhua Pharmaceutical Co. Ltd.

Shandong Xinhua Pharmaceutical Co. Ltd.(former Shandong Xinhua Pharmaceutical Factory) established in 1943, is a backbone enterprise of Shandong Xinhua Pharmaceutical Group Co. and one of national king signed enterprise. Mechanical plant with 3000 square meters in area and 15 million yuan of fixed assets is affiliated to the company.

The plant succeeded in making china's first glass-lined reactor in 1953, and was awarded with national creation & invention prize. It can manufacture 7 kinds of glass-lined equipments in 55 specifications, BRI class pressure vessels, preparartions facilities, and chinese medicine extractors with the licences for the manufacture of class I and class II pressure vessels and industry products.

The brand of product:Xinhua Legal representative:He Duanshi

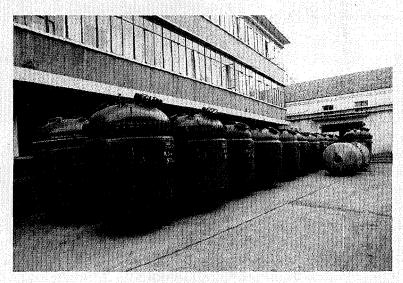

Authorized agent: Zheng Dezeng

Address: No.25 Xinhua street Zhangdian district Zibo city

Shandong province PRC.

Post code:255005 Telex:(0533)2185538

Tel:(0533)2184223-2340 or 2344



中国武汉制药机械厂

CHINA WUHAN PHARMACEUTICAL MACHINERY FACTORY



Wuhan Pharmacy Machinery Factory of China is one of the designated factories by the State Administration of Medicine for manufacturing equipment of traditional Chinese medicine.

The factory has been in service for pharmacy, chemical industry, foodstuff, health care, light industry and other trades for over 30 years. With powerful technical forces, advanced technology, fine processing equipment, reliable testing measures and complete quality assurance system, the factory integrates the design and production of pressure vessels, leading the tide in developing in research for traditional Chines medicine equipment.


Its "Jinguan" multi – functional extractor has been awarded with quality product titlis by the State Administration Bureau of Medicine in 1989, and its WC – 400 type pulverizer and PG – 26 spray drying unit have been awarded with prizes from the State Administration Bureau of Medicine and Wuhan City reapectively. Meanwhile , its "Dynamic extraction production line" filled the vacuum in China.

Quality first and client priority have been always the unique guiding principle of China Wuhan Pharmacy Machinery factory.


Add: 25Hangxiang (Red Lane), WuChang District, Wuhan Hubei Province, China.

Postcode: 430061 Cable: 1000


Tel: 027 - 8844815 8875662 8872022 8875663 Fax: 027 - 8874935

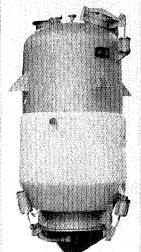





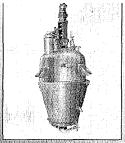

Multi-Function Extraction Tank

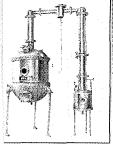


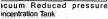

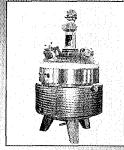


Preparing Tank



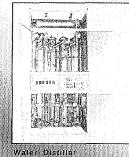

in Exchange Water Vessei

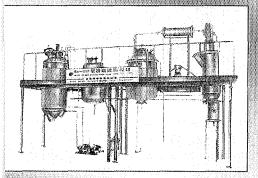

e Complete set of Equipment for Extracting and concentrating

Changshu Pharmaceutical Machinery Factory is one of the core enterprises which are under Jiangsu Changshu Pharmaceutical Group Co.,Ltd. It is an enterprise designed by the State Pharmaceutical Administration and the State Traditional Chinese Administration. The Factory has the licences of desibning,manufacturing and installing class A and class B pressure vessels. It mainly produces 4series of equipment, such as extractors, separators, concentrators, dryers, etc and over 100 varieties and specification in 10 categories of Pharmaceutical, food-making and chemical machinery.

Main products

| Product Name | Specifications |
|----------------------------------|-------------------------------------|
| 1 Multi-Function Extraction Tank | I~7M3 |
| 2 High-Effect Film Evaporator | 2.2M ² ~60M ² |

3 External Circulating Uni-Effect evaporator

350°2000kg/h

4 Coering-Type External Circulating Evaporator

700kg/h~1000kg/h

Tri-Effect Vaperizing Concentrator

500kg/h^2000kg/h

Vacuum Reduced pressure Concentration Tank

| | 100L 2000L |
|---------------------------------|----------------|
| 7 Spherical Concentration Tank | 250L*3000L |
| 8 Return-Flow Concentration Tan | k 5001,~10001. |
| 9 Tilting-Type Reaction Caldron | 501,75001, |

10 Sterilizing Dryer 11 Alcohol Recovery Tower Model 200-Model 800 12 The Complete set of Equipment for Extracting and

concentrating TN500-TN3000

13 Pipe-Standing-Type Multi-Effect Water Distiller

500kg/h,1000kg/h 14 Coil-Type Water Distiller 300kg/h,500kg/h 15 High-Purity Water Distiller 200kg/h,400kg/h 16 Ion Exchange Water Vessel 1000kg/h,2000kg/h 17 Injection Water Storage Tank 2.3°5M3

18 Preparing Tank 300L~2000L 19 Material Storage Tank $0.5^{\circ}10M^{3}$

20 Alcohol-Depositting Tank 300L, 2000L (Water-Depositting Tank)

21 Ultra-Sonic Cleaning Machine


Add: No.38 Yushan North Road, Changshu, Jiangsu Province,

P.R.China

Tel: 0520-2895982 2897666

Fax: 0520-2897000

P.C: 215500

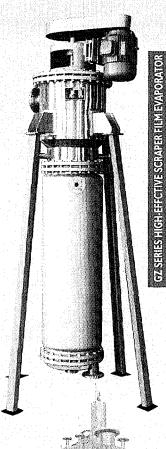


BRIEF INTRODUCTION OF ENTERPRISR

Wenzhou chemical machinery factory (another name: Wenzhot pharmaceutical machinery factory) was set up in 1966 and is the members of national Association of Pharmaceutical Machinery& Equipment Industry and the branch of National Special Pressure Container Techinical Committee with the qualification to manufacture BRI-grade pressure container(No.of manufacturing licence; RZZZhe039-99). It holds 16.86 million yuan of the fixed assets, 16, 170 sq.m of the factory land, 8,058 sq.m of the building area and 3.000T of the annual designing capacity for pressurecontainer; and in the past years, the products of pharmaceutica machinery were many times given with prizes such as the Secondary prize of scien tifical& technical Advance of Wenzhou City, the Silverbottle prize of Invention of Zhejiang province the Third prize of Scientifical&technical Advance of the National Medicine Bureau and the prize of Qualified New-product of the National Economica Committee.

with powerful technical force, perfect quality guarantee system and exquisite technology of processing, the product quality has been made good and stable for a long time, the equipments for raw material medicines, traditional Chinese medicines, preparations and pharmaceutical packings provided for pharmaceutical industry and the equipments of pressure containers provided for chemical, foodstuff, beer etc. industries are well sold to the 31 cities, provinces and autonomous regions of the nation and abroad, and also, the design, manufacture and installation for non-standard pressure containers are deeply in the users' good graces.

The factory always sticks to the motto of "paying much attention to quality and honouring all contracts and sincerely welcomes all users to give advice and to mutually create the brilliant future under the

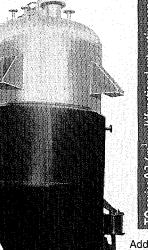



es Garostile, Riblet Painting Va

Add:No.106,Xueshan Rd,Wenzhou,China Tel:0577-8524295 9015271 9039421 Fax:0577-8525448 Cable:0209 Postcode:325005




遊

Ruian Pharmaceutic Machinery Factory Zhejiang China

Brief Introduction


Our plant was set up in1965 it's a member of associationin pharmaceuticeuticequip in China. It is entitled to bestrictly perform obligation bycontract's plantby province. Village&town's bureaus acknowledge 7 items of managment qualified. It has not only super production equip& perpect inspect means but abundant technique and powerful ability of research. We can design a full of equipment of China medicinematerial concentration, unusual equiment, pipe install for user.

SJZN-1000 Triparbite Energysaving Evaporating

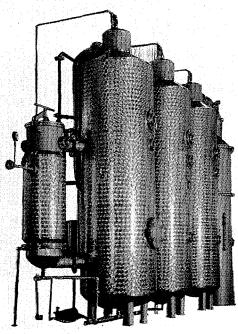

Other Products List


| name of article | model | remark |
|---|---|--------|
| Double-blade Mixer | WSH-150 | |
| Triparbite Energysaving Evaporating | SJZN-500 SJZN-2000 | Serial |
| Rolling style Automatic Bubble Packing Machine | DP-250 | |
| Flatbed Style Automatic Bubble Packing Machine | DPP-(140~250) | Serial |
| External Circulating Vaccum Evaporator | WZ-(500~1000) | Serial |
| Fast Setting Alcohol Recovery Tower | JH-(200~600) | Serial |
| Muliti — extractor | TQZO(3M3~6M3) | Serial |
| Boiling Box - desiccator | XFG ~ B | |
| Model – V High Efficiency Blender | GHV-(0.5 ³ ~2.5M ³) | Serial |
| Efficiency Scraper Film Evaporator | $0.8M^2 \sim 1.5M^2$ | Serial |
| Balling Concentrated Pot | NQ-(200~1000) | Serial |
| Double Up & Down Film Evaporator | 7M ² -25M ² | Serial |
| Efficiency Move Drawing Concentrated Machines | TQWZ-(1-200 ~3-500) | Serial |
| Boiling Water Rotary kill germ machine | 150 | |
| Concentrated Drawing Alcohol Machines | 300 ~ 1000 | Serial |

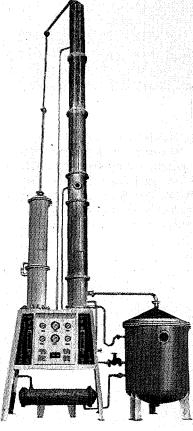

JH-Serial Alcohol Recovery Tower

Add:20 yanjiangxincun Ruian City Zhejiang China

Tel: 0577-5622267 5623309 5653288 Fax: 0577-5622267 Ptc: 325200

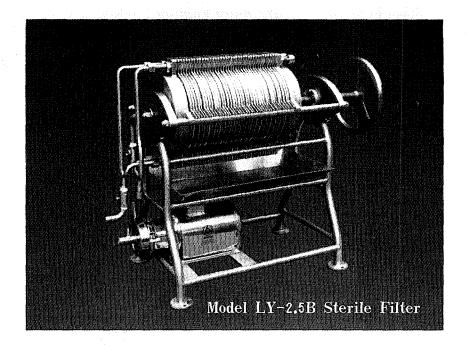

Brief Introduction

Located in Leshan City, Sichuan which is famous for its beautiful scenery and traditional culture, Long March Pharmaceutical Machinery – a subsidiary of Sichuan Long March Pharmaceutical Co., Ltd., is a manufacturer of different kinds of filters, bottle shakers, pressure containers and customized equipment. With the attractive scenic spots nearby including Leshan Grand Buddha and Mount E'mei, the plant is targeting to further develop and market its products supported by best quality and post-sales services, and to establish better reputation in both domestic and overseas market.

Main Products List

Model LY-2.5A(B) Sterle Filter Model ZP-96A(B) Bottle Shaker Model ZP-450 Bottle Shaker MF-B Serial Metal-Powder Air-Filter

Address:


Long March Pharmaceutical Machinery Sichuan Long March Pharmaceutical Co., Ltd. 120 Baiyang Road, Leshan Sichuan 614006, China

Tel: 0833 2440059 Fax: 0833 2440059 Legal Representative: Meng Dejun




WENZHOU ZHONGZHI MEDICINAL EQUIPEMENTS FACTORY

SJN-500B-8000B LATEST TYPE 3-EFFICIENT ENERGY-SAVING CONCENTRATOR. (ALCOHOL CA BE RETRIEVED)

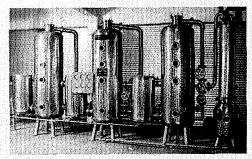

DG FOG-SPRAY DESICCATOR

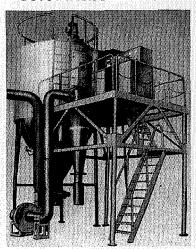

Our factory is presently by the largest manufacturing base for facilities for extracting, concentration, alcohol recovery, precipitation, filtration, bacterin killing and drying, etc. of herbal medicine and plant medicine. With brand new, advanced and energy conservative equipment, the factory owns high prestige, low price and immense customers. Especially Three Function Energy Conservative Concentrator, Multi-purpose Alcohol Recovering Concentrator, the best seller appeal to domestic customers who account for 70% of all customers. Equipment like dynamic extractor, single function concentrator, alcohol recovering tower, melangeur and alcohol precipit ator, spray dryer and filter, etc. have been exported to southeast Asian countries.

After the use of complete equipment by Thailand Shilefu International Co.,Ltd. they appraise our factory as low price excellent quality, high prestige and satisfying service.

Our main products are:DTOO. 5m³-6m³ dynamic extractor forced outer circulation extractor, SIN-300-800B double function, tri-function energy-concervative concentrator SIN-1000-10000 tri-function falling-film evaporator. DIN-500-2000A, B multipurpose alcohol recovering & concentrator. JC-500-6000 alcohol stir precipitator, PYG-3-150 one-step spray drying pelleter, etc.


Add:Shipu Ind. Zone, Yongzhong, 559 Airport Ave, Wenzhou
P.C.:325024 Factory Manager:Wang Chengdong
Tel:86-0577-6373357 86-0577-6371414 Fax:86-0577-6876357
Shift phone:86-01396636918 86-0577-9009087 86-0577-9100529

WUHAN MEDICAL APPARATUSES AND INSTRUMENTS PLANT

Wuhan Medical Apparatuses And Instruments Plant is one designated by the State Medical Administration Bureau to manufacture medical apparatuses and instruments and pharmaceutical equipment. It has a histry of more than 40 years.

The plant is well equipped with all kinds machines and testing apparatuses, as well as powerful technical team. It has been qualified to manufacture pressure containers. The medical apparatuses and instruments and pharmaceutical equipment are widely used in hospitals, pharmaceutical plants, beverage and food industry and scientific reseach organizations. It has gained good reputation among its customers all over the country. Its main products are medicine preparation machines, multifunctional extracting & concentractor machines, ointment batch pouring machines, steam sterilizer, distilled water machines in 5 kinds of more than 20 varieties.

Multi-Functional Extracting And Concentrate-Reclaiming
In Outer-Circulation Set

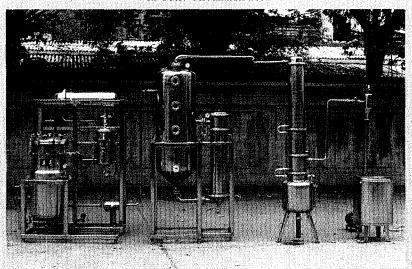

Add:NO.1 Xincun, Sanjiao Road, Wuchang, Wuhan, China. Tel:027-6811350 027-6811202 Fax:027-6811258


Cable: Wuhan6830 Post Code: 430062



PGL SPRAY DRYING GRANULATOR

(ONE-STEP GRANULATOR)

Patent No.ZL95240611.X

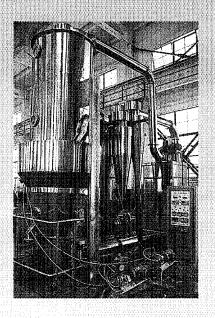

Mixing Spray drying Granulating Coating

Liquid ——Granules
All parts in contact with the product made of stainless steel
Good Manufacturing Practices

Data:


| Size | 3 | 5 | 15 | 30 | 50 | 80 | 120 |
|-------------------------|-----|-----|-------|-------|-------|--------|---------|
| Liquid kg/h | 2-4 | 4-6 | 10-25 | 20-40 | 40-75 | 65-105 | 105-150 |
| Steam consumption kg/h | 35 | 55 | 120 | 211 | 284 | 341 | 411 |
| consumpressd air m³/min | 0.5 | 0.6 | 0.6 | 0.9 | 0.9 | 1.1 | 1.3 |
| capacity kw | 6.7 | 10 | 10.7 | 13.5 | 21 | 32 | 41 |
| Electric heating kw | 6 | 9 . | 27 | 36 | 60 | 75 | 99 |

Chang Zhou One-Step Drying Equipment Factory

Address: JiaoXi, Chang Zhou, China, 213116

Tel: 0519-8902618 8902818 Fax: 0519-8902818 Cable: 7078



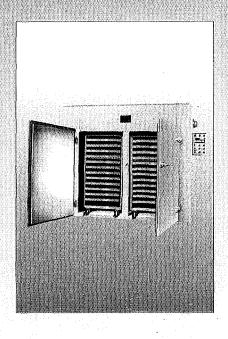

WON High Quality Prizes Province And ministry

Applications:

It suits to be used in pharmaceutical Food And chemical Industries For drying Wet Material

| Size | Rxh-14-b | Rxh-27-b | Rxh-41-b | Rxh-54-b | Rxh-5-c | Rxh-25-a |
|---------------------------|-----------------|-----------------|------------------|-------------------|------------------|------------------|
| capacity tkg/batch | 100 | 200 | 300 | 400 | 25 | 50 |
| steam consumptian kg/h | 20 | 40 | 60 | 80 | ,5. | 10 |
| wind volme m³/h | 1400 | 5200 | 9800 | 9800 | 3450 | 2500 |
| power kw | 111 | 1.1 | 2.2 | 2.2 | 0,45 | 0.45 |
| overoll dimensions mm | 2430×1200× 2375 | 2430×2200× 2433 | 3430× 2200× 2620 | 4380 ×2260 × 2620 | 1550× 1000× 2044 | 1200× 1000× 1600 |

ChangZhou One-Step Drying Equipment Factory

Address : JiaoXi, ChangZhou, China, 213116

Tel: 0519-8902618 8902818

Fax: 0519-8902818


Cable : 7078



GHL HIGH-SPEED MIXER GRANULATOR

Mixing: 2min/Granulating: 1-4min

Standard Features

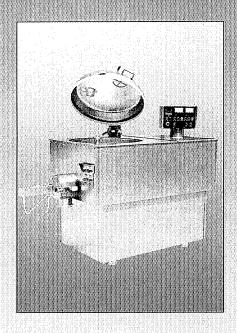

constant spiral product flow facilitated by the pear-shaped process
pan in vessel lid integrated
abservations window with wiper
chopper designed with stringed
crossscrew blades
additional stainless steel cover for
chopper motor protection
quick discharge by rotating blades
all product contacting parts are in
stainless steel, mirror polished
CGMP construction with air purge
sealing of drive shafts

Data:

| volume j | 10 | 150 | 200 | 250 | 400 |
|----------------------|----------|---------|-----------|---------|---------|
| capacaty kg/batch | 3 | 50 | 80 | 100 | 180 |
| agitator speed c.p.m | 300/600 | 240/360 | 180/270 | 180/270 | 106/155 |
| mator capacity kw | 1.5/2.2 | 6.5/8 | 9/11 | 9/11 | 18.5/22 |
| chopper speed c.p.m | | | 1500/3000 | | |
| mator capacity kw | 0.85/1.1 | 2.4/3 | 4.5/5.5 | 4.5/5.5 | 6.5/8 |


Chang Zhou One-Step Drying Equipment Factory

Address: JiaoXi, Chang Zhou, China, 213116

Tel: 0519-8902618 8902818 Fax: 0519-8902818 Cable: 7078



FL FLUIDIZED-GRANULATOR

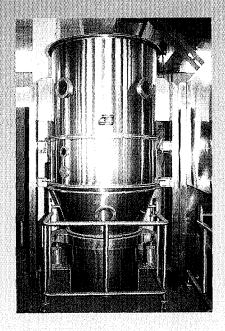

(ONE-STEP GRANULATOR)
Mixing
Granulating
Drying
Filmcoating

Export:

DaMei Medicine Industry USA ITK.CLT Japan etc.

Applications:

It is applicable in pharmaceutical, food and chemical industries to make granules and dry product. The granules can be used to press tablets or guick-soluble

Standard Design Features:

All part in contact with the product made of stainless steel.

Complete spraying device with coaxial air atomized spray nozzle and peristaltic pump. Separate fan unit
Observation windows
Pneumatic outlet air flap for control of fluidization height
Good Manufacturing Practices

| capacity:kg/bacth | 3 | 5 | 15 | 30 | 45 | 60 | 90 | 120 | 200 | 300 | 500 |
|-------------------|----|----|----|-----|-----|-----|-----|-----|-----|------|------|
| volume:l | 12 | 22 | 45 | 100 | 155 | 220 | 300 | 420 | 660 | 1000 | 1500 |

Chang Zhou One-Step Drying Equipment Factory


Address: JiaoXi, Chang Zhou, China, 213116

Tel: 0519-8902618 8902818

Fax: 0519-8902818

Cable: 7078

PRODUCTION LINE FOR MAKING TRADITIONAL CHINESE MEDICINE BY DYNAMIC WATER LEACHING PROCESS

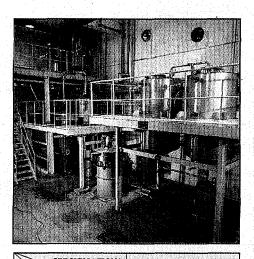

| SPECIFICATIONS | Raw material feed (kg/batch | | | | |
|---------------------------------------|--|--------------|--------------|------------|--|
| ITEMS PARAMETERS | 50 | 100 | 200 | 300 | |
| Method | Hot w | ater dyr | ımic lea | ching | |
| Raw material feed (kg/batch) | 50 | 100 | 200 | 300 | |
| Max.feed water leaching(L/min) | 500 | 1000 | 2000 | 3000 | |
| temperature for Leaching | 95℃±5℃ or 100℃±2℃ | | | | |
| Leaching time | 1.2 hrs./batch (time being calculated when water temperature meets requirement) | | | | |
| Solid-liquid separation Mode | Centrifugal | | | | |
| Mode for residue removal | Mechanical scraping | | | | |
| Mode for concentration | Vacuum co | oncentration | by natural c | irculation | |
| Max.concentration speed (kg/hr) water | 300 | 600 | 1000 | 1500 | |
| Mode for drying | Spray drying | | | | |
| Liquid feed for spray drying (kg/hr.) | 30 | 60 | 120 | 180 | |

STEPS:

- A.Dynamic water leaching:
- B. ThreeOstage Liquid—solid separation;
- C Low temperature vacuum concentration;
- D. Spray drying

Main functional parameters

ShangHai Jin Qiao Spray and Dry Complete Equipment Plant

Address: 60 Ying Jia Bang Yang Gao Rd Pu Dong ShangHai China

Tel:0086-021-58545772

Fax:0086-021-58995886

JIANGYIN FREEZE DRYING EQUIPMENT PLANT

We are specialized in manufacturing various kinds of freeze drying equipment, our products can be classified into three basic series, which contains of 16 detailed specifications.

 FREEZING DRYER FOR TESTING PURPOSE: 0.2m²~1m²


II. FREEZING DRYER FOR PRODUCTION PURPOSE: 2m²~10m²

III.LARGE-SCALE FREEZING DRYER, ESPECIALLY USED FOR FOODSTUFF INDUSTRY 20m²~100m²

We can also supply the following equipment:

- *Circulating hot air oven
- *Tunnel sterilizer
- *V-type mixer
- *Ultrasonic washing machine
- *Vacuum dryer
- *Alcohol recovering tower
- *ZS series vibrating screener
- *Mixing granulator

GDA Series Freeze Dryer Machine

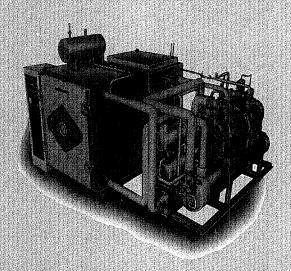

Add:West street, Shen Gang Zhen, Jiangyin,
Jiangsu, 214443, China
Tel:86-510-6621001,6621481
Attn.:Mr.Zhang Yao Fang

天利集团

北京天利医药设备有限公司

Pharmatech Group

Beijing Tianli Medicine Equipment Co., Ltd

Dear friends.


Beijing Tianli Medicine Equipment Co., Ltd., a subsidiary as a joint venture of Pharmatech Group. The TL-Lyoepic Freeze dryer which is produced first in China by introducing the leading foreign techniques, adopting imported parts with high quality, which receives high opinions in the related field and from our customers. So far, TL-Lyoepic $40(40m^2)$ freeze dryer which reaches the internationally leading level of the 1990s, has already been put into practical use in some of biggest national corporations e.g. NCPC, Harbin Pharmaceutical Factory and Harbin #2 Pharmaceutical Factory.

We are ready to offer good service to you.

Sincerely yours

甜放股

(Nie Qiming) General Manager

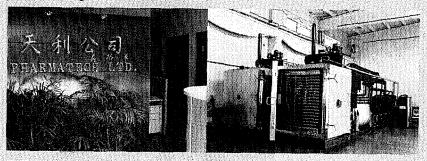



Pharmatech

Address:

香港: 九龙白加士街25号 庆云商业大厦12写字楼

Hong Kong: 12/F Hing Wan Commercial BLDG. No.25 Parkes Street, Kowloon

Hong Kong

TEL: 852-27303668 FAX: 852-27303668

Home Page: WWW.tl-pharmatech.com

北京: 朝阳区麦子店西路3号新恒基

国际大厦311室100016

Beijing: Rm. 311 Towercrest Plaza Maizidian West Rd., Beijing

P.R. China 100016


TEL: 86-10-64671989

FAX: 86-10-64671986


E-mail:sales.tl@tl-pharmatech.com



山东新华医疗器械集团有限公司 SHANDONG XINHUA MEDICAL INSTRUMENT GROUP CO.,LTD.

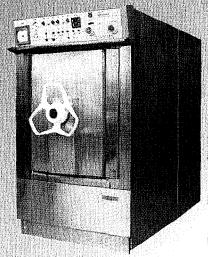

XG1 - P型脉动真空灭菌器 XG1 - P Pulse Vacnum Sterilizer Esterilizador vacio pulsatil de xg1.p

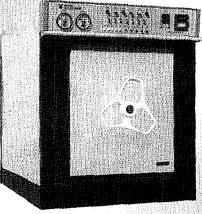

XG1·k型快速冷却灭菌器 XG1·K Quick Cooling Sterilizer Esterilizador con enfriamiento rapido de XGLK



公司地址: 山东淄博市张店区洪沟路2号

邮编: 255027 国际贸易部

电话: 0533-2184349 2183171-343

传真: 0533-2184349

Add: 2Honggou Road, Zhangdian District, Zibo. Shandong Proy, P. R. CHINA Post Code: 255027

INTERNATIONAL TRADE DEPARTMENT Tel:0533-2184349 2183171-343

Fax:0533-2184349



YANGCHUNHUI SCIENCE & TECHNOLOGY DEVELOPING Co. Lt.d.


This company is a hi-tech. industrial company, which specializes in developing and manufacturing sterilizer and drier with electromagnetic wave.

Main products are GM series sterilizers for traditional Chinese medicine (foods). Its patent No. is ZL932i6624.5, whose sterilizing source is electromagnetic field with low temperature, atmospheric pressure. When the wonderful sterilizing is finished, it is sure that the effective component isn't damaged. Specially, it has absolute sterilizing effect on pathogenic bacteria.

It is operated easily and runs smoothly and reliably. Effect is excellent, which is in line with the GMP standards.

This machine was awarded the golden prize of "National carefully-chosen patent products castle" in 1993 and the first grade "Golden bridge" achievement prize of S.& T. by Ji-lin province.

The manufacturing machine is supervised by CAPE.

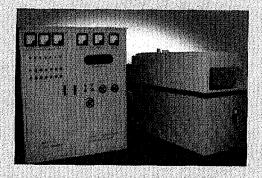

Address:No. 73 Shuang-Quan-Pu

Hai Dian district

Beijing, China

Liaison officer:Mr. Gao Xiao

P.C.:100085

Fax:0086-10-62329131

Tel.:0086-10-1391388610

COMBINED MACHINE FOR CLEANING RUBBER SHOPPERS OF PHIALS




Combined Machine For Cleaning Rubber Stoppers Of Phials. It finds its application in the pharmaceutical industry for producing various injections. While rinsing, sterilizing and drying rubber stoppers, it meets the requirements of the GMP. As a product recommended by the State Bureau of Medicine Administration, it is welcomed by our customers at home and abroad.

HUXI IRON WORKS

ADD:168/185 PINGTANG ROAD SHANGHAI TEL:021-62394624 FAX:021-62394624 P.C:200335

DOUBLE CONF SELVED ROTARY VACUUM DRYFR

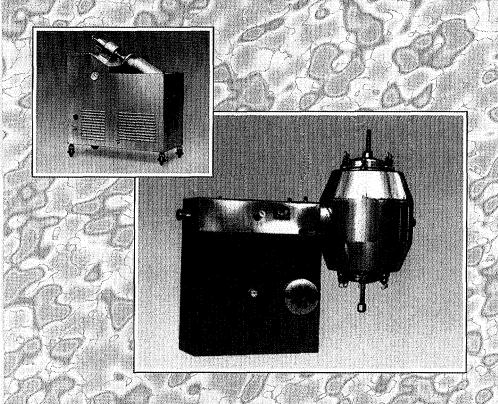

Double Cone Shaped Rotary Vacuum Dryer. It meets the requirements of bacteria-free operation in medicine, chemical and food industries and finds its application in refining small and roundish particles, e. g. antibiotics, crystals. In comparison to the conventional vacuum dryer, it raises efficiency by 200-300%.

HUXI IRON WORKS


ADD:168/185/PINGTANG ROAD SHANGHAI TEL:021-62394624 FAX:021-62394624 P.C:200335

MULTIFUNCTIONAL FULLY-AUTOMATIC ULTRASONIC RUBBER-CORK WASHER

DQX-1B

PATENT PRODUCT



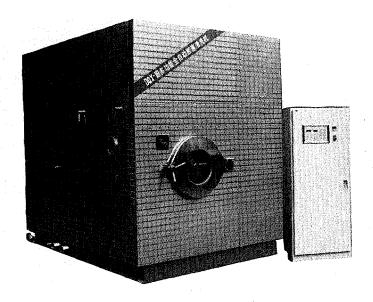


ZHENJING WENZHOU YAGUANG MACHINERY MANUFACTURING CO., LTD.

Add:99 Canghe Rd.(W). Yongzhong Town. Wenzhou, P.R. CHINA

Tel:0086-0577-6372550 Fax:0086-0577-6373389

P.C:325024 General Manager: Chen Guo Hua

BAOJI JIAN HUA PHMARMACHINERY CO.LTD

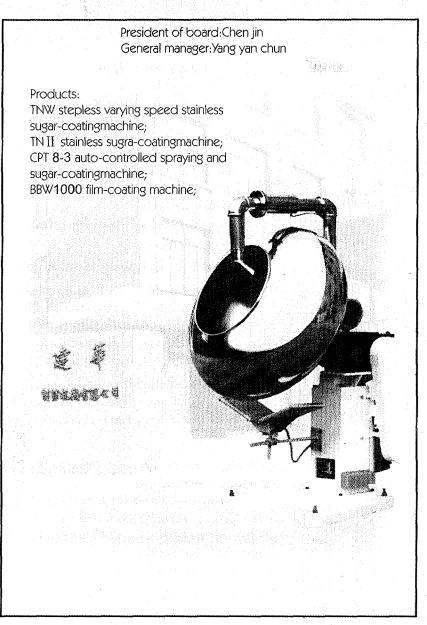

ADD: No.75 jianhua road,baoji,shaanxi TEL: (0917)3554670 (0917)3553721


FAX: (0917)3554670

P.C: 721001




BEIJING TONG-DA HI-TECH COMPANY


The 13th Institute,
China Academy of Launch-Vehicle Technology,
China Aerospace Corp.

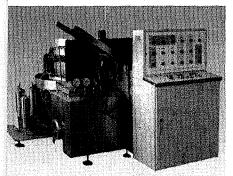

Beijing Tong -Da Hi-Tech Company is founded on t Launch -Vehicle technology, China Aerospace Corras its foundation and develops new and high technolog


The company develops, manufactures and manage pharmaceutical - machinery, medical apparatus ar Petroleum prospecting instrument, computer application develops and manufactures new products, innovates a train workers, process material for other consumers a are main products:

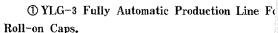

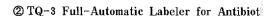


③ XP-3 Ultrasonic Cleaning Machine

The machine runs smoothly and reliably in steple bottles/minute in rated. The dirt and filth in the bottle ϵ smashed by high power ultrasonic device. The bottles ϵ and turned over two times automatically.

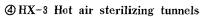

The heat area temperature can reach 350°C and a to sterilize. The clean degree in tunnels is as high as

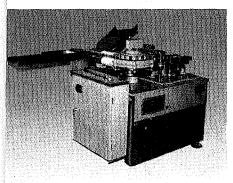

6

⑤GS-3 fully Automatic volumetric filling machine with seali machanism for roll-on caps. Advanced photo-electricity control and PLC contorl.

Accurate filling in fixed volume by six pumps.


Sealing the caps are finished by six 3-talon rolling heads.

Stop pouring for absent bottle.

The rated output is 8000 bottles/hour with rate of broken bottle ≤ 2

© Full-Automatic Labeler For Oral Liquid Bottle

It is suited to label and print production date and invalid date ampoules and other cylindrical containers with diameters from 10 min. to 23mm max. and with a cylindrical body height from 22 min. to 130mm max. It can stop labeling for absent bottle and has s protection set.

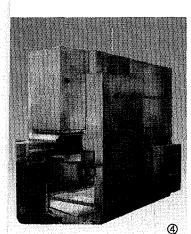

Address: No.52 yong-ding road Post address: P.O.Box 3913.Beijing.P.R.china Cable: 5405 (Beijing) Post code:100854 el: (8610) 68388768, 68388767

3th institute, China academy of he company takes space technology and products.

ax: (8610) 68388652

echano - electronic products, such as strument, Precision instrument, and so on. The company designs, ansfers technology. The company can tegrate with other units. The following


illing Oral Liquid Sealing With

ttles.


eed regulating and washes 133 hering to the bottles are shaked out and oured three times, blooded three times

ual, circulating hot air makes sure



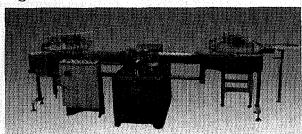




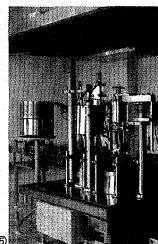


TD fully automatic production line for filling lyophilizing biologicals

This production line consists of the fully automatic volumetric filling machine ,rubber plug plugging machine and the sealing cap machine. The rated output is 200 bottles/min . Filling volume is adjustable from $1\sim 8 \mathrm{ml}$.

Screwing caps tightly. It is in line with the National Stand.. Automatic safe protection set.

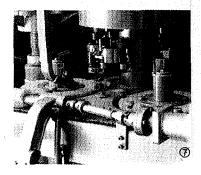

The 15th Institute of China Academy of Launch-Vehicle Technology Beijing Tian Min Hi-Tech Development Co.

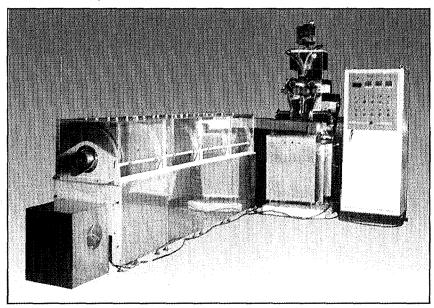

RJNJ-2 ROTARY DIE ENCAPSULATOR

We also supply the following equipment:

- BZJ-1000F GRANULATING AND COATING MACHINE
- BZJ-360M GRANULATING AND COATING MACHINE
- DWJ-1 AUTOMATIC SEAMLESS ENCAPSULATOR

Add: Donggaodi, Yongdingmen Wai, Beijing China


P.O.Box: 9200-71 Beijing, China

Fax: 86-10-68381851 Tel: 86-10-68382704

Zip: 100076


SHK Series High Speed Mixen/granulator

This is an equipment for tablet production in pharmaceutica industry, which can prepare ideal particle raw material for medicine. The particles produced by this equipment have the advantages of perfect sphericity, good fluency, and are homogeneous. The productivity of the equipment is high and the tablets produced by the equipment have excellent quality, smooth surface, better hardness, crumbliness and liberation ability. The equipment is optimized in design and has fine workmanship, carefully designed interlocking mechanism. In addition it is provided with cleaning function. Production cost can be saved greatly by using this equipment.

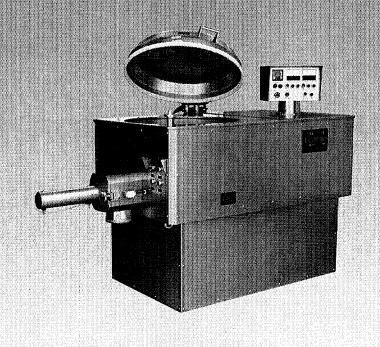



The advanced spraying mist system

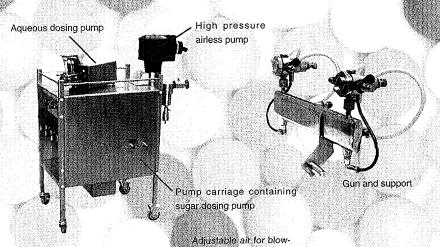


Spraying Gun

The spraying gun is controlled by three individual ways of compressive air: controlled air, adjustable atomization air and adjustable air for blowing atomization liquid into flat shape. The controlling the gun switch. When turning off the gun, the pin of the gun protrudes to prevent from the gun blockage.

The interval time brtween the gun turning off and another gun


ing the atomization liquid into flat shape Atomizing solvent

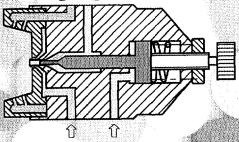

Adjustable atomiza- Controlled air tion air

turning off can be controlled by program or by hand. The adjustable atomization air can control atomization of the medium.the adjustable air for blowing atomization liquid into flat shape can control the size of the atomization sector angle. The spraying gun has thoroughly solved the problem of the gun blockage. The single way air can be individually controlled and adjustable, so the coating mediums with different consistencies, different viscosities and different solvents can gain ideal sector angles and atomization efficiencies.

北京航空工艺研究所

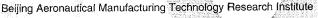



地 址: 北京市朝阳区八里桥北 信 箱: 北京市 340 信箱 201 部

多相:北京市340倍相2078 据:100024

电 话: (010) 65797893 传 真: (010) 65762306

电子信箱: bamtri @ public 3.bta.net.cn


Location: Baliqiao, Chaoyang District,

Beijing, 100024, P.R. China Address: P.O.Box 863 Beijing, P.R. China

Tel: (010)65797893 Fax: (010) 65762306

E-mail: bamtri @ public 3.bta.net.cn

BGB Series High Efficiency Coating Machine

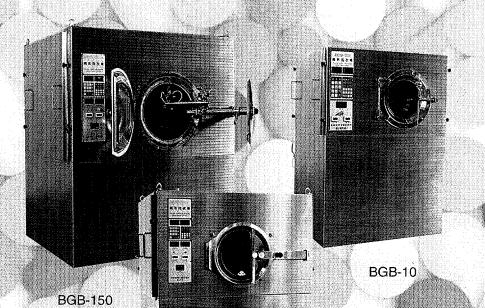

BGB-75






BEIJING HONG HUA ELECTROMECHANICAL NEW TECHNOLOGY COMPANY

Brief Introduction


Beijing Hong Hua Electromechanical New Tech. Company is sponsored by 206 Institute of China Aerospace Corporation and has registered in the Beijing New Technology Development about the Company which combining technology, industry, production and sale together, has the qualification of enterprise legal person.


In recent years, the company had sold 1300 sets of pharmaceutical equipment to 300 factories in 28 provinces of China, including TaiWan, and has exported its equipment to Indonesia, Thailand, Malaysia, the Philippines, the United States, Huangary, and HongKong.

The company has been elected as advanced unit for development of pharmaceutical machine, "50 outstanding enterprise" by the Beijing New Technology Development and Experimental Zone, and has also been honored the title of " Aero-Space Star ".

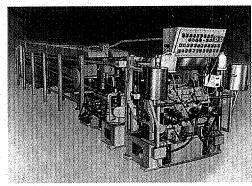

JSZ- III Production Line for Hard Capsule

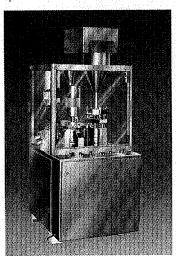

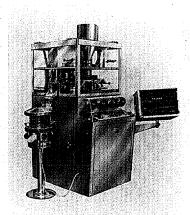


Series Products of ZTJ-400, 700, 800A Fully Automatic Capsule Filling Machine


Series Products of GZP-28, 28A, 28B Fully Automatic High-Speed Tablet Press

Company's name: Beijing Hong Hua Electromechanical New Technology Company Address: 50 Yongding Road, Haidian District, Beijing Postal address: P.O.Box 142-208, Beijing, China

Fax: 010-68386763 Tele: 010-68387242 010-68385234 Zip code: 100854



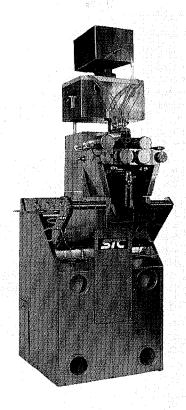

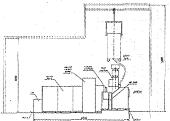

RG2 - I ENCAPSULATOR

RG2 - I a well designed rotary die encapsulator, with the cooperation of injector, spreading chamber and die to from all kinds of oil and liquid into variform soft capsules in pharmaceutical food and cosmetics fields.

Adopting a frequency inverter to regulate the die rollers steplessly, a special high precision temperature control system. It can complish successive capsules producing as the skin mading and forming at the same time. It provides high production efficiency and steady resional quacity to meet GMP standward and reached advanced level in the word.

Main Technical Specifications

| Stepless Speed Regulation | 0 ~ 5 r/min | | | | |
|---------------------------|-----------------------|--|--|--|--|
| Rated Rotation Speed | 3.5 r/min | | | | |
| Die Dimension | ø103×152mm | | | | |
| Material Supplied | 0 ~ 2ml | | | | |
| Dosage Filling Tolerance | ±2% | | | | |
| Motor | 1. 1kw 380V 1400r/min | | | | |
| Electric Heater | 150W × 6pc. 220V | | | | |
| Electric Heater Sleeve | 55W × 2pc. 220V · | | | | |
| Dimension | 750 × 900 × 1900mm | | | | |
| Weight | 900 (kg) | | | | |

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd.


Add: 1380 Kang Ding Road, Shanghai, China.

Tel: 021 - 62583414 021 - 62560800

Fax: 021 - 62582067 Cable:

Cable: 2771 Post Code: 200042



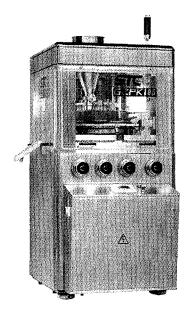

GZPK100

HIGH SPEED ROTARY TABLET PRESS APPLICATIONS:

This kind of machine can be used not only in pharmaceutical industry but also in chemical, food, electonic, daily necessities, is controlled automatically with computer, etc industry, the machine is designed for pressing tablet from various kinds of granular materials, It can press not only round tablets but also different geometrical shape tablets, Annular tablets. The tablets may have the impressed letter on both sides. We will supply special design for customers.

The design idea and the quality of manufactoure, STC company have achieved a new level The GZPK100 series high – speed rotary table press has adopted the newest technological achievement on mechanism and electricity that results into high producing capacity perfect quality of tablets, accurate weight of tabelts, lower waste, meeting GMP standrad.

Main Technical Specifications

| Model | | GZPK138 | GZPK132 | GZPK126 |
|-------------------------|-----------|-----------------|-----------------|-----------------|
| Number of stations | | 38 | 32 | 26 |
| Tool type | | IPT19 (EU19) | IPT19 (EU19) | IPT1" (EU1") |
| Max. table diameter | mm | 13 | 18 | 25 |
| Max. depth of fill | mm | 18 | 18 | 18 |
| Turret speed | r/min | 100 | 100 | 75 |
| Max. output | tablets/h | 228, 000 | 192, 000 | 117, 000 |
| Max. operating pressure | KN | 80 | 80 | 80 |
| Max. pre - pressure | KN | 20 | 20 | 20 |
| Power | ĸw | 7. 5 | 7. 5 | 7.5 |
| Weight | kg | | 2300 | |

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd.

Add: 1380 Kang Ding Road, Shanghai, China.


Tel: 021 - 62583414

021 - 62560800


Fax: 021 - 62582067

Cable: 2771



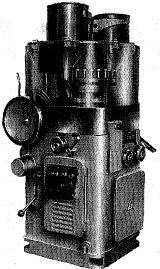

ZPW21A **ROTARY TABLET PRESS**

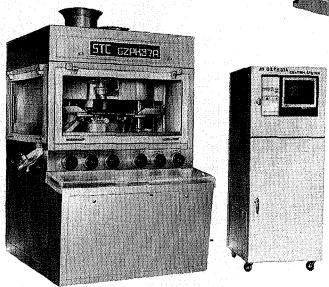

GZPK37A HIGH SPEED ROTARY TABLET PRESS

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd.

Add: 1380 Kang Ding Road, Shanghai, China.

Tel: 021 - 62583414

021 - 62560800

Fax: 021 - 62582067

Cable: 2771 Post Code: 200042



ZP100 SERIES ROTARY TABLET PRESS

ZP100 Rotary Tablet Press is a serial product widely used for pharmaceutical, food industry to produce all kinds of round tablets and irregular tablets.


This series meets the GMP standard, Pressing chamber is utterly separated with the gear system There is no any paint in the chamber The parts contacting with tablets are made of stainless stee or the surface of those parts are specially handled. The parts are poisonousless and rotten - resistant, The chamber is also sealed with or ganic glass to prevent pollution.

The gear system, the trach and the major friction parts are fully lubricated by the continual and concentrated lubrication device of the ZP100.

Being accordance with it can be equipped with forced feeding device to ensure feeding of powdery material customer's demand.

It's very convenient to observe, operate the pressure, thickness feeding, speed adjustable hand wheels and all kinds of buttons, indicating lights and indicators are at the same side.

The design of ZP100 series is reasonable. Most of parts are universal in this series. So it's easy to manage and change the spare parts The punch and dies meet the IPT standard.

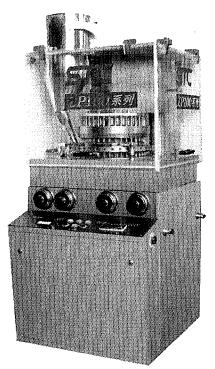

Main Technical Specifications

| Model | | ZPY136 | ZPY129 | ZPY124 | ZPY134 | |
|----------------------|---------------------------|-----------|--------|--------|--------|--|
| Station | | 36 | 29 | 24 | 34 | |
| Max. Pressure | (kN) | 80 | 80 | 80 | 80 | |
| Max. Depth of f | Fill (mm) | 18 | 18 | 18 | 18 | |
| Max. Dia. of Ta | ablet (mm) | 13 | 18 | 22 | 13 | |
| Production Cap | acity (ten thousand/hour) | 9. 7 | 7.8 | 6.4 | 10 | |
| Motor (kw) | | 4 | | | | |
| High (mm) | Without Material Vessel | 1,580 | | | | |
| With Material Vessel | | 1,800 | | | | |
| Cover Space | (mm) | 930 × 950 | | | | |
| Weight | (kg) | | 1, 4 | 400 | | |

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd.

Add: 1380 Kang Ding Road, Shanghai, China.

Tel: 021 - 62583414

021 - 62560800

Fax: 021 ~ 62582067

Cable: 2771 Post Code: 200042









HLSG250 MIXING GRANULATOR


Main Range of Application:

The process of Mixing, Wetting and granulating can be completed in single working period. It is one of the advanced facilities which meets the GMP requirement and has the characteristics of high efficiency, best quality, low consumer, none pollution as well as safety, etc.

- 1) Adoped wholly closed case, all the parts on machine are made of stainless steel.
- 2) The application of the conical barrel. Which the upper mouth is smaller than the lower botton and the two mixing arm guaranted the sufficient revolving of material. When the material is sufficiently mixed, the high speed revolving of smash cutter made the material cutted to be even and fine franule
- 3) The Consumption for adhesive is 15 20% less than the traditional technology.

Main Technical Specifications

| Bowl Volume | 250 L |
|----------------------|-----------------------|
| Working Volume | 220 L |
| Batch Size - Wet mas | 100 kg |
| Time per batch | 8 – 10 min |
| Mixing Power | 11/14 kw |
| Mixing Speed | 180/270 r/min |
| Chopper Power | 4,5/5.5 kw |
| Chopper Speed | 1500/3000 r/min |
| Granule Size | 20 - 80 mesh |
| Working Air Pressure | 0.5Mpa |
| Dimension | 2272 × 1050 × 1515 mm |
| Net weight | 1200 kg |

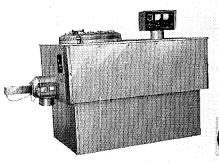



GFG80 120 EFFICIENT BOILING DRYER

This dryer is composed of dryer air filter, air heater, draught fan, and automatic temperature controller.

Its distinguishing features are: high drying speed and uniform temperature; hermetic and unpolluted draught after being filtered; and simple operation and easy cleaning. It is therefore, an advanced drying equipment conforming to G. M. P. requirements.

Main Technical Specifications

| ìtem | | Type | GFG120 | GFG80 | | | |
|---------------------------|-------------|----------------|-------------------------|--------------------|--|--|--|
| Product | ivity | | 120 Kg | 60 Kg | | | |
| Drying Intake Temperature | | | 50 ~ | T001 | | | |
| Compre | ssed Air | Pressure | ≥0. | 4Mpa | | | |
| Draught | Fan Vol | ume | 3400 m ³ /hr | | | | |
| Draught Fan Power | | | 11 Kw 380V | | | | |
| Mixer Motor Power | | | 0.55Kw 380V | | | | |
| Filler - 8 | ag Vibrati | on Motor Power | 0. 12Kw 380V | | | | |
| Mixer R | otating S | peed | 11 rpm | | | | |
| | Steam | Pressure | ≥0. | 4mPa | | | |
| Heater | Steam | Consumption | 160 | Kg/h | | | |
| | Electricity | | 40 Kw 380V | 30Kw 380V | | | |
| Overall Dimensions | | | 3020 × 1600 × 2950 | 2735 × 1390 × 2395 | | | |
| Weight | | | . 1750 Kg | 1200 Kg | | | |

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd.


Add: 1380 Kang Ding Road, Shanghai, China.

021 - 62560800 Tel: 021 - 62583414

Fax: 021 - 62582067

Cable: 2771



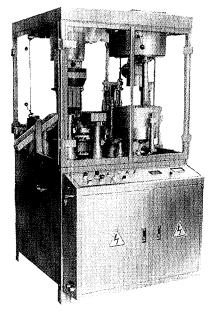

NJA 330 AUTOMATIC CAPSULEFILING MACHINE

It is designed with automatic intermittent rotating motion, every cycle of motion there are 8 times of dwell during which the various operation (feeding, opening, filling, faulty capsule ejection, closing, full capsule ejection, and moules cleaning) take contemporarily place.

It is suits to fill chinese and western medicines and the filling dosage can be adjusted correctly.

It is equipped with A.C frequency converter and LCD.

The work domain of the machine is sealed with transparent polymethyl methacrylate. The parts contact the capsules and madicines are made of stainless steel.

It is an ideal equipment for all pharmaceutical factories to produce hard capsules.

Main Technical Specifications

| No. 1000 To 1000 | ical Specifications | | | | | |
|------------------|---------------------|--|--------------------|--|--|--|
| Model | | NJA330 | NJP1000 | | | |
| Output | | 20000n/h | 60000n/h | | | |
| Suitable Cap | sule | No 0、1 | .2.3.4 | | | |
| Capacity Diff | erence | ± 5% | | | | |
| Main Power | | 1. 28kw | 4. 05kw | | | |
| Dimensions (| w×D×H) | 1015 × 1000 × 1850 | 1270 × 1060 × 1960 | | | |
| Weight | | 1000kg | 1500kg | | | |
| Annandia | Vacuum Pump | XD - 040 1. 1kw 30m ³ /h 0. 4 ~ 0. 6MPa | | | | |
| Appendix | Vacuum cleaner | XCJ 36 1.1kw | | | | |

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd.

Add: 1380 Kang Ding Road, Shanghai, China.

Tel: 021 - 62583414

021 - 62560800

Fax: 021 - 62582067

Cable: 2771



ZP35A ROTARY TABLET PRESS

Application:

This machine is a double-press type, continuous automatic tablet press for pressing granular raw materials into tablets. It is mainly used in pharmaceutical industry and also the chemical, food, electronic industries.

Characteristics:

The outside part of the machine is fully-enclosed, The material of the cover and the inside table face is stainless steel. The turret face is specially disposed keeping the lustre of the face and prevent the cross pollution, meeting the G. M. P. requirement

It has transparent windows so that press condition can be observed clearly and the window can be opened cleaning and maintenance is easier.

The lay-out of all the controllers and the operate parts is reasonable.

The speed adjustment of the driver is controlled by the frequency conversion device.


All the drive devices are located inside the machine so keep the machine clean.

The hydraulic system can assure the stable of the system.

A overload protection unit is included in the system to avoid the damage of the punches and apparatus, when overload occurs, the machine stops automatically.

This machine also has electronic magic driver and other safety protection device which can be adjusted and operated when running.

The punches and dies are also applicabel to Model ZP19 ZP33 ZP33B and ZP35.

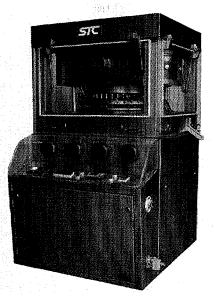

Main Technical Specifications

| Dies | 35sets | | |
|--------------------------|------------------------|--|--|
| Max. Pressure | 60KN | | |
| Max. Dia. of Tablet | 13mm | | |
| Max. Depth of fill | 15mm | | |
| Max. Thickness of Tablet | 6mm | | |
| Turret Speed | 14 – 36r/min | | |
| Production Capacity | 150000 pc/h | | |
| Motor | 3KW960r/min/380/50v/hz | | |
| Overall Size | 910 x 1120 x 1650mm | | |
| Net Weight | 2300kg | | |

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd.

Add: 1380 Kang Ding Road, Shanghai, China.

Tel: 021 - 62583414

021 - 62560800

Fax: 021 - 62582067

Cable: 2771



The Brief

Shanghai Tianxiang & Chentai

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd. is a joint venture between Shanghai No. 1 Pharmaceutical Machinery Factory belonging to Shanghai Light Industry Machine Inc. and Singapore Chentai Machinery Co. Ltd. The registered capital is five million one hundric eighty thousand U. S. dollar. The total investment is nine million fifty-eight thousand U. S. dollar.

Our company is specifized in producing pharmaceutical machinery equipment and releval chemical machinery equipment not only in pharmaceutical industry but also in chemical, food, elektronic, battery, metallurgical industry. The trade mark"天祥""STC" are generally famous brand both home and abroad. Up to now the company has long history of 60 years and become one of the biggest profession enterprises in China. Our products cover twenty-nine provinces of China, the sands of customers and have good reputation in the market. The products have been exporting sin 1950's, selling well in South Asia and openning American and Eurpoean Market. The company is the best enterprise making foregin currency among the same industry.

The company has had more than 30 kinds of products and has developed a seres of GN meeted products, such as ZP35A, GZPK37A ROTARY TABLEYT PRESS, GZPK100 SERII TABLET PRESS, ZP100 SERIES TABLET PRESS, DPP205 BLISTER PACKING MACHIN

| Series No. | Process of Craft | Available Equipment | Series No. | Process of Craft | Available Equipment |
|---------------|----------------------------------|--|---------------|----------------------------------|--|
| 1 | Pulverizer | GF300 Efficient universal pulverizer | 3 | Mixer Machine & Granulator | ZL250A Rotary granulator |
| | | GF300A Efficient universal pulverizer | | | HLSG250 Mixing granulator |
| | | GF300A2 Efficient universal pulverizer | | | HLY10 Mixing granualtor |
| | | SGF130A High - speed pulverizer | | | HLY100 Mixing granualtor |
| | | DF20 Efficient universal pulverizer | | | HLY200 Mixing granualtor |
| 2 | Vibrating Sreens | XZS400Vibrating sreens | 4 | Drier | GFG80 Efficient boiling dryer |
| | | XZS500Vibrating sreens | | | GFG120 Efficient boiling dryer |
| | | XZS800Vibrating sreens | 5 | Granulate | ZLY140 Rotary granulator |
| 3 | Mixer Machine & Granulator | CH10 Trough type mixing machine | 5 | Regulator | ZLY180 Rotary granulator |
| | | CH50 Trough type mixing machine | 6 | V – type | 0.05M - 3MV Mixer |
| | | CH150 Trough type mixing machine | | Mixer | Vacuum feeding device of v - type mixe |
| | | CH150A Trough type mixing machine | 7 | Tablet Press | GZPK37AHigh speed rotary tablet press |
| | | CH200 Trough type mixing machine | | | GZPK100High speed rotary tablet press |
| | | CH200A Trough type mixing machine | | | ZP100 Series rotary tablet press |
| | | LYK60 Oscillating granulator | | | ZP35A Rotary tablet press |
| | | LYK160 Oscillating granulator | | | ZP35 Rotary tablet press |
| | | ZL250Rotary granulator | | | ZPW21A Rotary tablet press |
| | | YK160A Oscillating granulator | | | ZPW21 Rotary tablet press |

OUR DEMAND OF PRESENT AND FUTURE

ntroduction f

harmaceutical Machinery Co., Ltd.

ILSG250 MIXING GRANUATLOR, BGC1000 FILMCOATING MACHINE, Now the company have had ablity to produce completely line in solid medicine line from pulverizing, screening, mixing, granulating, lrying, mixing, pressing to filmcoating, or sugarcoating, packing.

In 1997 the company is the first enterprise passed examination of ISO9001, meanwhile, it is also omfirmed as an advanced technology enterprise and a model enterprise praised by Shanghai muicipal administration using CAD.

By the spirit of "Unity and Struggle, Striving for perfection, Pioneering. Spirit, writing the new istory", we welcome all the customers, business patners to visit our company.

Address: 1380 kand Ding Rd. Shanghai P. R. of China

Zip: 200042

Tel: 021-62560800 62583414 Fax: 021-62582067 62583648

· Manager: Mr. Ma Ling Geng Mr. Yan Guo Jing

| ocess Craft | Available Equipment | Series No. | Process of Craft | Available Equipment |
|----------------|---|---------------|---------------------|---|
| et press | ZPY15 Rotary tablet press | 9 | Coating | BGC1000High efficient film coating machine |
| | ZP15 Rotary tablet press |] " | Machine | BYG1000/800/400 High efficient film coating machine |
| | ZP25 Rotary tablet press | | | DPT140Roller blister type packing machine |
| | ZP33B Rotary tablet press | | | DPH140Roller blister type packing machine |
| | ZP33 Rotary tablet press | | | DPP250 Plate blister packing machine |
| | ZP19 Rotary tablet press | | Packing Machine | DPP200 Plate blister packing machine |
| | THP Flower basket type tablet press | | | DPP140 Plate blister packing machine |
| | TDP Single punch tablet press | | | DPP100 Plate blister packing machine |
| | | | | DPP205 Plate blister packing machine |
| | CBN20Semi - automatic capsule filling machine | 10 | | DLB240 Twin - foil packing machine |
| | NSA330 Automatic capsule filling machine | | | DLB120 Twin - foil packing machine |
| osule | YDL60 Capsule printer | | | SP Capsule and tablet counting machine |
| | NJP1000 Automatic capsule filling machine | | | FSZ Filling paper machine |
| er & sher | ZWS137 Horizontal rotary tablet sifter | | | FXG Lids turner |
| | ZS100 Round brush tablet sifter | | | FR Ailfoil sealer |
| | MJP2 Capsules tablets polisher | | | DCF200 Pillets filling and packing machine |
| ating | BY300/400 Sugarcoating machine | | | DCF200 Powder filling and packing machine |
| nine | BQ1000/800 Sugarcoating machine | 11 | Encapsulator | RG2 - I Soft capsule encapsulating machine |





DPP205 Plate Blister Packing Machine

Application:

This machine is used to pack all kinds of tablets, sugar-coated tablets and capsule etc, with PVC and alumi-nium material in pharmaceutical, food industries. Characteristics:

This machine adopted advanced technique, owns compact structure, reliable function, highlevel quality.

Shaping, sealing, batch number printing, transfering, indenling, cutting, etc, are controlled by intensive position transducer and programmable controller, so that positioning is accurate.

It has overload protection of the motors and overtemperature protection function. The emergency brake in case of lacking materials are controlled by eletropneumatic control, and the mechanism will stop at a fixed and safe position which will ensure the machine keeping the best condition

The jog-run function will facilitate the adjustment and the testment.

The plastic blister formed by compressed air, so the wall-thickness of the blister is uniform and the surface is smooth.

The die can be made according to various customer's requirement.

This machine is drived by machanical infinitely variable-speed driver, so action is reliable and accurate, operation issimple.

Pneumatic system is smoothed by no mist lubrication.

Various feeding structure can be provided and matched with the main part to suit various medicine to be packed.

Main Technical Specifications

| Packing Speed | 20 – 50 1/min | |
|-------------------|----------------------|--|
| Max. Forming Area | 190 × 98mm | |
| Blister Depth | 1 ~ 12mm | |
| Max. Power | 5kw | |
| Overall Size | 3000 × 1030 × 1420mm | |
| Weight | 1000kg | |

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd.

Add: 1380 Kang Ding Road, Shanghai, China.

Tel: 021 - 62583414 021 - 62560800



BRIEF INTRODUCTION OF SHANGHAI YUAN DONG PHARMACEUTICAL MACHINERY GENERAL FACTORY

Being a professional equipment manufacturer for the pharmaceutical industry in China, this factory is provided with a strict production control and tight quality assurance system. It also holds the State-issued "Design Certificate for Pressure Vessel", "Production Certificate for Grade 1 & Grade 2 Pressure Vessel" and "Design & Production Certificate for Enamel & Glass Reaction Tank."

With the strong technical force and consummate experience about pharmaceutical machinery fabrication, this factory can manufacture the equipment for the injection, tablet, powder, medicinal material, chemical and water processing. The products of this factory are of 8 categories and over 100 specifications such as the medicinal packing machine, medicinal crushing machine, liquid filling machine, solid preparation machine, antibiotic separate packager, vacuum freeze drier, water distillator, water processor with ultrafiltration membrane, chemical equipment for medicinal material and customized equipment. This factory enjoys a high reputation because of its reliable product quality.

Factory

Add: 777, JinYe Road, Jlaoqiao Industrial Park, Shanghai, China.

P.Code: 201206 Fax: 86-21-58991994

Tel: 58995186(to different extensions) 58991992(direct line)

Sales Department

Add: 359, Kang Ding Road, Shanghai, China

P.Code: 200040 Fax: 86-21-62564739 Tel: 62170175.62580660

SHANGHALYUAN DONG PHARMAGEUTICAL MAGHINERY GENERAL FACTORY








FG150 Fluid-Bed Drier



ACSD1~20 Ampoule Washing, Drying, Filling & Sealing Combine

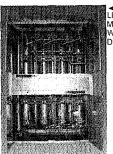


◀ LD Tubular Multi-Effect Water Distillator



GZL Series Vacuum Freeze Drier

SXZ20,SXZ50 Rotary Evaporator

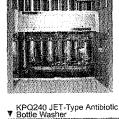



kMZ600/35 Tunnel -Type Penicillin Bottle Sterilizing Drier



KG**Ğ**300 Rolling Capper





Shanghai Tian he Pharm Aceutical Machinery





Director:zhao zhi nena

Shanghai Tianhe Pharmacy Machinery Factory

Address:Chengjia Town, Chongming County, Shanghai

Tel: 86-021-59431285
Fax: 86-021-59432024
Business Dept. Address:Room 706,31

Hengfeng Lu, Shanghai

Tel: 86-021-63811920 63813712

Zip code:200070

General Introduction, Shanghai Tianhe Pharmaceutical Machinery Factory(here inafater named as STPMF) being of national pharmaceutical machinery industry association and of Shanghai medicine industry association, is formed by reorgan izing the former Shanghai No.1 pharmaceutical machinery factory being of Shanghai light industry machinery stock company except its joint venture.

STPMF is not only the largest professional factory producing tablet press equipment, but also most rapidly develops new products of tablet press equipment. All tasks from casting to metal machining, and general assembly are aiconplished by itself. ZP 37 tablet press, SLB high speed auto. Alfoil packing machine and PZ10/55 table bottling machine etc. all meet GMP requirement, which are on behalf of tablet press product developing direction and win honour certificates for Shanghai excellent science techiology products. It has many service offices located in every privince and city of our country. Every custoner used TIANHE products can enjoy the best service after sale.

STPMF has obtained customer's trust and satisfaction depending of its good quality, complete varieties, good service and cheep price. It has a teshnology troop with international level and rich experence in manu facturing tablet press equipment. STPMF kindly prorides excellent services for the customers all over the world, and welcomes to make business.

DGF Elextromagnetic Induction Sealing Machine with Composite Aluminium-foil

- Application: The machine is used for sealing the plastic and the glass bottle full of medicine, food, cosmetics and chemical products, efficiently to prevent them from overflowing when being squeezed and getting damp.
 - 2. Features:
 - 2.1 Heating by the process of contactless electromagnetic induction is adopted for the machine, to have the aluminium-foil on the bottle-rim generate heat and then melt on the bottle-rim, and the sealing will be realized.
 - 2.2 The sealing rate of the machine is 100%. Even if the remainer on the bottle-rim, the bottled products can not be affected.
 - 2.3 Auto, counting display and auto, detecting unit for taking-off the bottle-over without aluminium-foil are used to assure the pass percentage 100%.
- 2.4 It is easy for the machine to be installed and operated. Its running is stable and without noise. When replacing the bottle(different sizes), it is easy to adjust the machine.
- 2.5 The electromagnetic radiation strength is up to Class I standard of state safety zone without injurying human body.
- Main Technical Data:
- 3.1 Working capacity:40~90 bottles per min.
- 3.2 Bottle sizes:Height 52~150mm(it can be more when specially required)Dia. of bottle mouth 25° 70mm(it can be more when specially required)
- 3.3 The time of induction heating:5~15 sec.
- 3.4 The total power of the machine:Type DGF-1-2kw single-phase 220V/15A Type DGF-2-3kw single-phase 220V/20A Type DGF-3-4kw 3- phase 220V/380V/15A(the power can be increased to 5~20kw when specially required)
- 3.5 Induction heating frequency:about 400KHz
- 3.6 The weight of a complete set of the machine:about 265kg(bottle transmitting unit=50kg,bottle arranging unit=35kg, induction generator=180kg)
- 3.7 Overall dimension, transmitting unit:1600 × 700 × 1500mm arranging unit:820 × 500 × 900mm induction generator:1000 × 600 × 1100mm



Shanghai Tian he PharmAceutical

® Machinery Factory

Address: Chengjia Town, Chongming County, Shanghai Tel: 86-021-59431285 Fax: 86-021-59432024 Business Dept. Address:Room 706,31 Hengfeng Lu, Shanghai Tel: 86-021-63811920 63813712 Zip code:200070

ZP33 ROTARY TABLET PRESS




This machine is designed for pressing round tablets from various kinds of granular materials. It is primary machine to fit batch production of tablets. The mechanism of the machine belongs to doublepress type. It consists of two sets of feed assemblage and press-roller unit. The turret carries 33 sets of punches and dies, so 66 tablets will be pressed for each revolition of the turret. During operation, the speed of the turret rotation, the depth of the filling material and the thickness of the tablets can be adjusted. The mechanical buffer unit will prevent punches and apparatus from overload damage. A powder suction unit is attached to take in the residual powder during the operation, so that the blockage by the residual powder material can be reused. This machine has won the Fine Quality Product Prize awarded by China National Medical Management Bureau.

Main Technical Specifications:

Dies 33 sets Max.Pressure 40kN Max.Dia of Tablet 12mm Max Depth of Fill 15mm Max Thickness of Tablet 6mm Turret Speed 11-28r/min Production Capacity 43.000-110.000pc/h 2.2KW960r/min380/50V/Hz Motor Overall Size $930 \times 900 \times 1600$ mm

Net Weight

SLB HIGH SPEED **AUTOMATIC ALUMINIUM FILM** PACKING MACHINE

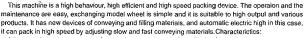

It has good structure, operaton is easy, has no noises and good behaviours.

With imported frequency conversion device, heating sealing speed can be adjusted with no step so that installation can be carried out at low sped, production can be carried out at high speed, and raw materials can be sowed.

It is automatically controlled, packing is accurate and beautiful.


Based on packing demand, cutting products packing of differents sizes can be adjusted. It needs no hands during working, and meets the G.M.P. requirement.

Use of aluminum film paper is suitable to pactablets, capsules and all kinds of blocks.

The volume is small, the capacity is high and it has many functions.

Conveying materials, filling, sealing, cutting are all automatically operated and with cooling wind pipe. It is air pressure type cutting, stable and noiseless.

| Main Technica | I Specifications: | | |
|--------------------|-----------------------|------------------------|------------------------|
| Model | SQY-APM-A | SQY-APM-B | SQY-APM-C |
| Capucity | 400-600pcs/min(4line) | 600-1200pcs/min(6line) | 800-1600pcs/min(8line) |
| Max.Wrapping Width | 160mm | 180mm | 280mm |
| Motor power | 0.4kw | 0.4kw | 0.4kw |
| Heater power | 3kw | 3kw | 3kw |
| No.of cutting | 1-5paekage | 1-5paekage | 1-5paekage |
| Machine size | 840 × 1050 × 1640(mm) | 840 × 1350 × 1640(mm) | 940 × 1450 × 1640(mm) |
| Packing size | 950 × 1200 × 1680(mm) | 100 × 1700 × 1900(mm) | 1050 × 1800 × 1900(mm) |
| Net weight | 370kg | 430kg | 510kg |
| Gross weight | 430kg | 510kg | 650kg |

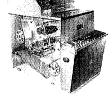

THP THP-1 THP-2 FLOWER BASKET TYPE TABLET **PRESS**

This machine is a medium continuous automatic tablet press which is one of the essential equipments for pressing granular raw materials into tablets or medicines taken with water in pharmaceutical/chemical. food, electronics industries. It can be used to press medicinal tablets, catalysts, candies, calcium tablets, coffee pills, powder metallurgical parts electronics components, various tablets of pesticide and fertilizer etc, in small-batch production, laboratory and hospital. The machine can press special-shaped ringlike tablets of different kinds, and can also process the tablets of pesticide both sides are impressed trademarks, letters and simple patterns.

| main recimical spe | CHICAHOHS: | |
|----------------------|--------------------|--------------------|
| THP | THP-1 | THP-2 |
| 40kn | 100km | 100km |
| 30mm | 40mm | 50mm |
| 25mm | 45mm | 45mm |
| 22mm | 35mm | 35mm |
| 5.10mm | 6-14mm | 6-14mm |
| 25mm | 45mm | 45mm |
| 34.50pcs/min | 35.46pcs/min | 46pcs/min |
| 1.1kw | 3.7kw | 5.5kw |
| 685 × 565 × 1610(mm) | 666 × 860 × 1824mm | 666 × 860 × 1824mm |
| 520ka | 750kg | 770ko |




Shanghai Tian he PharmAceutical

® Machinery Factory

Address: Chengila Town, Chongming County, Shanghai Tel: 86-021-59431285 Fax: 86-021-59432024 Business Dept. Address:Room 706.31 Hengfeng Lu, Shanghai Tel: 86-021-63811920 63813712 Zip code:200070

ZP19 ROTARY TABLET PRESS

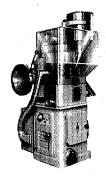

This machine is suitable for pressing round tablets from various kinds of granular materials, specially for small-scale batch production. The turret carries 19 sets of punches and dies, so 19 tablets will be pressed for each revolution of the turret. During operation, the speed of the turret rotation, the depth of the filling mateerial and the thickness of the tablets can be adjusted. The mechanical buffer unit will prevent punches and apparatus from overload damage. A powder suction unit is attached to take in the residual powder during operation, so that the blockage by the residual powder can be eliminated and the collected powder material can be reused.

This machine has won the Fine Quality Product Prize awarede by China National Medical Management


MAIN TECHNICAL SPECIFICATION:

Net Weight

Dies Max.Pressure 40KN Max.Dia.of Tablet 12mm 15mm Max.Depth of Fill Max.Thickness of Tablet 6mm 20-40r/min Turret Speed 25,000-45,000pc/h **Production Capacity** Motor 2.2kw960r/min380/50V/Hz Overall Size 880 × 750 × 1600mm

560kg

ZPW21A ROTARY TABLET PRESS

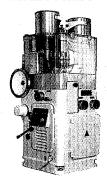

This machine is designed for pressing tablets from various kinds of granular materials. As its feature, the machine can press not only round tablets but also different geometrical shape tablets, double-layered and annular tablets, the tablets may have the impressed letters on both sides. The turret carries 21 sets of movable-punches and dies, so 21 tablets will be pressed for each revolution of the turret. During operation, the speed of the turret rotation, the depth of the filling material and the thickness of the tablets can be adjusted. The pressure for tablet pressing is set by the hydraulic drive system, a pressure gauge is attached to the system to show the pressure for tablet pressing is set by the hydraulic drive system, a pressure gauge is attached to the system to


show the pressure value. A overload protection unit is in cluded in the system to avoid the damage of the punches and apparatus. When overload occurs, the machine stops automa tically. A independent powder suitin unit is attached to the machine. It takes in the residual powder through its nozzles and hoses during operation. So that the blockage by clustered powder and residue can be eliminated and

collected powder material can be reused. This machine is developed recently in 1980's by our factory. Output sells abroad.

Main Technical Specifications:

Dies 21sets Max.Pressure ลดหพ Max.Dia.of Tablet 20mm Max Deoth of Fill 18mm Max. Thickness of Tablet 9mm Turret Speed 12-36r/min Production Capacity 15.000-45.000pc/h 3KW1440r/min380/50VHz Motor Over all Size 900 × 840 × 1640mm Net Weight

ZP25 ROTARY TABLET PRESS

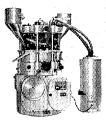

This machine is designed for pressing tablets from various kinds of granular materials. Besides round tablets the machine can press the tablets in different geometrial shape and impressed letters on both sides of the tables. The mechanism of the machine belongs to double-press type, including two sets of feed assemblage and press-roller unit. The turret carries 25 sets of movable-punches and dies, so 50 tablets are produced by each revolution of the turret rotation. When round tablets are pressed, there are four sets of punches and dies to choose:1-tablet mode of punch and die through 4-tablet mode of punch and die. The choice is dependent on the diameters of the tablets. So this machine is more suitable for mass production. When operated, the speed of the turret rotation, the depth of the filling material and the thickness of tablets are adjustable. A independent suction unit is attached to the machine. The residual powder during operation can be taken in by this suction unit through its nozzles and hoses. Thus the blockage by clustered powder is eliminated and the collected powder material can be reused.

Main Technical Specifications: Dies 25 sets Max.Pressure **60KM** Max.Dia.of Tablet 20mm Max Depth of Fill 15mm Max.Thickness of Tablet Turret Speed 6-10r/min Motor 3KW950r/nin380/50V/Hz Overall Size

1,800kg

Net Weight

1-tablet punch and die set 18,000- 30,000pc/h(up to 20mm Dia.) 2-tablet punch and die set 36,000-60,000pc/h(up to 12mm Dia.) 3-tablet punch and die set 54,000-90,000pc/h(up to 9mm Dia.) 1100 × 1100 × 1700mm 4-tablet punch and die set 72,000-120,000pc/h(up to 6mm Dia.)

Production Capacity:




Shanghai Tian he PharmAceutical

® Machinery Factory

Address: Chengjia Town, Chongming County, Shanghai Tel: 86-021-59431285 Fax: 86-021-59432024 Business Dept. Address:Room 706,31 Hengfeng Lu, Shanghai Tel: 86-021-63811920 63813712 Zip code:200070

LT1400 Vertical **Labelling Machine**

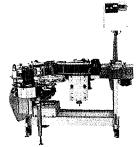

Main Technical Specifications:

1.Speed of conveyer belt:7.5M/min

2.Range of labelling: {(+) 20mm- +> 70mm full circumference}~+> 150mm}

3.Range of working table's height:(750-830)mm

4.Adjusting range of electric elevator:(0-100)mm 5.Adjusting range of webbing apparatus:(0-100)mm

6. Width of conveyer belt: 190mm

7.Specification:L × 8 × H 1.4m × 1.2m × 1.5m

8.Voltage,power;220v,250w

9.Labelling precision: ± 0.3mm(left-right direction) ± 0.8mm(turning direction)

10.Working ability:15-200bottles/min

11.Label:special webbing sticker

12.External diameter of webbing sticker:300mm(near 240m long)

13.Internal diameter of reel:76mm

V-TYPE EFFICIENT

MIXER




V-TYPE EFFICIENT MIXER This Mixer has the following distinguishing features, no blind angle in the cylinders, no accumulation of materials, high speed, and short mixing time. With the body made of stainless steel and both inner and outer walls polished, the cylinders have elegant appearance. The charging of materials can be made either by manual operation or through vacuum suction. For discharging, rotating valves and rapid swivel heads are adopted.

This Maxer is applicable for mixing powdery or granular materials in pharmaceutical and other industries. It consists of two v-type cylinders of different attitudes, and through mechanical transmission the materials in these cylinders are made to move back and forht and turn over again and again, thereby achieving the goal of uniform mixing.

Main Technical Specifications:

| mon secin | mear specin | cations, | | |
|--------------------|--------------|-----------|-----------------------|----------------|
| Specification | Capacity(kg) | Power(kw) | Overalldimensions(mm) | Net Weight(kg) |
| 0.05m ³ | 12.5 | 0.37 | 1150 × 800 × 1200 | 200 |
| 0.15m ³ | 37.5 | 3.0 | 1600 × 1100 × 1700 | 250 |
| 0.3m ³ | 75 | 1.5 | 2100 × 1600 × 2400 | 480 |
| 0.5m ³ | 125 | 1.5 | 2200 × 2000 × 2700 | 500 |
| 1m ³ | 250 | 2.2 | 3100 × 2300 × 3100 | 700 |
| 1.5m ³ | 375 | 3 | 3300 × 2600 × 3550 | 900 |
| 3m3 | 750 | 4 | 4200 × 2700 × 3600 | 1100 |

ZP37 ROTARY TABLET PRESS

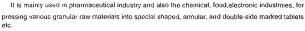

Characteristics:

The machine is a double-press type, continuous automatic lablet press for pressing granular raw materials

The outside part of the machine is fully enclosed. The material of the cover and the inside table face is stainless steel. The turret face is specially disposed keeping the lustne of the face and prevent the cross pollution, meeting the G.M.P. requirement.

The speed adjatment of the driver is controlled by the frequeny conversion device.

It has transparent windows so that press condition can be observed clearly and the window can be

opened cleaning and maintenance is easier.

The lay-out of all the controllers and the operate parts is reasonable.

A overload protection unit is included in the system to avoid the damage of the punches and apparatus, when overload occurs, the machine stops auto maticaly.

The hydraulie system can assure the stability

The punches and dies are-also applicable to ZP19, ZP33, ZP33B, and ZP35 type.

Main Technical Specifications:

37sets

| - 1 | vlax. Pressure | 80KN | | | |
|-----|---|-------------------------------|--|--|--|
| 1 | Vlax.Dia of tablet | titer-no-sums tablet) | | | |
| 1 | Max.Depth of fill | térentapopal snapirlong side) | | | |
| | ax.Depth of till ax.Thickness of tablet irret speed | 15mm | | | |
| | | 6mm | | | |
| | | 12-40r/min | | | |
| | Max.Capacity | 170000pcs/h | | | |
| | Viotor | 4kw,1030r/min.380v/50Hz | | | |
| • | Overall size | 1170 × 1150 × 1750mm | | | |
| 1 | Net weight | 2300ka | | | |



Company brief:

Yong he machinery CO.LTD is the member of China medicine equipment union. Our company is an important enterprise which is engaged in designing and manufacturing wrapping machines and medicine machines. We have strong technology power and enough equipment. The technology of productions is advanced, the construction is reasonable and we have met the normal requirement of GMP medicine producing.


Our company is located in Hairnen city jiiangsu province. The distance between Shanghai and our company is only 60 kilometies. The traffic is convenient and surroundings is beautiful. The companity holds the land of 22650 square meters. Main products: (1) Vials power filling and wrapping production line (2) Vials liquid (frozen biological products) filling production line (3) easy-opening cap (steal-protecting cover) oral liquid production line. Meanwhile we can design and manufacture abnormal products for customers. The quality of products is credible and the price is reasonable. Our good service behind selling is also praised by customers.

Our company would like to cooperate and develope with customers home and abroad.


> Legal representive: Yu yong-ping. Add: Haimen cotton plant Jiangsu province (North of haimen narbour)

> > P.C.: 226114


Fax: (0513)-2292308 Handy phone: (0513)-9622697, 01396279890

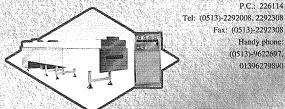

Model MSH-800 High Temperature Sterillizing Tunnel

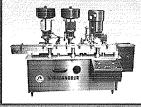


ZFJ-120 Straight Line Power Metering and Filling Machine for Vials

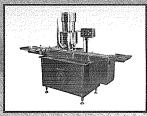

Model FGL Multiple Function Capping Machine



CXP-5-10 Lying Supersovmdwave Bottle Cleaning Machine

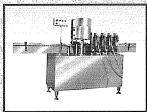


GCS Series Linear Frozen Product Liquid Filling Machine



ZYG-80 Easyopening cap Liquid Filling Machine



Fzl2-120 Computer-controlled Screw Type Powder Bottling Machine

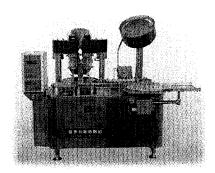

Main Technical Data

Production Capacity:60-120 bottles/min(stepless speed regulation)

Antibiolic Bottle Size:5-20ml(with some parts changed)

Conforming Rate of Capping:99.7%

Powder Charging Error:no more than +/-4%

Powder Charging Containment: 0.1-4 gram

Powder of Corollary Motor:0.8KW

Outline Dimentions:800 × 1200 × 1550mm

The main products include complete sets of powder-injection equipment. The factory was established in 1971, with more than twenty years of experience in manufacturing pharmaceutical machines. As one of the major manufacturers of powder-injection equipment in China, it owns rich technical resources and advanced processing equipment. Its products have reached advanced level in the domestic market. We enjoy high reputation for our management, which features "Credit First, Customers Paramount".



Title of Enterprise:Shandong Penglai Pharmaceutical Machinery Factory Address:Penglai City, Shandong Province, China Post Code:265609 Phone:0086-535-5971223 Fax:0086-535-5971223



400, 600, 800. 系列全自动 胶囊充填机



CFM •••• 400.600.800. Serial automatic capsule filling machine.

High filling percentage and dosing is the rule for design of CFMserial filling machine.

Small volume .low energy cost. high efficacy. automatic-standard and serial is the best consideration for the customer.




Jin Wan Co. LTD

ADDRESS: 559 HAMATANG RD. DANDONG LIADNING P. R. CHINA TELEPHONE: 86 – 415 – 4151322 FAX: 86 – 415 – 4151380

FAX: 86 - 415 - 4151380 POSTCODE: 118003

E - mail: jiayi @ public2. ddptt. ln. cn.

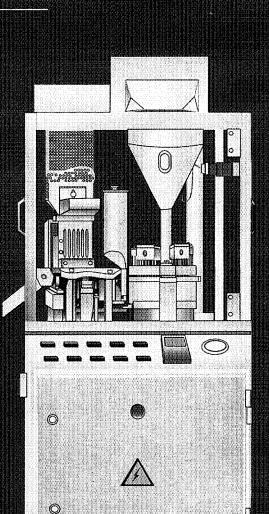



COMMODITY: FLUID - BED CYCLONE CATCHING DRYER MODEL: FGX - 120

FEPM

Application:

It is used in pharmaceutical, food and chemical industries for drying powder and granular materials.

Features:

Based on the principle of fluid – bed dryer, cyclone separator is adopted in this machine to catch powder, where, powder trapped by the first stage cyclone separator returns to the main unit through vortex gas – seal valve so as to avoid problem of simultaneously drying in a batch the wet powder tapped in early stage of operation and overcome the trouble of cleaning powder in the filtering bag of fluid – bed dryer during operation and shortcoming of fiber shedding.

Technical data:

Productivity(kg/time): 80-160Air - in temperature(°C): 20-120Moisture in finished product(%): 0.5-3

Obtainable ratio(%): 99
Pressure of gas source(Mpa): 0.6

Compressed air consumption(m³/min): 0.6

Power(kw): 11

Overall Dimensions(mm):


1800x 1800x 2800 Other Specification:

FGX3, FGX10, FGX30, FGX60, FGX300, FGX1000



400039 Tel: China

Fax: (023)68604229 (023)68604229

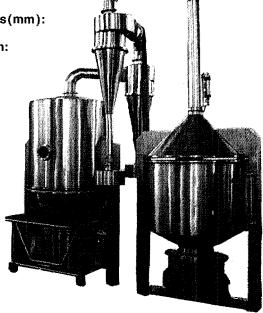


CHONGQING FINE AND ELABORATE PHARMACEUTICAL MACHINERY COMPANY, LTD.

HARBIN AIRCRAFT MANUFACTURING CORPORATION

Founded in 1952, it is a large backbone enterprise of AVIC. Besides aircrafts and helicopters, it has developed a large variety of civil products, pharmaceutical machinery is also included.

I:Hard Gelatin Capsule
Filling Machine
YJF15

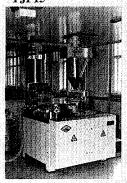

III: Horizontal Breakable Neck Ampoule Maker WABI-2



II: Glass Tube Drawing
Production Line
BLGX10



IV: Antibiotic Vial Antibiotic Powder

Automatic Production Line

HFZ 300




ADD:15 Youxie Street, Pingfang District, Harbin, Heilongjiang, 150060, China

TEL:86-451-6501726 FAX:86-451-6502223



GUANGDONG PROV.

HUIYANG MACHINERY FACTORY

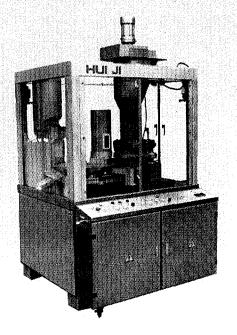

AUTOMATIC CAPSULE FILLING MACHINE SERIES

The Automatic Capsule Filling Machine produced by our factory, have speciality such as: controlled by computers, having central rotary unit, intermittently rotating around the turning table. It can automatically finish: feeding capsules, separating capsules, filling medicine, rejecting faulty capsules, closing capsules, ejecting full capsules and cleaning die. The advantages are reliable in capsule's right position, adjustable in filling weight, and accurate in dosage. The error of weight dosage less than ±5%. The rate of finished products is more than 98%. These machines completely corresponds to the pharmaceutical G. M. P. standards due to them use sealable device.

TECHNICAL DATA

| Туре | ZJT-20A | ZJT-40 | ZJT-40A | NJP-500 | NJP-700 |
|------------------------|--|-----------------------------|-----------------------------|-------------------------------|-----------------------------|
| Output Capacity | 20000 CAPSULES/H | 48000 CAPSULES/H | 56000 CAPSULES/H | 500 CAPSULESMIN | 700 CAPSULES/MIN |
| Capsules to be used | 00"~ 4" | 00*~ 4* 00*~ 4* | | 0*~ 4* | 0*~ 4* |
| Types of Fillings | POWDER GRANULE PELLET | POWDER GRANULE PELLET | POWDER GRANULE PELLET | POWDER GRANULE PELLET | POWDER GRANULE PELLET |
| Power | 3.88kW | 4.35kW | 4.35kW | 3.4kW | 3.4kW |
| Total Weight | 997kg | 1597kg | 1597kg | 750kg | 750kg |
| Dimension | Dimension 1015x1000x1850 1234x1200x2 LxWxH LxWxH | | 1234x1200x2268 LxWxH | 784x930x2260 LxWx H | 784x930x2260 LxWxH |
| Accessories | Vacuum pump 1 set, Cleaner 1 set | | | | |




GUANGDONG PROV.

HUTYANG MACHINERY FACTORY

GUANGDONG PROV. HUIYANG MACHINERY FACTORY was set up in 1956. We have all the facitilies at hand for machinery buliding. After ten years of hard work, now, we have been able to supply hard capsule filling machine series, traditional Chinese medicine pill machine, super micro grain pulverizer and program-control rubber stopper washer, etc. more than twenty kinds of products. Many of them have obtained the Scientific and Technoglogical Progress Reward from the Ministry of Machinery.

Our products have been exported to USA. Canada, Ukraine, Turkey, Pakistan and Southeast Asian countries, we have set up the agency at USA.

We welcome your inquiry and cooperation for a greater expectation.

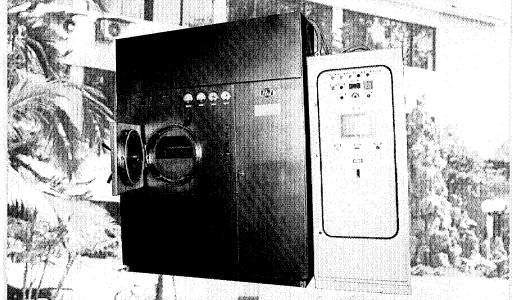

Address: 26. Xia jiao middle road, Hui Zhou

Guang Dong Prov. P . R . China

Tel: 86-752-2227000

Fax: 86-752-2241409

P.C: 516001



Shanghai Ying Hwa Pharmacy Machine Industrial Company

An Excellent Science-Technique Enterprise

offer you the novel series products of the GNP powder and ampoule assembly line.

Since the founding of our company in Nov of 1990, it has successively set up Ying Hwa pharmacy machine research bureau and Ying Hwa pharmacy installation works. It possesses higher technique strength and speciliges in the research and manufacture of pharmacy powder and ampoule production installation. It is a productive real body concentrated technique, industry and trade.

The main prodelectronic numerical screw dividing ashonourably won Price" awarded by "Nation Pharmacy ucts Exhibition;"PSR washing machine; channel oven, from minute to 800 amtomatic tagging maand ampoule as-



control double head sembly machine, "Golden Sheep the first session Machine New Prod-1-2 vertical voil YHF steriliging 300 ampoules/poules/minute; auchine; GNP Powder sembly line novel


ucts are:PLF1-4

series products; all kinds of conveger set series. The quality of the products attains to the standard.


Service is best, upto now it has made 19 assembly lines for each famous group and pharmacy works. Besides selling well the products 188 sets to 25 provinces and cities, it exports them far to Indonisia, Pakistan, Thailand and other southeaster countries.

Our company has been appraised as "Excellent science Technique Enterprise" by Shanghai municipality and yangpu district science technique committee in successive many years since 1991. Zhang yue ying, legal person representative, chaiman of board of company and higher engineer, has won the title "Nation Advanced Producer and Excellent Enterpriser" many times.

Zhang Yue Ying, legal person representative and chairman of board of the company, and general manager Xu Bao Min welcome warmly broad masses of cuslomfers and people of all ranks to give advice.

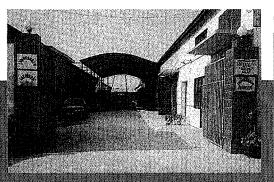



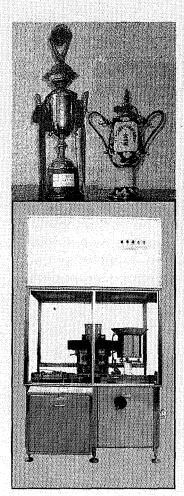

Company Address: No 8 Passage 510 San Men Lu, Shanghai Fax: (021) 65444114 Tel: (021) 65366089 Post Code: 200434

Works Address: Shanghai 318 National Highway QingPu County ZhaoHang own Xin Chen Economical District No 4—No 8

Tel: 021-59753538

Post Code: 201703

(1)Characteristics:

1. The filling machine is fully computerized.

2.The machine applies Man-Machine talk, screen shows.

3. Eight menus for easy selection.

4.No filling while no vial on the filling station.

5.Easy to set and amend the parameter of the stepmotor, and chosse

to add opposite roatation if you need.

6.If the auger touches the surface of the funnel the computer will show

A or B auger is malfunction and alarms.

(2)Applicable Range:

7 C.C.(25mm,Dim)mould vials,(22mm,Dim)rubing aials

10 C.C.(28mm,Dim)mould vials,100

C.C.(53mm,Dim)mould vials

(3)Capacity:

0.06-2.0gng/vial-120vial/min

2.0-3.0gng/vial-70-80vial/min

5.0gng/vial-30-40vial/min

(4)Tolerance:

0.5-2.0gm---±2-4%

150mg---±3-6%

(5)Powder Range:

Crystalline powder, spray-drying powder, freezedrying powder and

powder with poor fluidity all can be filled.

(6)Dimension:L \times W \times H 1450 \times 900 \times 1700mm.


(7)Driving pattern: Applying numerical controlling stepmotor.

(8)Stoppering:Stopper is picked up by vacuum and then insert to the

vial mouth, Successful rate of stoppering operation 99.9%

(9)Power:

2KW,380V,50Hz (10)Weight:650Kg.

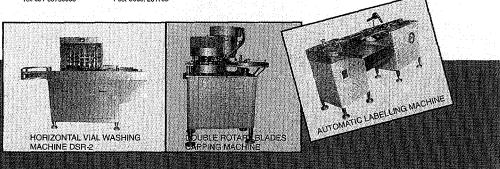


CHONGQING GUANGXIA DRYING EQUIPMENT & ENGINEERING CO.

GENERAL INTRODUCTION:

It's a member of China Association for Pharmaceutical Equipment, a first grade school run enterprise of Chongqing. The company specialized in reaserching and development of drying and granulating equipment.

For the coming century we will keep on winning our user's confidence by our powerful talent system, distinguished quality and satisfying service. Our products have been widely sold to more than twenty provinces in our country and also exported to Pakistan, Malaysia, Indonesia and Tazania with high praise from users. With our effort and widely cooperation with more friends, we wish to be more helpful to parmaceutical industry, particularly to the solving of the manufacture problem of the form of medicine.

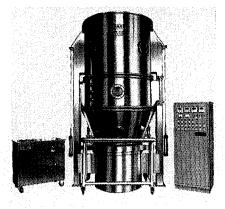

MAIN PRODUCTS:


Model PGL spraying dryer granulator Model FL-A fluid-bed dryer granulator Model FG-A fluid-bed dryer

Legal person: Mr. Xu Zemin

ADD:128# Shuangfeng Road, Yubei Area Chongqing China

TEL:86-23-67823244,67822860 FAX:86-23-67822860

TELEX:0633 POST CODE:401120

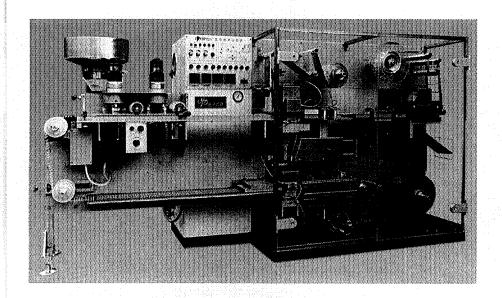

Beijing Parmaceutical Mechanization Equipment Factory is a appointed factory by State Pharmaceutical Administration of China. It has history more than twenty years and has more than fifty technicians. The main products as follow: Plate Aluminium-Plastic Packing Machine Type DPP170, Plate Aluminium-Aluminium Packing Machine Type DPP170-S, Plate Aluminiun-Plastic Packing Machine Type DPT250, Plastic Tube Ointment Machine Type ZGC25. Our factory is also appoint factory who can produce I, II container with pressure used in pharmaceutical and chemical industry.

BEIJING PHARMACEUTICAL MACHINE & EQUIPMENT FACTORY

Address:Dajiaoting Guang Qu Street, Beijing China Tel: (010)67716779 67722573 Fax: (010)67721941

Cable: 0009

Post: 100022


Contactor: Li Qi Yun Wang You Jian






ENTERPRISE BRIEF INTRODUCTION

Our enterprise is a member of Pharmacy Industry Equipment Association of China specialized in manufacturing pharmaceutical equipment. We have supplied one thousand blister packaging and filling equipments for customers of country and overseas. It is well-received by the user for the superb quality and good service after-sale. The main products have gained the prizes from nation and concerning departments, and have also passed quality certification by technology inspection centre of National Medicine Administration.

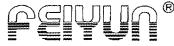



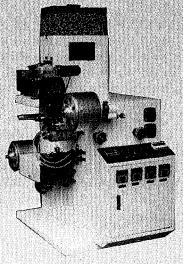

ZGL-160 Automatic High-Speed Double Aluminium Foil Packer


Main Application: Applied to double aluminium packing of tablet and capsule, etc. in pharmaceutical industry.

Operating Principle and Structural Characteristics: Adopting frequency control of motor speed, punching the array feeding material into required plates after it has been heat sealed by the moulding wheel, and has automatic temperature regulation, automatic counting for finished product and other functions.





DPT-140 Automatic Aluminium-Plastics Packer

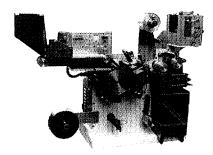


Main Application:Used for packing tablet, capsule, suppository and medical utensil in pharmaceutical industry.


Operating Principle and Structrual Characteristics:PVC after having been heated will form bubble cap under negative pressure, and after the medicine material is fed in, the bubble cap will be heat sealed through knurled roll and aluminium foil, and then it will be punched into finished plates.

NCB-200B Semi-Automatic Capsule Filling Machine

Mina application:Used for producing hard capsules in the pharmaceutical industry.

Operating Principle and Structural Charcteristics: Fulfilling the capsulearranging, setting, separating, filling, locking and other processes of capsules through mechanical drive.

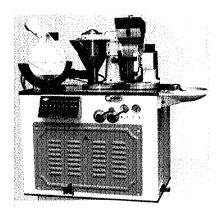

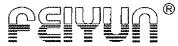

ADD:Dongshan Development Zone, Ruian, Zhejiang, P.R.China (Opposite to Dongshan Post Office)
Tel:0577-5603082 5603317 5603081 5603316 5603315
Cable:7845 Fax:0577-5603083 P.C.:325200

JINZHOU OLDTREES PACKING MACHINERY CO.,LTD.

MODEL DPB210 BLISTER PACKING MACHINE



JINZHOU OLDTREES PACKING MACHINERY CO.,LTD.is a joint venture for manufacturing medicine and food packing machines, medical PTP aluminum packing materials.


Main Products:

DPA250,DPA130,DPB210,DPT220,DPT130 Blister Packing Machine;


DQ325,DB420 Multifunction Packing Machines:

DSPI,DSPII Plastic Ampoule Forming and Filling-Sealing Packing Machines; DXDK200,DXDK320,DXDKD660 Shading Packing Machines.

MODEL DPT130 BLISTER PACKING MACHINE

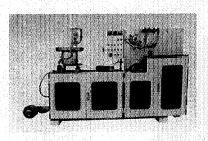

MODEL DPA250 BLISTER PACKING MACHINE

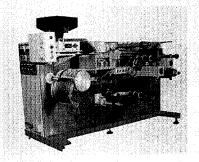

GENERAL MANEGER: SHI LIZHI

ADDRESS: NO.63 CHONGQING ROAD, GUTA DISTRICT JINZHOU CITY, LIAONING PROVINCE CHINA
TEL: 86-0416-4167750 4165752 CABLE: 2623 FAX: 86-0416-4168359 POST CODE: 121001

JIANGNAN PACKING MACHINERY CO.,LTD. RUIAN ZHEJIANG

Our company, a member under the Chinese Trade Association of pharmaceutic Equipment, which passed the ministerial level acceptance check of total quality control standard in 1993, is accredited as the Model Enterprise of Zhejiang Sparkling program, patent Top loo Enterprise of Zhejiang Province, Provincial- level Green Advanced Enterprise, Key Enterprise of Wenzhou City and Top lo Enterprise of Ruian City etc.


The corporation is producing annually over 300 blister packing machines, which have found a ready maket in all over the country (Hong kong, Taiwan), Europe, Middle-East Area, South-A frica and whole Southeast Asea. It has over 1000 domestic and foreign users. Its quantity on sale has achieved the first among the product of same-type over the world.

Series BG (150 · 80 · 40)

AUTOMATIC TABLET COATING


96 National Grade New Product




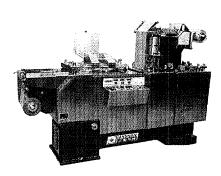




YH-420 Paper Box Marking Machine

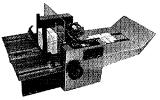

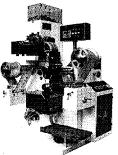

DSL-160B Full Automatic Doubleside Aluminium- foil Heat-sealing Packing Machine



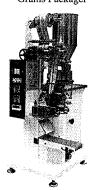

地址: 浙江省瑞安市紅旗工業區 電話(Tel):5662604 5663701 電挂:(Cable):7738 傳真(Fax):5663983 Add:Hongqi Industrial Area Ruian Zhejiang 手機(M.T): 1386802121

general manager: YIJUN YANG


HEILONGJIANG DEER PHARMACEUTICAL MACHINERY CO.,LTD



Automatic Pill Making Machine ZW Series:processing capacity from 20kg~120kg/h

Established in 1970, Heilongjiang Deer Pharmaceutical Machinery Co.,Ltd is specialized in manufacturing Chinese traditional medicine processing machinery.

Owing to our dedication and hard work, we have become one of the biggest suppliers in this field, our machines see a ready market in mainland China, Hongkong & Taiwan and exports to Southeast Asia countries. Your satisfaction to our machine is what we pursue.

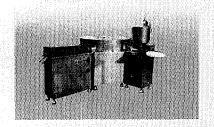

Suppository Drug Filling
Machine Sets
Packing strip for suppository



Plastic Bottle Filling, Counting & Sealing Assembly Line Consist of :Strenghtening, Sealing, Cap Screwing, Selfadhesive Labeling Machines

COMPANY: HEILONGJIANG DEER PHARMACEUTICAL MACHINERY CO.,LTD

ADD:No.37,Linhai Road,Jiamusi,Heilongjiang Province, 154007.China

TEL:0086-454-8782182 ATN:Mr.Gao Chuan FAX:0086-454-8782365



BRIEF INTRODUCTION

OF

SHANGHAI XIN-XING INDUSTRIAL AN TRADIMG MACHINERG CO.LTD SHANGHAI XIN-XING PHARMACEUTICAL MACHINERE PLANT

Our factory is a modern packing machinery plant with strong technical capacity. More than 90% of the engineering staff come from pharmaceutical, cosmetics and food factories. So the designed and manufactured auto packing machines are especially suitable for the specialised production of hospitals, pharmaceutical, cosmetics and food factories.

We are good at designing and manufacturing different types of auto liquid filler, oral medicament bottle sealer, auto ointment filler and the automatic machines for special packing materials. We are able to perform pesticide filling, plugging, inner revolving plug, sealing assembly line and large infusion combination line. We also design and install purifying shop and other equipment.

We undertake various work of non-standard designing, processing and manufacturing.

Our feature is to solve the customers' problems and to serve the customers with heart and soul.


Zhu Shi-hai General Manager Factory Director




Add; Roome • 501' House • 8' Lane 561 Yan Chang Rd • (Middle) Shanghai China P.Gode: 200072

P.Code; 200072 Fax: 86-21-56770079 Tel: 86-21-56625768 86-21-56770079



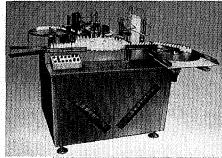

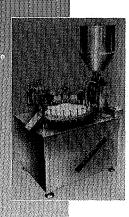

Description:

The machine employs computer-manipulating pneumatic equipment, suitable

for the liquid filling of plastic bottles 5-30ml. It has the function of filling while there are bottles, pausing while no bottles and there is no leaking. Filling brecision ≥ 99%. Output is 60 bottles per minute, quality rate ≥ 98%. There is also the equipment of not plugging and holding the outer cock while no bottles. The machine can form an assembly line with the label sticker.

Packing material requirements:

Inner plug edge > 3.0m/m, discernable of weight to facilitate oscillation. In order to improve the end product rate, the bottle neck of the inner plug and the tolerance with the mouth should be > 1-1.5m/m.

The outer cock should not have burr. It is to prevent the obstruction to actions.

Power:2kw

Voltage:38OV, 50H

Size:1.5M*1.3M*1.4M

Weight:520kg

Compressed air pressure: >4kg/cm(disposed by the customer)





ZHONGNAN NO. 1 PHARMACEUTICAL MACHINERY

ZTZT20/1000 DRUM LABEL-STICKING MACHINE



AUTOMATIC PACKAGING LINES FOR PLASTIC BOTTLES




ZHONG NAN PHARMACEUTICAL MACHINERY FACTORY NO.1 Add:QIYANG HUNAN CHINA P.C.:426100 Cable:2894 Tel:(0746)3222388 3223738 3223608 Fax:(0746)3222617 Zhongnan No. 1 pharmaceutical Machinery Factory (formerly: Zhongnan Pharmaceutical Machinery Factory) is one of the major production bases of pharmaceutical machinery in China, and the key backbone enterprise under the SPAC. It has over 1200 employees, and over 100 of them are technicians. The works has established a pharmaceutical machinery institute with senior engineers as its main staff. It now manufactures over 60 varieties of products in three types, i.e.pharmaceutical machinery, packaging machines and pressure vessels. Pharmaceutical machinery includes mixing machine, granulator, tablet press and sugarcoating machine and bubble-dispelling filling machine, etc. Packaging machines cover automatic production line for filling tablets and pills into bottles, automatic plastic bottle-filling production line, automatic filling production line for liquids, compact produciton line for large volume parenteral, blister packaging machine, and table type filling machine and bubble-dispelling filling machine, etc. Pressure vessel range contains multifunction extraction tank, Chinese medicine extraction tank (by means of filling water or alcohol), three-function evaporator and other pressure containers and equipment, as well as various nonstandard pressure vessels.

Encouraged by the policy of reform & open-door and the idea of developing science and technology to make our country strong, Shanghai Shan Li Pharmaceutical Equipment Technical corp. (herein after referred to as SL) was born in 1992. It consist of Ms. Wang Guan (Legal representative) who used to be the director and the senior engineer of Shanghai Pharmaceutical Equipment Corp., specialist, professor and senior technicians. All of them have rich experience in pharmaceutical, chemical and elector-mechanical industry form colleges, instituting and enterprises.

The purpose of SL is: Quality is First, Reputation is important.

The series of Electromagnetic Induction Sealing Cap is designed and

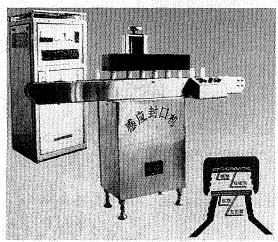

manufactured by SL. According to several testing, it is an excellent product with superior function, safety and reliability. This serials has received many prize. Especially, Model II, one product of this series, has many excellent characteristic, such as advantage tectonic, small shape, light

weight, simple and artistic and easily carried. It has caught the foreign advanced level.

上海三利制藥裝備技術公司

SHANGHAI SANLI PHARMACEUTICAL EQUIPMENT TECHNICAL CORPORATION


地址:中國上海延安西路900號(江蘇路口) Add: 900 Yan An Rd., Shanghai China 電話: (086)021-62131812 Tel: (086)021-62131812 電柱: 7785 Cable: 7785

郵編:200052 Zip: 200052



NUZHOE INDUSTRIAL PARK. AS MACHINERY & ELECTRONIC APPLIANCE CO., LTD.


Suzhou Industrial Park Hua Yuan machinery Co., Ltd. is a joint venture established Park. It is a high - tech enterprise mainly combining machinery electronics and optical instruments. It is a specialized company in researching, developing and producing capsule & tablet printing machinery, the products of AP - axial, DP round & axial, RP directional capsule printing machines have been applied for three patent rights from China National Patent Bureau. Our products Are sold very well all over China. Furthermore, they are exported to the Southeast Asia.

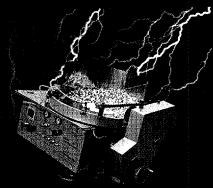


On abjects are Developing nationality quality-greating word - class level, footholding in Suzhou and moving towards the world. Our aim is: best quality and good service for customers. Our promises are good quality, all specifications of now products and best service.

We'd fike all friends to join us and bring about a good furnite.


ADD:No. & Nan Yuan Road. Suzhou ZIP:215006 TEL:0512-5308737 FAX:0512-5107028

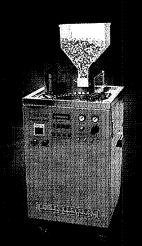



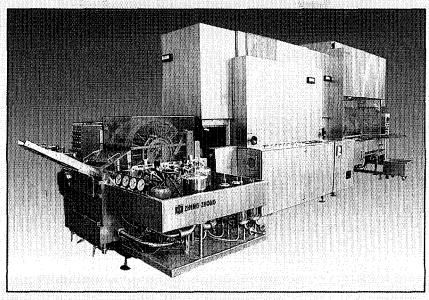

BXSZ1/20 LINKAGE SET FOR AMPOULE CLEANING STERILIZING FILLING AND SEALING

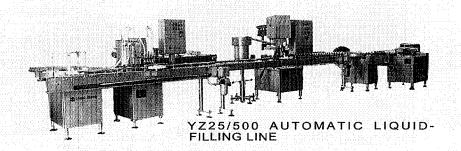

ZHONGNAN PHARMACEUTICAL MACHINERY FACTORY NO.2

Add:Xianjia Hu Road Changsha Hunan China

Tel:(0731)8809343 8809250 Fax: (0731)8809250 Post code: 410006 Factory director:ZHANG XIYAN

BRIEF INTRDUCTION OF ENTERPRISE

Aiming at offering excellent pharmacy for the health of mankind, according to standard of pharmaceutical quality(GMP) issued by state our factory rsearches, develops and produces pharmaceutical equipment with internation level, which sell well at home and abroad owing to the advaned process technology stable performance and exellent quality.

Our main products are:


BGB series high efficient sugarcoating equipment;

HD series multidirection movement mixer;

HL series material bucket mixing machine;

FZ series graining and granulating..ect


Our products have high marking share over the country and far into the markets in Indonesia, Pakistan, Philippine, Vie Nam, Cambodia, Hong King...

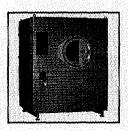

BGB series



HD series

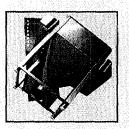


HL series

WENZHOU PHARMACEUTICAL EQUIPMENT FACTORY

Director: WangXiaolun
Add:Shanhou lane Liming Road
Wenzhou city,china
TEL:0577-8317378 8378153 9004282

Fax:0577-8362420 Zipcode: 325003

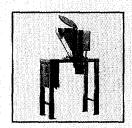

FZ series


DASONG ENTERPRISE GROUP

DASSNG




With rich experience in the production of machine and numbers of advanced engineers and technical personnel Beijing Dasong Machinery and Electronics Co., a high-new technology enterprise in the New Technology Development and Experiment Zone of Beijing Municipality, works mainly on the development, production and management of various packing machine, printing machine and relating consumers' goods.

Among about thirty types of the packing machine in over ten variety, SB series small sized packing machine was entitled as national-grade excellent new product in 1992 and got the measure licence issued by National Measure Bureau, as well as the DXD1000B automatic packing machine was appraised as national-grade torch item. For its wonderful reputation, Dasong packing machine has been sold through out our country and exported to more than ten foreign countries and regions.

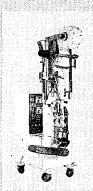


SB-80 GRAIN AUTOMATIC PACKAGING MACHINE

> SB-60F POWDER AUTOMATIC PACKAGING MACHINE

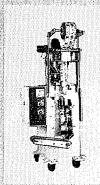

SB-60PJ
TABLET
AUTOMATIC
PACKAGING
MACHINE

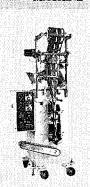

SB series automatic packaging machines, composing of powder, tablet, grain and liquid varieties, are the latest small-sized and vertical bag-packaging machines with the intermittent gearing and the three sides sealing, which were developed by Beijing Dosong Machinery & Electronics Co. It can automatically finish bag-making, filling, sealing and cutting.


Beijing Dasong Machinery & Electanics Co. Add: No. 118 Yimuyuan Haidian, Beijing(100091) Tel: (010)62882351 62881351 Fax: (010) 62881498

Liaison department: Packaging Machine Dept. ZhangXiaoPing



WENZHOU FANGZHENG ENTERPRISE GROUP(WENZHOU PACKAGNG MACHNERY GENERAL FACTORY)

Automatic Strapping Machine

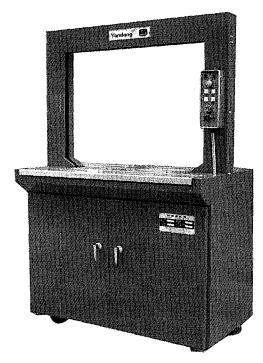

SK series Automatic Strapping Machine has been used for the strapping of packages in a wide range of all kinds of fields. It can be operated separately or combined into packaging line.

ADD:33 Feixia South Rb, Wenzhou, Zhejiang, China.

(TEL):0577-8631349 8635981

(Cable):0001

(Fax):0086-577-8632224

Postal Code:325003



Brief Introduction

Jiangsu GUI BAO Group is located at ZHUTANG Town of Jiangyin, famous for its beautiful scenary on the south bank of Yangtze River. It lies with Shanghai to its east, Suzhou and Wuxi to its south, and with Yangtze River Bridge and Shanghai-Nanjing Express Way called "Golden Passageway" nearby. It has a unique geographical advantages with convenient transportation both on land and by water. So, it is a best place for investors to settle in.

This group company has 185 million RMB yuan of fixed assets. It covers 181500 square meters of land space. And it has employed more than 1580 staff and workers, among them, 283 are senior and middle grade scientific and technical personnel. So, it possesses a strong scientific and technical force and it has owned advanced manufacturing and processing machines and a complete set of testing and inspecting equipments.

Its main products are medicines, food, chemical machinery, Chinese medicine machinery, plastic soft-packing machinery of which, the zs series of rotating-vibrative screens, GFSI series of high-efficiency crushers. WFI series of micro-crushers, GK-70 Type dry Grain-formers, Powder-Engineer Series' equipments, Finished Chinese Medicine Cutting gears, sugar-coating machines, twin-screwstem extruders, etc., which adopt sophisticated technology of the world. All the above-mentioned products won the "Golden Sheep Prize" from National Pharmacy Mnagement Bureau, the "Excellent Product Prize" from the Chinese Medicine Management Bureau, "Golden Ox Prize" from Jiangsu province, ane Gold and Silver Prize of state Star-Light Planning, It also won the Excellent Prize for new-and-advanced technique and newly created products etc., "GUI BAO" has been ranked as the key trunk enterprise by the State Pharmacy Management Bureau and Chinese Medicine Management Bureau, and evaluated as Jiangsu provincial New and Advanced Enterprise and a key close and congregated Enterprise for "Star-Light Tech.". It has been evaluated as an advanced or a first-class enterprise for its good reputation and always-obeying contracts for 8 yesrs in succession by WUXI and Jiangyin governments. Also, it has gained AAA Certificate from Jiangsu provincial System. It has been evaluated as Wuxi City Star Enterprise for 3 successive years and as Jiangyin City Star Enterprise for 5 years in succession.

GUI BAO products have enjoyed good reputation, possesses an important position in market and being popolar with customers both at home and abroad.

Mr Zhou Yuan Long, Chairman of Jiangsu GUI BAO Group, and Mr Kong Pei Du, the General Manager, on behalf of all the staff members and workers, deeply appreciate all their customers, old and new, for their sincere cooperation and support, devote themselves to providing qualified and high-efficiency products and offer their optimal services to develop and expand international trade.

Warmly welcome to GUI BAO for business talks and create a better and brighter future! Jiangsu GUI BAO Group Corp. LTD.

LEGAL REP.: ZHOU YUAN LONG

ADD: NO.58, RENMIN ROAD, ZHUTANG TOWNSHIP, Jiangyin, Jiangsu Province TEL:0510-6381883 6381305 FAX:0510-6381242 PC:214415

Main uses: The machine is designed according to the theory of mixing spraying powder. It is widely used in the pulverizing materials in pharmacheutical foodstuff and chemical industries. Structural Features: The machine consists of the pulverizer and the fixed quantity adder it has a simpel struture, easy to detach and clean. Some main parts inside the chamber that bouches the materiats to be pulverized are all made of stainless stell thus the machine is acid-resiting corrosion-resiting.

The machine works properly, smoothly with low noise low cosumption of power.

Working Principle: It is high speed revolving machine, it uses double sharp blades, it cuts so fast that the materials in the chamber is well pulverized by the revloing blades To meet the differ-

Model GFSJ-16 High-Efficient Pulveriger




ent demands. Vatious kind of blades with different shape and size can be applied, the screen meshes can be substituted to get the right granules or adjust the speed.

ULTRA-CRYOGENIC MILL DFJ-2





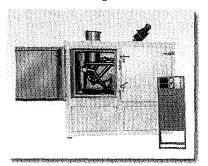

Main uses: It's applied widely to grind material such as pharmaceuticals, foodstuffs, chemical industry, national defence and science research etc. Especially, it's possible to grind thermo-sensitive material, such as rubber, plastics, medicine, animals and plants etc. Which are tenacious,


not easy to grind under normal temperature. Low temperature grinding is that material is cooled by water and grinded to demanded level without destroying original ingredient. It's more efficient to combustibles by means of low temperature grinding.

Struture: Pre-cooling barrel, the spiral coveyor, lower temperature air-blower, whirl-separator, powder filter, electric control machine etc.

Model YG Filling and Sealing Machine for Nourishing Bottles

Main uses: This machine is specially designed for the filling and sealing of the straight tube bottles which contains oral nourishing liquid.


Properties: The machine preforms the following procedures: first converying bottles with a rod, then filling, drop-

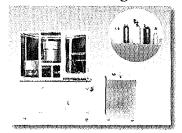

ping bottles covers, then fastening the covers and finally sending bottles out. Themachine applies rod-feeding bottles to avoid lacking of bottles. While fillingthe bottles. A flexible device is fitted in the fastening place to avoid bottle-crashing. A mini decelerating box is equipped to control the speed. The various actions in the filling and sealing is controlled by the convex theory. The machine is compactly designed, easy to handle and maintain.



MODEL WPJ-15 MINI-EFFICIENT PULVERIZER

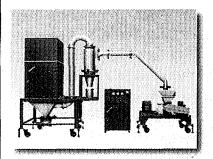


Main uses: The machine is made up of main engine, auxiliary engine and electric box, featuring winnowing without sieve or mesh, even granule size, continuous operation, etc.

The machine reaches advanced world level and is widely used for

material grinding in pharmaceutical, chemical, foodstuff industries, etc. Working Principle: Material goes through the funnel and is sent to the grinding chamber by the spiral conveyor, where it is cut into small bits by the whirlwind separator and bag-type dust remover by merns of negative pressure wind transporation.

ZS-Syst Vibrating Screens


Main uses: The machines are suitable for flow process, and they are perfect equipments used for granulas to screen. Characterstics: The machine cosist of hopper, vibrating chamber, connecting shaft, and electricmotor. There are eccentric wheel, rubber softwate softwate,

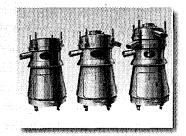

nain shaft and shaft-bearing in the vibrating chamber. When the adjustable eccentric weight is passed the centre of main shaft by motor driven, it may be out of balace and causes eccentricity, and so that the materials form regular eddy. The amplitude of the weight adjuster may regulate by changing the different material and the meshes of screens. It is small in volume, thightness, dutness, noise-free, high output, low in energy consumption, easy to move and easy to maintain.



RXH-Type Oven Series

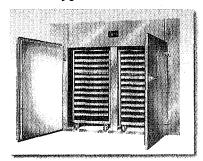


Main uses: It's suitable for phRMcy, chemical, foodstuff, light industry, heavy industry ect. Such as original medicines, raw medicine doses, packing bottles, drying vegetables, foodstaffs, plastic resin, electric pressure, original workpieces, heating lacquer etc.

GFGQ-100Multiple Function

Churning Drier

Main uses: The machine is suitable for drying granulated material of pharmacy, Chemical, foodstuff, light industry and viscous material. With the multiple drying function, it is the advanced equipment of 1990's.

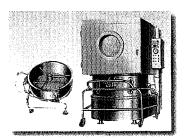

Characteristics:

- 1. The machine makes material dry by hot air churning bed. It feature rapid drying, well distributed temperature. The hot air meet the GMP standard.
- 2. The main drier is supplied in sets of air heater, pilot air pump and automatic temperature controller. The material tank of drier is jacked up by cylinder, It has the characteristic of simple structure, easy operation and sealing guarantee.
- 3. The material tank is equiped with agitator and special air scattering plate, so there is no blind angle in drying bed. It is very convenient to unload wash since adopting the system of bottom turn board unloading.

PURIFICATION EQUIPMENT (SUZHOU)CO,LTD.



Expert in Micro-Pollation Controlling



*Specially design center,data center,judge center of national purification engineering
*Responsible conduct,operating,trainning,maintainace of pure water & purification engineering
*Using modern sophisticated purificate technology & U.S.FDA standard
*Unique whole CAD system of purification engineering
*A lot of famous comsumers & key projects

Provide following service in pharmaceutical, electric, chemical & food area

Pure water & waste water treatment system Clean room engineering system Serving system for purified projects

Add:No.3 Jinshan Road the National New & Hi-Tech Industrial Development Zone Suzhou Tel:0086-512-8256112.8255652 Fax:0086-512-8256112


Post card:215011 General Manager:Ke Ming-ming





BEIJING AEROSPACE CITY MACHINERY ELECTRICEQUIPMENT PLANT BRIEF INTRODUCTION

As the specialised manufacturer of pharmaceutical machinery and machinery electric equipment, Beijing Aerospace City Machinery Electric Equipment plant comprises technique, matured process equipment and reliable quality. Products are sold to over 20 provinces, cities and areas, Furthermore, exported to South East Asia and America regions. Any enquiries and order are cordially welcom

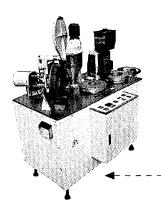


JTJ-III Semi Automatic Capsule Filling Machine

JTJ-III Semi Automatic Capsule Filling Machine, combining electrial and pneumatic control, fitted AC frequency conversion stepless speed converter, and equipped with and electronic automatic counter and trigger button, is reliable to use, sensitive to control, and easy to operate automatically performs capsule positioning, the finished products rate of 97% or more is obtained.

| productivity size of capsule granule/h | Cilling type | power supply kw | weight kg | dinension L *W*H (mm) | subsidiary installation |
|--|-----------------|--------------------|--------------|--------------------------|---|
| 6000~13000 0 ~ 3 | powder, granule | 2.61 | 400 | 1240 × 860 × 1480 | air compresser, circulating vacuum system |

MG-1 Polishing Machine

Advanced technology, reliable product Make use of polish method, the capsule ! and tablet surface shall be clean and smooth Stepless variable speed With the change of main axle inclination angle, the polish time can be adjusted

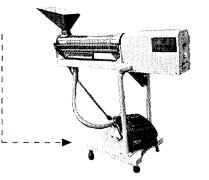

| Technical Data | | | | |
|---------------------------|--------------|------------|-------------------|--------------|
| productivity guanule/h | power supply | power W | size L × W × H | weight Kg |
| 3000 ~ 7000 | AC200V, 50Hz | 180 | 1150 ×1250× 400 | 620 |

Address: No.56A Wu Ke Song Rd. Haidian District, Beijing, 100039, China.


Director: Wang Yong Qi

Tei:(8610)68182941, (8610)68211217 Fax:(8610)68182941


ZHOUKOU PHARMACEUTICAL MACHINERY MANUFACTURER

We are specialized in manufacturing various kinds of slicing machines for the pre-processing of traditional Chinese medicine. What we pursue is fine quality & excellent service. We can supply various kinds of pre-processing machinery for herbal medicine, to meet the requirements of different clients. High quality, high reliability, high efficiency, high safety, & high flexibility have become the ear-make of our machines.

RECIPROCATING SLICING MACHINE WQY 240-2

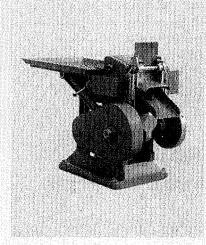

Add:No.180,Bayi Road, Zhoukou, Henan 466001,China ROTATING TYPE SLICING MACHINE OYSERIES

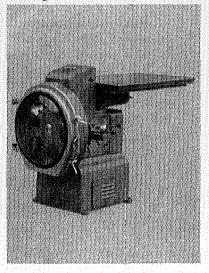


Tel:0086-394-8223481 0086-394-8223480 Fax:0086-394-8232737

ZHOUKOU PHARMACEUTICAL MACHINERY MANUFACTURER


HERBAL MEDICINE WASHING MACHINE

XY 720

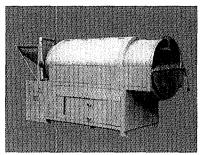

It is suited to wash fruits and other similar goods, widely used in the industries of Pharmacy, Food, Tinned Food, Vegetable Processing, etc. The processing capacity varies from 200~750kg/h. We can also meet the special requirements of the clients.

FRYING MACHINE CY SERIES

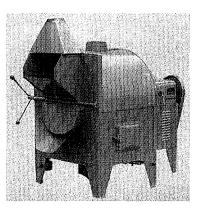

Mainly used in concocting of traditional Chinese medicine, such as: clear frying, sand frying, bran frying, vinegar frying, frying with charcoal, and honey moxibustion etc, as well as in processing and frying of medicinal herbals and fruits.

We are the experts in this circle. Reliable quality is what we offer you.

Add:No.180, Bayi Road, Zhoukou, Henan 466001, China Tel:0086-394-8223481 0086-394-8223480 Fax:0086-394-8232737


TIANJIN CHINESE MEDICINE MACHINERY WORKS

We are the major manufacturer of Chinese medicine machinery in China, and we can supply you with the following machine series. If you have any special requirement, please contact with us, we are always at your service. We adhere to the principle of "Quality first, customer first".


I:CRUSHING UNIT TF-700




III:CHOPPER TYPE MEDICINE CUTTER D74-10

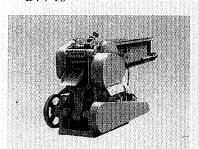

II:THREE ROLLER HONEY PILLS FORMING MACHINE ZTM 80-3

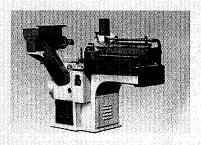

IV:TRIPLE EFFECT CONCENTRATOR TSNQ

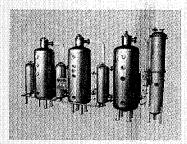

ADD:No.278, Xi Qing Dao, Xi Qing District, Tianjin,300380,China


TEL:86-22-27391201

FAX:86-22-27391201


Add: No.92, Weijin Road TianJin Tel:(022)27405545 Fax:(022) 23350851 P.C.:300072




Director Gong Li Qing

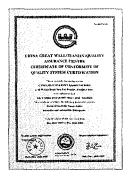

The radio factory and the precision instruments of TianJin University possesses the comprehensive techniue consisting of microdectronics, computer, precision machine, optics and so on. It is an enterprise equipped with high science and technique that produces the drug testing instruments and was awarded the centificat for Quality System of GB/T 19001-ISO 9001:1994 in 1997.

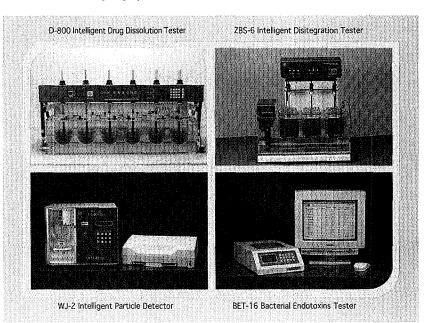


TIANJIN MEDICAL INSTRUMENT INSTITUTE

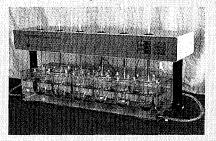



Intellectual Pharmacal Dissolution Testor(RC-610).

Intellectual Disintegration Testor(BJ) Series

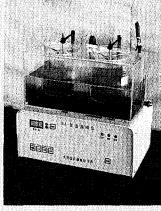

Tablet Hardness Testor (YD-II)

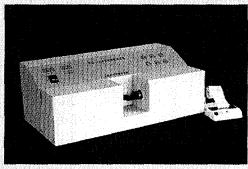

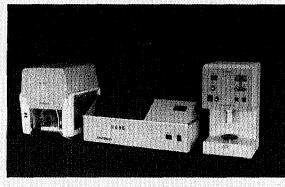

Gel Strength Testor (JS-I)


Add:No.237 Hong Qi South Road Nankai Tianjin China

Tel:(022)23363018 23369102

Fax:(022)23365845

Post Code:300191

Ry Thermocompression Distiller

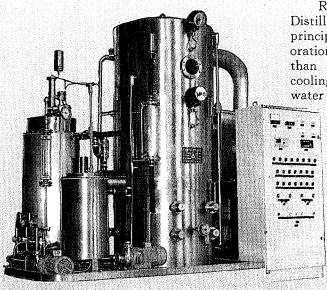

Ry Thermocompression Distiller is based on the principle of heat pump evaporation. It's saved more than 90% energy without cooling water. The distiller water is high quality, low

cost and stable water quality. It conforms to the standard set by the People's Republic of China for use in medical injection. It's also suitable for use in chemical, petroleum, food, desalination of sea water, medical treatment and light industry.

Main Products

- 1)Air Lift Ferment Reactor
- 2)Thermocopression Distiller
- 3)Horizontal Closed Peeler Centrifuge 4)Doubble—cone Rotary Vacuum Dry-
- er 5)Oscillating Granulator
- 6)Multi-effect Water Still
- 7) Glass lined Closed Reactor
- 8)Mixer

- 9)Powder Forming Machine
- 10)Sugar Coatiog Machine
- 11) Honey Refining Unit
- 12) Non standard Chemical Equipment
- 13)Motion Crystallizing Machine
- 14) Spraying Drying Tower
- 15) Auto Ion Exchange Water Treatment Machine

Brief Introduction

Baoji Pharmaceutical Machinery Works is one of the key manufactures of the State Medicine Administrative Bureau, specializing in producing pharmaceutical and chemical equipment. It has 30 years kindly service experience for state and abroad user in pharmaceutical industry. It also gets good name from user.

The works has chief production shops of riveting and welding glass lining forging machine and heat treatment. It has a backbone of rich experience, skilled business and excellent technology. It has perfect system of quality administration and advanced test means. The works research institute is in contact with the enterprises of the pharmaceutical, chemical and college of the whole country. It devotes itself to the development of pharmary chemical machines and equipments new products and technique. It has put a lot of high technology products on the market.

Director: Li Guitang Post Code: 721006 Add: 21 Chuan Shaan Road, Baoji. Shaanxi. Tel: 0917—3315195 Fax: Cable: 2750 Beijing Pharmaceutical Machinery and Equipment Factory Add: Dajiaoting Guang Qu Street, Beijing, China

Tel:0086-10-67716779

Fax:67721941

P.C.:10022

Beijing Chinese Medicine Machinery Works

Add: 33 Totiao Chang Xiang, Chong Wen Qu, Beijing, China Tel: 0086-10-65115345 Fax: 65115343 P.C.: 100051

Beijing Tong-Da Hi-Tech Company

Add: No. 52 Yongding Road, Haidian District, Beijing, China

Tel: 0086-10-68388767 Fax: 68388652 P.C.: 100854

Beijing Tian Min Hi-Tech Development Co.

Add: P.O. Box: 9200-71 Beijing, China.

Tel: 0086-10-68382704 Fax: 68381851 P.C.: 100076

Beijing Hong Hua Electro—mechanical New Technology Company

Add: 50 Yongding road, Haidian District, Beijing China

Tel:0086-10-68385234 Fax:68386763 P.C.:100854

Beijing Aeronautical Manufacturing Technology Research Institute

Add: P.O. Box: 863 Bejing, China

Tel: 0086-10-65797893 Fax: 65762306 P.C.: 100024

Beijing Dasong Machinery & Electronics Techology Co.

Add: No. 118 Yimuyuan, Haidian District, Beijing, China

Tel:0086-10-62882351 Fax:62881498 P.C.:100091

Yangchunhui Science & Technology Developing Co., Ltd Add: No 73 Shuang—Quan—Pu, Hai Dian district, Beijing, China

Tel:0086-10-62329131

Fax: 62329131

P.C.:100085

Beijing Aerospace City Machinery Electric Equipment Plant Add: No. 56A Wu Ke Song Rd. Haidian District, Beijing, China.

Tel: 0086-10-68182941

Fax: 68182941 P. C.: 100039

Beijing Tianli Medicine Equipment Co., Ltd

Add: Rm. 311 Towercrest Plaza Maizidian West Rd., Beijing, China

Tel: 0086-10-64671989

Fax: 64671986

P.C.:100016

Tianjin Chinese Medicine Machinery Works

Add: No. 278, Xi Qing Dao, Xi Qing District, Tianjin, China

Tel:0086-22-27391201

Fax: 27391201

P.C.:300380

Tianjin Medical Instrument Institute

Add: No. 237 Hong Qi South Road, Nankai District, Tianjin, China

Tel:0086-22-23363018

Fax:23365845

P.C.:300191

Tianjin University Radio Factory

Tianjin University Precision Instruments Factory

Add: No. 92, Weijin Road, Tianjin, China

Tel:0086-22-27405545

Fax: 23350851

P.C.:300072

Liaoyang Pharmaceutical Machinery Co., Ltd

Add: No. Sec1, Shengli Rd, Liaoyang, Liaoning, China

Tel: 0086-419-2126405 Fax: 2158106

P.C.:111000

Jinzhou Oldtrees Packing Machinery Co., Ltd

Add: No. 63 Chongqing Rd, Guta District, Jinzhou, Liaoning, China

Tel:0086-416-4167750

Fax:4168359

P. C. :121001

Jin Wan Co., LTD

Add: 559 Hamatang Rd, Dandong, Liaoning, China

Tel: 0086 — 415 — 4151322

Fax:4151380

P. C. :118003

Harbin Aircraft Manufacturing Corporation

Add: 15 Youxie Street, Pingfang District, Harbin, Heilongjiang, China

Tel: 0086 - 451 - 6501726

Fax:6502223

P. C. 150060

Heilongjiang Deer Pharmaceutical Machinery Co., Ltd

Add: No. 37, Linhai Rd, Jiamusi, Heilongjiang, China

Tel: 0086 - 454 - 8782182

Fax:8782365

P. C: 154007

Shanghai Yuan Dong Pharmaceutical Machinery General Factory

Add: 359, Kang Ding Rd, Shanghai, China

Tel:0086-21-62170175

Fax: 62564739

P. C. :200040

Shanghai Tianxiang & Chentai Pharmaceutical Machinery Co., Ltd

Add: 1380, Kang Ding Rd, Shanghai, China

Tel:0086-21-62583414

Fax: 62582067

P.C.:200042

Shanghai Tian He Pharmaceutical Machinery Factory Add: Chengjia Town, Chongming County, Shanghai, China

Tel:0086-21-59431285

Fax:59432024

P. C. :202162

Shanghai Ying Hua Pharmaceutical Machine Industrid Company Add: No. 8 Passage 510 San Men Lu, Shanghai, China

Tel:0086-21-65366089

Fax: 65444114 P.C.: 200434

Shanghai Jin Qiao Spray and Dry Complete Equipment Plant Add: 60 Ying Jia Bang Yang Gao Rd, Pu Dong, Shanghai, China Tel: 0086-21-58545772 Fax: 58995886 P.C.: 201206

Shanghai Sanli Pharmaceutical Equipment Technical Corporation Add: 900 Yan An Rd., Shanghai, China

Tel: 0086-21-62131812 Fax: 62131812 P.C.: 200052

Shanghai Xin—Xing Pharmaceutical Machinery Plant Add:Roome • 501' House • 8' Lane 561 Yan Chang Rd • (Middle) Shanghai, China

Tel:0086-21-56770079 Fax:56770079 P.C.:200072

Huxi Iron Works

Add: 168/185 Pingtang Rd, Shanghai, China Tel: 0086-21-62394624 Fax: 62394624 P.C.: 200335

Jiangsu Guibao Group Corp. Ltd.

Add: No. 58 Renmin Rd, Zhutang Township, Jiangyin, Jiangsu, China

Tel: 0080-510-6381271 Fax: 6381242 P.C.: 214415

Suzhou Purification Equipment Co. Ltd.

Add: No. 3 Jinshan Rd, The National New & Hi—Tech Industrial Development Zone, Suzhou, China

Tel: 0086-512-8255652 Fax: 8256112 P.C.: 215011

Suzhou Industrial PARK Huayuan Machinery & Electronic Appliance Co., Ltd.

Add: No. 8 Nan Yuan Rd, Suzhou, China

Tel: 0086-512-5308737 Fax: 5107028 P.C.: 215006

Haimen Yonghe Machinery Co., Ltd.

Add: Haimen Cotton Plant (North of Haimen Harbour), Jiangsu, China

Tel:0086-513-2292008

Fax: 2292308

P.C.:226114

Changzhou One-Step Drying Equipment Factory

Add: Jiaoxi, Changzhou, Jiangsu, China

Tel:0086-451-8902618

Fax:8902818

P.C.:213116

Changshu Pharmaceutical Machinery Factory

Add No. 38 Yushan North Rd, Changshu, Jiangsu, China

Tel:0086-520-2895982

Fax:2897000

P.C. :215500

Jiangyin Freeze Drying Equipment Plant

Add: West Street Shen Gang Zhen, Jiangyin, Jiangsu, China

Tel:0086-510-6621481

Fax:6621001

P. C. :214443

Wenzhou Fangzheng Enterprise Group Co.

(Wenzhou Packaging Machinery General Factory)

Add: 33 Feixia South Rd, Wenzhou, Zhejiang, China.

Tel:0086-577-8631349

Fax:8632224

P.C.:325003

Wenzhou Pharmaceutical Machinery Factory

Add, No. 106, Xueshan Rd, Wenzhou, Zhejiang, China

Tel:0086-577-8524295

Fax: 8525448

P.C.:325005

Ruian Pharmaceutical Machinery Factory Zhejiang China Add: 20 yanjiangxincun, Ruian, Zhejiang, China

Tel: 0086 - 577 - 5623309

Fax: 5622267

P.C.:325200

Wenzhou Pharmaceutical Equipment Factory

Add: Shanxiahou Lane Liming Road East, Wenzhou, Zhejiang,

China

Tel:0086-577-8317378

Fax:8362420

P.C.:325003

Jiangnan Packing Machinery Co., Ltd. Ruian Zhejiang Add: Hongqi Industrial Area, Ruian, Zhejiang, China Tel: 0086-577-5662604 Fax: 5663983 P.C.: 325200

Ruian Feiyun Machinery Factory(Co)

Add: Dongshan Development Zone, Ruian, Zhejiang, China
Tel: 0086-577-5603082 Fax: 5603083 P.C.: 325200

Wenzhou Zhongzhi Medical Equipment Factory Add:Shipu Ind. Zone, Yongzhong, 559 Airport Ave, Wenzhou, China

Tel:0086-577-6373357

Fax: 6876357

P.C.:325024

Zhejiang Wenzhou Yaguang Machinery Manufacturing Co., Ltd. Add: 99 Canghe Rd. (W). Yongzhong Town, Wenzhou, Zhejiang, China

Tel:0086-577-6372550

Fax:6373389

P.C.:325024

Shandong Xinhua Pharmaceutucal Co., Ltd.(Former Shandong Xinhua Pharmaceutucal Factory)Add: No. 25 Xinhua Rd, Zhangdian district, Zibo, Shandong, China

Tel:0086-533-2185538

Fax:3286157

P.C.:255005

Shandong Xinhua Medical Instrument Group Co., Ltd Add: Honggou Road, Zhang Dian District, Zibo, Shandong, China Tel: 0086-533-2184349 Fax: 2184349 P.C.: 255027 Shandong Peng Lai Pharmaceutical Machiney Factory

Add: Peng Lai, Shandong, China

Tel:0086-535-5971223 Fax:5971223

P.C.:265609

Zhoukou Pharmaceutical Machinery Manufacturer

Add: No. 180 Bayi Rd, Zhoukou, Henan, China

Tel:0086-394-8223481

Fax:8232737

P.C.:466001

China Wuhan Pharmaceutical Machinery Factory

Add: 25 Hangxiang (Red Lane), Wuchang District, Wuhan, China.

Tel:0086-27-8844815

Fax:8874935

P.C.:430061

Wuhan Medical Apparatuses and Instruments Plant

Add: NO. 1 Xincun, Sanjiao Rd, Wuchang District, Wuhan, China

Tel:0086-27-6811350

Fax:6811258

P. C. 430062

Zhongnan No. 1 Pharmaceutical Machinery

Add Qiyang, Hunan, China

Tel:0086-746-3222388

Fax: 3222617

P.C.:426110

Zhongnan Pharmaceutical Machinery Factory No. 2

Add: Xianjia Hu Road, Changsha, Hunan, China

Tel:0086-731-8809343

Fax:8809250

P.C.:410006

Guangdong Prov. Huiyang Machinery Factory

Add: 26. Xia jiao middle road, Hui Zhou, Guang Dong, China

Tel:0086-752-2227000

Fax:2241409

P. C. 516001

Chongqing Pharmaceutical Machinery Factory
Add: No. 35 Shapingba Main Street, Chongqing, China
Tel: 0086-23-65312352 Fax: 65315646 P. C. 400030

Long March Pharmacertical Machinery
Sichuan Long March Pharmacetical Co. Ltd.
Add:120 Baiyang Road, Leshan, Sichuan China
Tel:0086-833-2440059 Fax:2440059 P.C.:614006

Chongqing Fine and Elaborate Pharmaceutical Machinery Company, Ltd

Add:Rm 6-2 SCI-Tech Building High/New Tech Devel,
Zone50 Yuzhou Rd, Chongqing, China
Tel:0086-23-68604229 Fax:68604229 P.C.:400039

Chongqing Guangxia Drying Equipment & Engineering Co. Add:128# Shuangfeng Rd, Yubei Area, Chongqing, China Tel:0086-23-67823244, Fax:67822860 P.C.:401120

Baoji Pharmaceutical Machinery Works
Add: 21 Chuan Shan Rd. Baoji. Shaanxi, China
Tel: 0086-917-3312691 Fax: 3311653 P.C.: 721006

Baoji Jian Hua Phmarceutical Machinery Co. Ltd Add: No. 75 Jianhua Rd, Baoji, Shaanxi, China Tel: 0086-917-3554670 Fax: 3554670 P.C.: 721001

| 单位名称 | 通讯地址 | 邮编 | 电话 | 传真 |
|---------------|---------------------------------|--------|--------------|----------|
| 北京制药机械设备厂 | 北京广渠门外大郊亭北站 | 100022 | 010-67716779 | 67721941 |
| 北京市中药机械厂 | 北京崇文区长巷头条 33 号 | 100051 | 010-65115345 | 65115343 |
| 北京市通大高技术公司 | 北京 142 信箱 403 分箱 | 100854 | 010-68388767 | 68388652 |
| 北京天民高科技开发公司 | 北京 9200 信箱 71 分箱 | 100076 | 010-68382704 | 68381851 |
| 北京弘华机电新技术公司 | 北京市 142 信箱 208 分箱 | 100854 | 010-68385234 | 68213632 |
| 北京航空工艺研究所 | 北京市 340 信箱技术开发处 | 100024 | 010-65797893 | 65762306 |
| 北京大松机械电子总公司 | 北京海淀区一亩园 118 号 | 100091 | 010-62882361 | 62881498 |
| 北京阳春汇科技发展有限公司 | 北京海淀区双泉堡 73 号 | 100085 | 010-62329131 | 62329131 |
| 北京市航天成机电设备厂 | 北京海淀区五棵松路 甲 56 号 | 100039 | 010-68182941 | 68182941 |
| 北京天利医药设备有限公司 | 北京朝阳区麦子店西路 3 号 新恒基国际大厦 311 室 | 100016 | 010-64671989 | 64671986 |
| 天津市中药机械厂 | 天津市西青区西青道 278 号 | 300380 | 022-27391201 | 27391201 |
| 天津市医疗器械研究所 | 天津市南开区红旗南路 237 号 | 300191 | 022-23363018 | 23365845 |
| 天津大学无线电厂 | 天津市南开区卫津路 92 号 | 300072 | 022-27404928 | 23350851 |
| 辽阳制药机械股份有限公司 | 辽宁省辽阳市胜利路 一段二号 | 111000 | 0419-2126337 | 2126413 |
| 锦州欧仕包装机械有限公司 | 辽宁省锦州市重庆路 63 号 | 121001 | 0416-4167750 | 4168359 |
| 丹东金丸公司 | 丹东市振安区哈蟆塘大街 559 号 | 118003 | 0415-4151322 | 4151380 |

| 単 位 名 称 | 通讯地址 | 邮编 | 电话 | 传真 |
|----------------------|--------------------------|--------|----------------|----------|
| 哈尔滨飞机制造公司 | 哈尔滨市平房区友协街 15 号 | 150066 | 0451 — 6501726 | 6502223 |
| 黑龙江迪尔制药机械有限公司 | 黑龙江省佳木斯市林海路 37 号 | 154007 | 0454 — 8782182 | 8782365 |
| 上海远东制药机械总厂 | 上海康定路 359 号 | 200040 | 021-62170175 | 62564739 |
| 上海天祥. 健台制药机械 有限公司 | 上海康定路 1380 号 | 200042 | 021 - 62583414 | 62583648 |
| 上海天和制药机械厂 | 上海崇明县陈家镇南首 | 202162 | 021-59431285 | 59432024 |
| 上海英华实业有限公司 | 上海市三门路 510 弄 8 号 | 200434 | 021 - 65366089 | 65444114 |
| 上海金桥喷雾干燥成套设备厂 | 上海浦东扬高路硬家浜 60 号 | 201206 | 021 — 58545772 | 58995886 |
| 上海三利制药装备技术公司 | 上海市延安西路 900 号 301 室 | 200052 | 021 - 62131812 | 62131812 |
| 上海新星制药机械厂 | 上海延昌路 561 寿 8 号 501 室 | 200072 | 021-56770079 | 56770079 |
| 上海房屋设备总公司 沪西工程处 | 上海市平塘路 168/185 号 | 200335 | 021 - 62394624 | 62394624 |
| 江苏瑰宝集团公司 | 江苏江阴祝塘镇人民路 58号 | 214415 | 0510 - 6381271 | 6381242 |
| 苏州普滤得净化有限公司 | 苏州高技术开发区金山路 3号 | 215011 | 0510-8255652 | 8256112 |
| 苏州工业园区华园机电 有限公司 | 江苏苏州市南圆路 8 号 | 215006 | 05125308737 | 5107028 |
| 海门永和机械有限公司 | 江苏海门棉种场 (海门港北首) | 226114 | 0513-2292008 | 2292308 |

| 单位名称 | 通讯地址 | 邮编 | 电 话 | 传真 |
|--------------------|---------------------|--------|-------------------------------|---------|
| 常州市一步干燥设备厂 | 江苏常州市焦溪镇 | 213116 | 0519-8902618 | 8902818 |
| 国营常熟制药机械厂 | 江苏省常熟市虞山北路 38 号 | 215500 | 0520 — 2895982 | 2897000 |
| 江阴市冷冻干燥设备厂 | 江苏江阴市申港镇西街 | 214443 | 0510-6621481 | 6621001 |
| 温州方正企业集团公司 | 浙江省温州市飞霞南路 33号 | 325003 | 0577 — 8631349 | 8632224 |
| 温州市制药机械厂 | 温州市雪山路 106 号 | 325005 | 0577 — 8524295 | 8522123 |
| 浙江瑞安市制药机械厂 | 浙江省瑞安市沿江新村 20号 | 325200 | 0577 - 5623309 | 5622267 |
| 温州市制药设备厂 | 温州市黎明东路山下后巷 | 325003 | 0577 — 8317378 | 8362420 |
| 江南包装机械有限公司 | 浙江瑞安市红旗工业区 | 325200 | 05 77 — 5663701 | 5663983 |
| 瑞安飞云机械实业公司 | 瑞安东山开发区 (东山邮局对面) | 325200 | 0577 5603082 | 5603083 |
| 温州市中药机械厂 | 温州永中中河路河东新村 15 号 | 325024 | 0577 — 6371414 | 6876357 |
| 温州亚光机械制造有限公司 | 温州市永中镇沧河西路 99 号 | 325024 | 0577-6372550 | 6373389 |
| 新华制药股份有限公司机械厂 | 山东淄博张店新华街 25 号 | 255005 | 0533-2185538 | 3286157 |
| 山东新华医疗器械集团 有限公司 | 山东淄博市张店区洪沟路 2号 | 255027 | 0533-2184349 | 2184349 |
| 山东省蓬莱制药机械厂 | 山东蓬莱市马格庄镇 | 265609 | 0535-5971223 | 5971223 |
| 河南省周口制药机械厂 | 河南周口市八一路北段 180号 | 466001 | 0394-8223481 | 8232737 |

| 单位名称 | 通讯地址 | 邮编 | 电 话 | 传真 |
|--------------------|------------------------|--------|------------------------|----------|
| 武汉制药机械厂 | 武汉市武昌区红巷9号 | 430061 | 027 - 8844815 | 8874935 |
| 武汉医疗器械厂 | 武汉市武昌区三角路新村 1号 | 430062 | 027-6811250 | 6811258 |
| 中南制药机械一厂 | 湖南祁阳县浯溪镇天马路 1号 | 426100 | 0746 — 3222388 | 3222617 |
| 中南制药机械二厂 | 湖南长沙市河西咸嘉湖 | 410013 | 0731 — 8809536 | 8809250 |
| 广东省惠阳机械厂 | 广东省惠州市下角中路 26 号 | 516001 | 0752 — 2227 000 | 2241409 |
| 重庆制药机械厂 | 重庆市沙坪坝正街 35 号 | 400030 | 023-65312352 | 65315646 |
| 四川长征制药股份有限公司制药机械分厂 | 四川乐山市柏杨路 120 号 | 614006 | 0833-2440059 | 2440059 |
| 重庆精工制药机械 有限责任公司 | 重庆输州路 50 号 科创楼 6-2# | 630039 | 023-68604229 | 68604229 |
| 重庆广厦干燥设备工程公司 | 重庆渝北区双凤路 128 号 | 401120 | 023 - 67823244 | 67822860 |
| 宝鸡制药机械厂 | 陕西省宝鸡市川陕路 21 号 . | 721006 | 0917 - 3312691 | 3311653 |
| 宝鸡建华制药机械有限公司 | 陕西省宝鸡市建华路 75 号 | 721001 | 0917 — 3554670 | 3554670 |
| | | | | |
| | | | | |
| | | | | |
| | | - | | |
| | <u> </u> | | j | |

POSTSCRIPT

«A Handbook of Chinese Pharmaceutical Machinery & Equipment for Export is organized by Department of International Cooperation, the State Pharmaceutical Administration of China (SPAC) and edited by the China Association for Pharmaceutical Equipment under the requirement of "Promotion of Cooperation Partnership Among Selected Developing Countries in Pharmaceutical Industry" (US/UT/RAS/97/138), a joint project between UNIDO and SPAC, to make people a better understanding of Chinese pharmaceutical machinery and equipment and enhance the cooperation between developing countries.

The Handbook is in the format of product specification to provide easy and quick search for product name, model, application, working mechanism, main technical data, fitted motor, overall dimensions etc.

All products listed within over one hundred pages of the Handbook are classified and listed in line with "Terms of Pharmaceutical Machinery" under China National Standard GB/T15692. 1~9-95, including machinery and equipment for bulk pharmaceuticals, preparation machinery, pharmaceutical milling machinery, pharmaceutical pulverizing machinery, water treatment equipment for pharmaceutical use, pharmaceutical packaging machinery, medicine detecting instrument and others.

The Handbook gathered products with the features of Chinese characteristics, high efficiency and energy—saving, technical advanced and incorporated electricity, optics and mechanics either single unit or automatic line. The unit system applied in the Handbook complies with China National Standard GB $3100\sim3102-92$.

The models listed in the Handbook are only typical specifications, other more specifications are set in the last item: "Other Specifications" owing to space limitation. The color pictures in the handbook are product photos, which is good reference to help understanding the content in the specifications and to see the appearance of products.

The appendix I, for contact needs, is a directory of major manufacturers of pharmaceutical machinery & equipment.

Many thanks should go to all those who have helped us to build this handbook. Particulally, thanks to major editors: Zhang Rihua, Shi Qing, Gao Chuan; editors: Fu Xinan, Wu Chunlin, Zhu Dalin, Qian Changsheng, Li Yurong, Luo Yuling and Yuan Mingya, and thanks to Huang Jiye, Chang Wenzuo, Chen Xingyu Liu Ai who made contributions to this handbook; thanks to Zhu Yingfa, Yuan Junan Xue Yiming Zhang Bin, Yuan Mingya and Ma Jingbo reviewed the English version of the handbook.

Department of International Cooperation, SPAC China Association for Pharmaceutical Equipment Beijing May, 1998

编者的话

《中国制药装备出口产品手册》(以下简称手册),是根据联合国工业发展组织(UNIDO)与国家医药管理局的合作项目一"推动发展中国家制药工业及装备合作研讨会"(US/UT/RAS/97/138)的要求,由国家医药管理局国际合作司组织,中国制药装备行业协会承办编辑而成,以便于人们更好了解中国制药装备的品种、水平,增进发展中国家的合作。

"手册"以产品介绍表的形式,便捷查找所需要的产品名称、型号规格、用途、工作原理、主要技术参数、配套电机、外型尺寸等。

按照中国国家标准 GB/T15692. 1~9--1995《制药机械名词术语》,对 100 多页的产品介绍表,加以分类排序。有原料药设备及机械、制剂机械、药用粉碎机械、饮片(中药)机械、制药用水设备、药品包装机械、药物检测设备和其他制药机械及设备等共八大类产品。

"手册"中,收集了具有中国特色的产品、高效节能的产品、国际先进水平的产品、机电光一体化的产品、符合 GMP 要求的产品等等,既有先进的单机,也有全自动化生产线。

产品介绍表,尽量采用中国国家 GB3100~3102-92 标准《量和单位》规定的计量单位。

由于版面、篇幅所限,产品介绍表的产品型号中,所列出的规格,是典型举例的规格。还有大量的规格,录入该表的最后一项"其他规格"栏内。

为替代产品样本,彩色插页附有该公司产品照片,可对照(配

合)产品介绍表中的产品型号规格,达到既方便看出表中内容,又可明快分析插页中设备外观照片;

为便于需方(制药厂)与供方(中国的部分制造厂)联系,"手册"提供了附录一:中国部分制药机械厂通讯录。

"手册"得到了有关领导、会员单位的鼎力相助,并请张日华、 石青、高川先生担任主要编审;傅新安、吴春林、朱大林、钱长生、李 玉荣、骆玉伶和袁明亚等女士/先生担任编辑;特邀朱英发、原俊 安、薛以明、张斌、袁明亚和马京波女士/先生为英文核稿,他们在 百忙中,拔冗惠于支持,借"手册"面世之际,致以衷心谢意。

> 国家医药管理局国际合作司 中国制药装备行业协会 一九九八年五月 北京